

# **Basic Photoshop For Electron Microscopy**

**(Mostly) Simple Techniques to Tune-Up Your Electron  
Micrographs with Adobe Photoshop**

**Authored by Eric Jay Miller**

© Eric Jay Miller 2013-2022

## **2022 Edition**

I offer my heartfelt thanks to Professor Vinayak Dravid, Ph.D., and Ben Meyers, Ph.D., for encouraging my work on this book while employed at Northwestern University. It would likely not have been completed otherwise. Most of the images in this book were taken on instruments at the *NUANCE* center at Northwestern University.

## Contents

	Introduction	2
<b>I.</b>	Levels	3
<b>II.</b>	Resizing Images	7
<b>III.</b>	Simple Pseudo-Color – Color Table	10
<b>IV.</b>	Simple Pseudo-Color – Gradient Map	14
<b>V.</b>	Complex Pseudo-Color – Color Table	20
*	Addendum to Section IV	34
<b>VI.</b>	Pseudo-Color – The Easy Way	42
<b>VII.</b>	Calculations	50
<b>VIII.</b>	Combining EDS Color Maps – Cut and Paste	54
<b>IX.</b>	Multi-Detector Color	61
*	Tips & Tricks for Multi-Detector Color	73
<b>X.</b>	Finishing Moves for Your Color Images	86
<b>XI.</b>	Drawing a Mag Line Marker (Pre-Existing Mag Line)	92
<b>XII.</b>	Drawing a Mag Line Marker (The <i>SUPER</i> Cheating Way)	97
<b>Appendix I</b>	File Formats	101
<b>Appendix II</b>	The Short Versions	102

## **Introduction**

This paper will discuss several simple procedures that can help to punch your micrographs up a notch. These instructions are written in such a way that just about anybody should be able to pick them and start using them without any previous experience with Photoshop (although that certainly won't hurt either).

This is not intended to be (nor does the author claim it to be) the end-all and be-all of Photoshop instructions. There are often multiple ways of doing the same operation in Photoshop. These are just a few that work pretty well for this particular application of electron micrographs. There might be other procedures that work better for specific applications (or work just plain better). The reason for the omission of any of those techniques is due to the author not having been smart enough to have thought of them prior to writing this.

The program used for this paper was Photoshop CS4 (aka the **full** version), though most older versions of Photoshop will work perfectly fine with most of these techniques. Lighter versions of Photoshop, like Photoshop Elements, will be unable to perform some of these operations.

In electron microscopy there are often concerns about the authenticity of an image once it has been processed through Photoshop. These can be valid concerns as there are a lot of really crazy things one can do with this high-powered program. The main purpose with the techniques discussed here is to either a) make the information in the image easier to see, thus yielding more usable information, or b) to make the image more aesthetically pleasing. These instructions are geared towards image processing and **not image analysis**.


## I. Levels

*The most basic thing* one can do to improve the overall look of a micrograph is to adjust the Brightness and Contrast levels. With Adobe Photoshop, this is a simple, quick and very accurate procedure.

**1.** Open your image in Photoshop.  
Use:

**File > Open**

or **Ctrl + O**

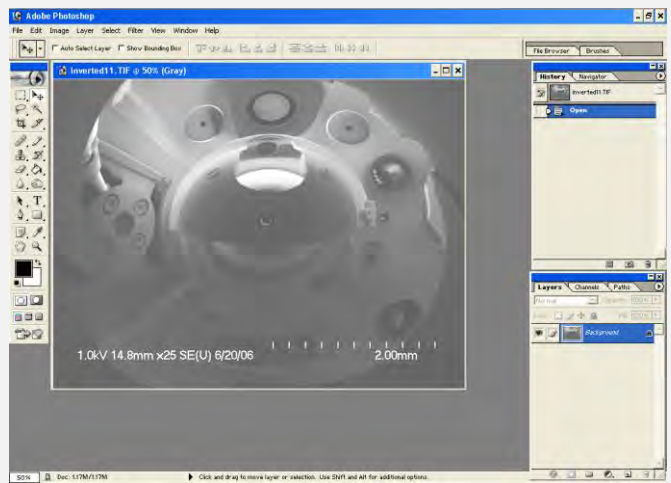

**2.** Make sure that your image is an 8 Bit Grayscale image by clicking on

**Image > Mode**

**Grayscale** and **8 Bits/Channel** should have checks next to them.




**3.** An ideal grayscale image has an area that is absolute black, an area that is absolute white and has the maximum shades of gray in between as possible. Even if the image we took on the microscope is not ideal, we can usually even it out, at least a little.


Open the Levels dialog box by clicking on:

**Image > Adjustments > Levels**


or **Ctrl + L**

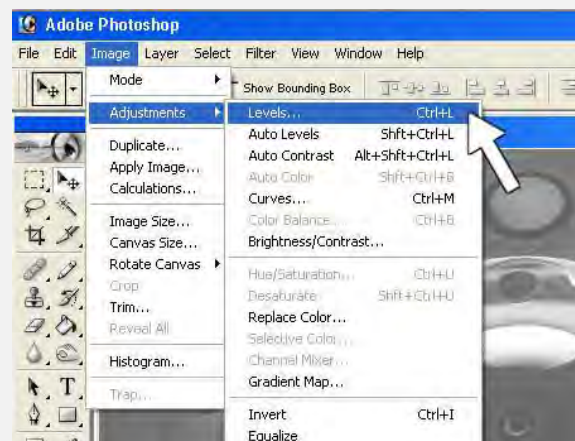


**4.** In the middle of the Levels box there will be a histogram that shows how the 255 shades of gray in your grayscale image are being used. The extreme left of the histogram represents black while the right side represents white. Underneath the histogram are 3 small arrows. These are set points for absolute black, absolute white and the shade of gray in between the two.

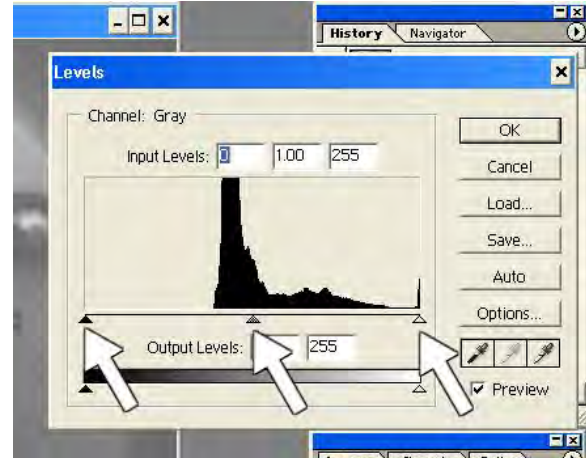


**5.** What you will want to do is bring the black triangle and the white triangle in so that there is no 'empty space'. In the example, there is an empty space in between the main body of the histogram and a peak on the right end. Sometimes a peak on the fringes can be omitted by bringing the triangle past the peak to the main body. But in this case, we will leave it.

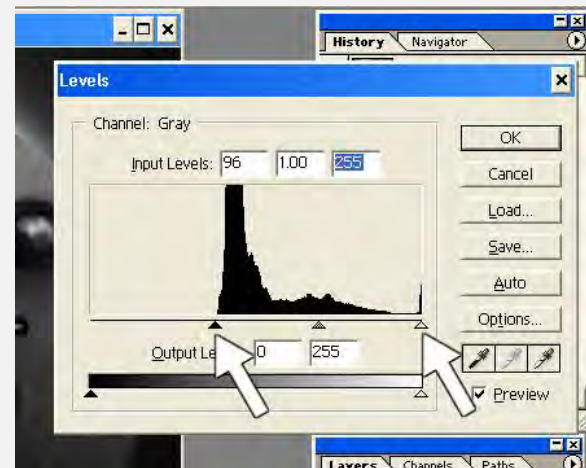

**6.** Next we will set the midpoint. Make sure that the **Preview** box (which is in the lower right corner of the Levels box) is checked. Drag the gray triangle back and forth until there is a good balance between light and dark in your image and click **OK**.

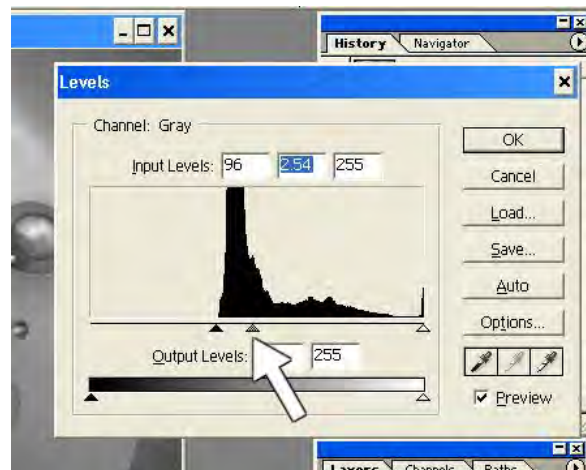

**7.** If we go back and take a look at the Levels box now we can see that Photoshop has stretched the histogram so that it now uses almost every shade of gray, whereas before it was only using about half. Also, the darkest color in the unprocessed image was almost halfway to absolute white and now it is now close to absolute black. This will add a lot of extra definition and visual punch to the image.

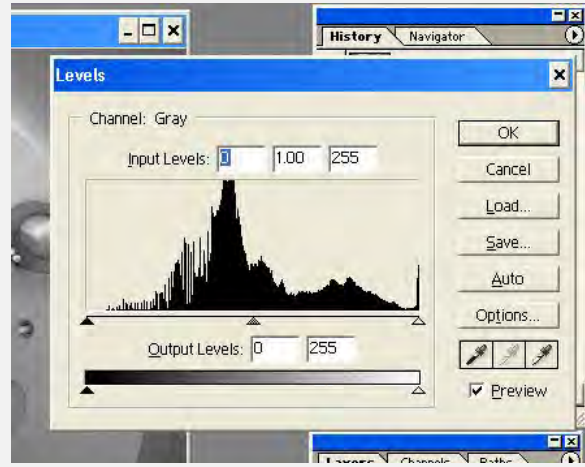


**NOTE:** By saving your image, you'll make these changes we made **PERMANENT!** You *might* not want to do that. A good idea is to always keep your original, unprocessed image, in case you want to do something different to it later on. So once we have processed an image, perform **Save As** to save your file as a new file.


**File > Save As...**

Or **Shift + Ctrl + S**

Then when you close your original file, PS will ask you if you want to save it, SAY **NO!** You should end up with 2 files, one is the original, the second is the one you did stuff to.

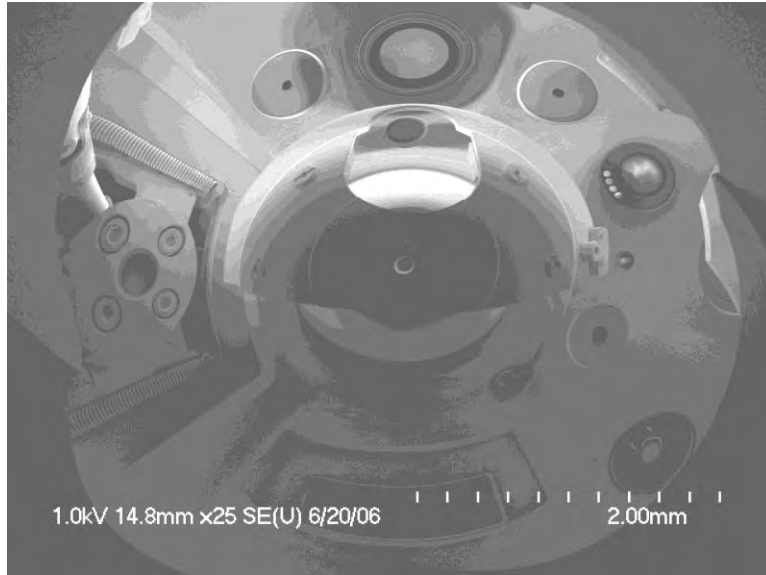

Here is the image before processing. It is not horrible, but it is a little dull.

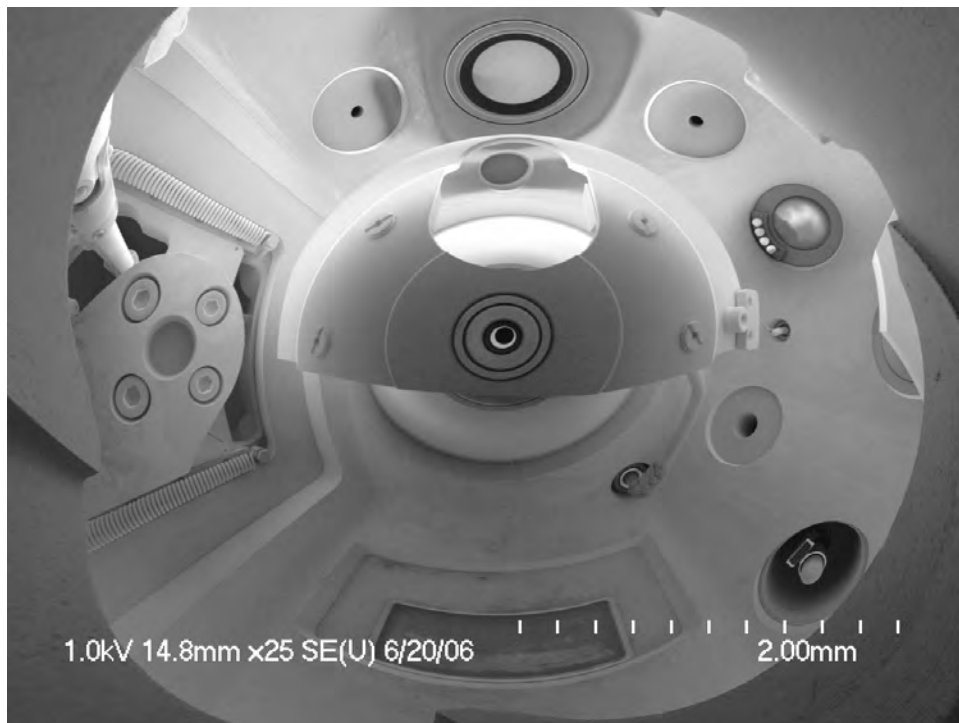


This is the image after the simple Levels procedure. It now has much more definition and is much snappier and appealing to the eye.



## II. Resizing Images

*Do not think* that clicking the corner of your image and dragging it around in Powerpoint or Word to make your image bigger or smaller is the same thing as resizing the image in Photoshop. It is *totally* not. Here we will discuss the right way of resizing an image.

**1.** Do all processing to your image *first*. Levels, colorizing, whatever. Resizing the image will be **THE VERY LAST** thing you will do.

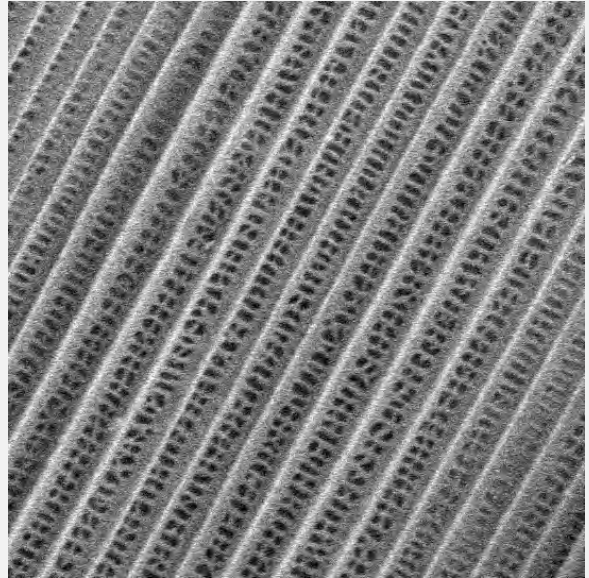

### **IMPORTANT: Do's and Don'ts of Resizing an Image**

**Rule 1:** When we resize images, we want to do it evenly if we can. Ex. Take an image that is 1000x1000 pixels and enlarge it to 2000x2000 or reduce it to 500x500. It's best if you can cut things evenly in halves (or quarters, eighths, etc.) or double them (or quadruple, etc.).


However, the final image size we need is not always going to correspond to such sizes. At some point we'll have to enlarge or reduce an image to some weird size, like 81% of its original size, or enlarge a 1000x1000 pixel image to 2365x2365. So if we cannot abide by Rule 1, we move to:

**Rule 2:** Don't ever resize an image twice. Do it once and you're done.

**2.** To resize our image, we'll open the **Image Size** window by clicking:


**Image > Image Size**

or **Alt + Ctrl + I**

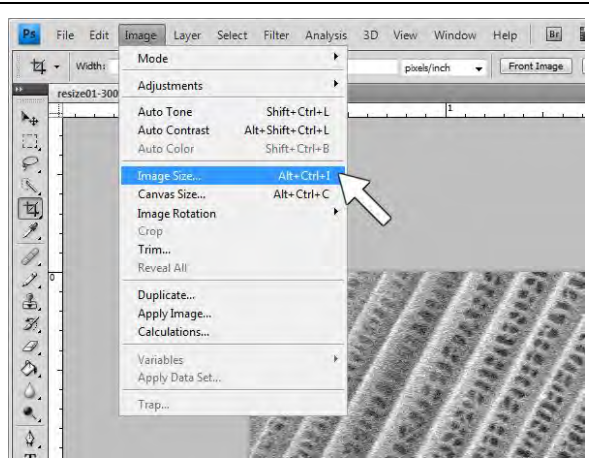

**3.** Once in the Image Size window, we have 2 main options we want to look at. First is the **amount** we will resize the image and also **how** we will resample the image.


We should know the amount we're going to resize the image, and it's usually some set amount based on what we're going to do with the image. But the method we will use to resample the image will depend on the image itself and whether we're reducing or enlarging it.

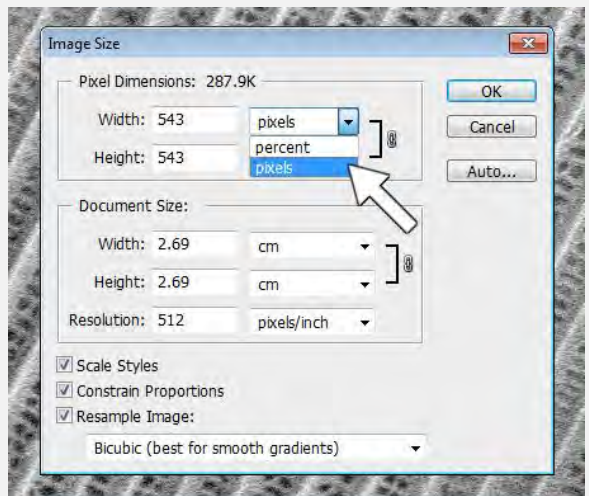







**4.** When we click the resampling dropdown menu, we have 5 options. Some specify they're good for enlarging or reducing. While that may be generally true, it might not be true for your image.

This can be complicated explaining what all these do and how they do it, so for simplicity's sake, this is what you should do:

Make 4 duplicates (**Image > Duplicate**) of your image (which should leave you with 5 identical images) and then resize each image using a different resample method. Then sit back and see which one did the job you wanted done best.

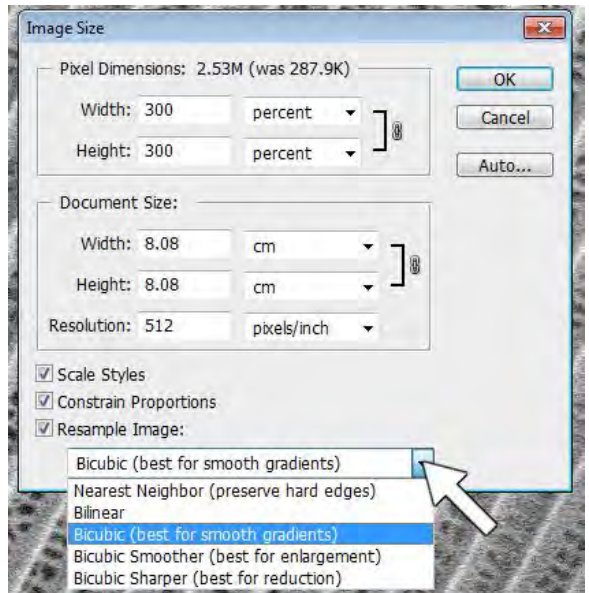

**5.** Depending on your image, different resampling methods might make your image look blurry, while reducing the noise in the image, while others might make your image look sharper, while increasing the noise in the image. Some methods will help, while others will hurt. It depends on your image, how much you are enlarging or reducing, etc.

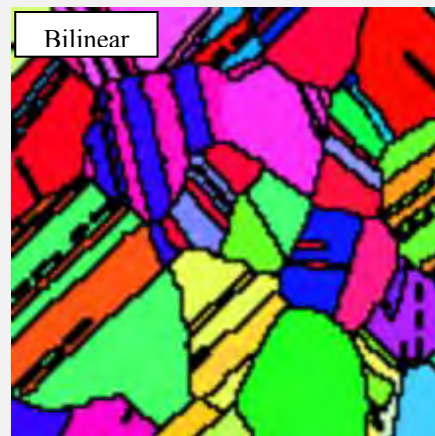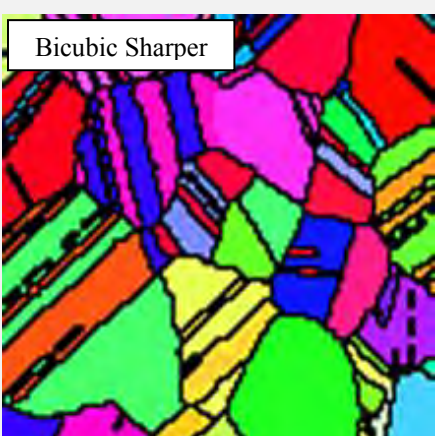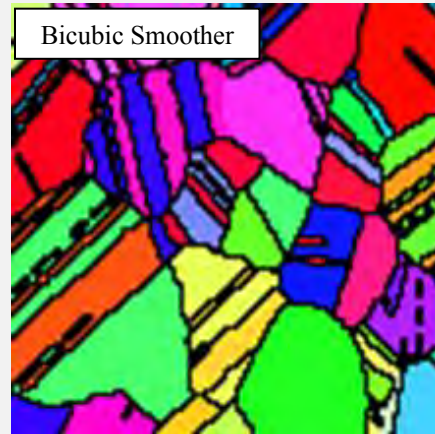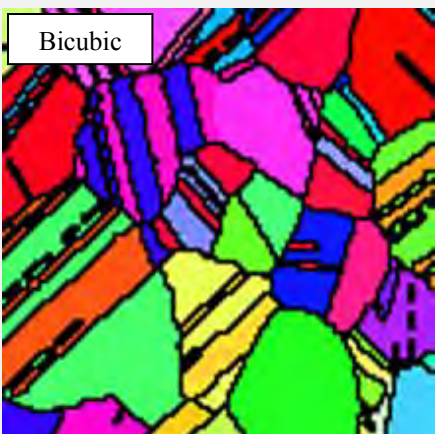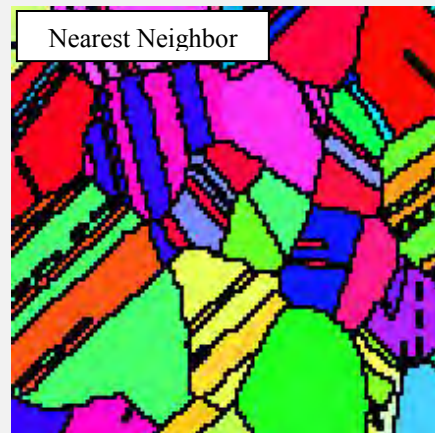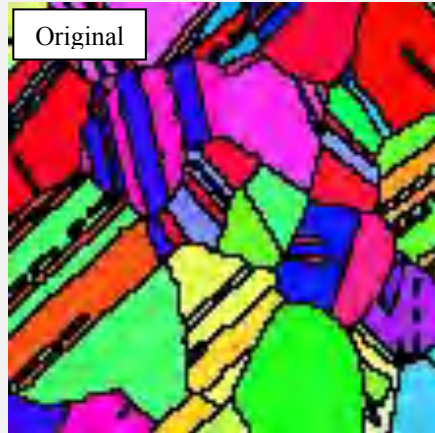

The differences may appear subtle here, but if you plan to blow up an image to print on a huge poster, small differences become very noticeable.

### III. Simple Pseudo Color – Color Table

Let's say we wanted to add a little color to our black and white image. This is a super quick and easy way to get a color image.


**1.** We'll start by opening our image and making a copy of it by clicking:

#### **Image > Duplicate**

Close your original image.

We don't *really* need to this, but it will prevent anything accidentally happening to our original image.


Once that's done we'll want to perform the Levels adjustment that we did in **Section I**.

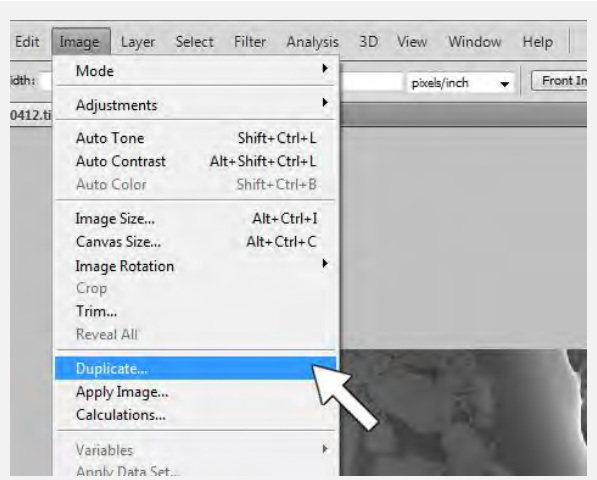

**2.** With this technique we are going to give a grayscale image fake color. In the example there is a grayscale SEM micrograph of gold dust that we will make gold colored.

We want to make sure our image is in **Indexed Color** mode:


#### **Image > Mode > Indexed Color**

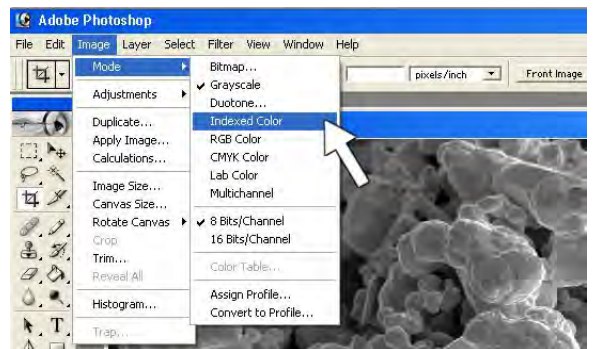

**3.** Next click on:


#### **Image > Mode > Color Table**

This will bring up the Color Table dialog box.

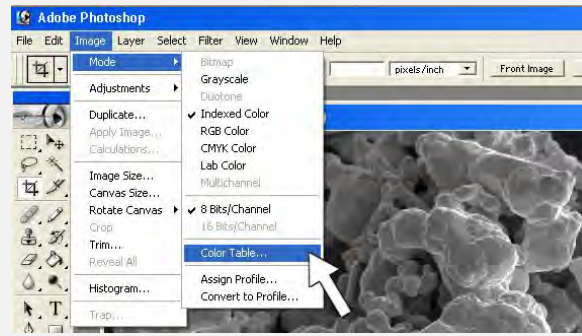


**4.** The **Color Table** box shows what color is being applied to each of the 256 possible colors in the image. If you want, you can click on each box and specify each color. *That would be insane.* Instead, we'll do something easier.

Click the top left box (which should be absolute black) and hold the mouse button down and drag it down to the box in the bottom right corner (absolute white) and release the mouse button.

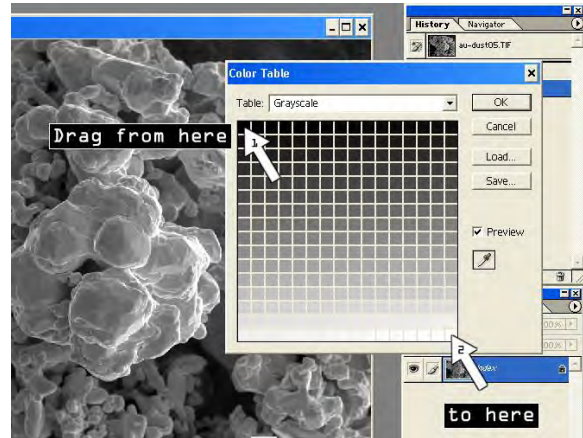


**5.** When you release the mouse button the **Color Picker** box will appear. Since black was the color at the top of the Color Table in the last step, we'll need to start off this step by picking the color you would like to represent the darker areas of your image. You can do that just by clicking on it. In this case, we will go with black. Click **OK**.

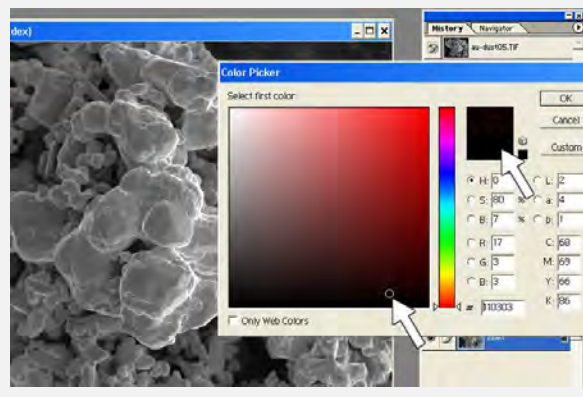


**NOTE:** If you picked the light color first when you should have picked the dark color first [or vice versa] it will be obvious pretty quickly as your image will look like a negative (with areas that should be dark being light, etc.). If this happens, just hit cancel on the Color Table and start over from Step 4.


**6.** Now choose which color you would like to represent the lighter areas of your image. In this case, we want gold to be colored gold, so a darkish yellow is selected. Click **OK**.

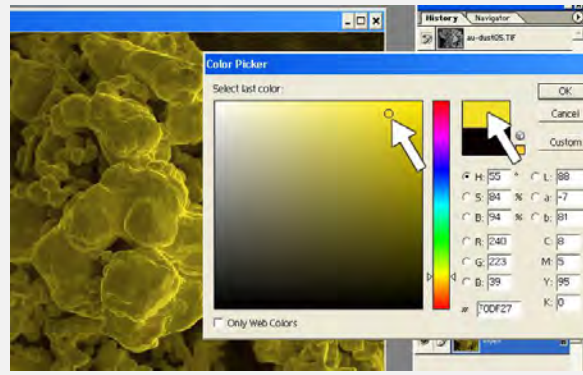


**7.** The Color Table will come up again and show you how your 2 colors have been stretched to fill all 256 colors. Make sure the Preview box is checked and if you like how it turned out, click **OK**.

You should be done! **Save** it!

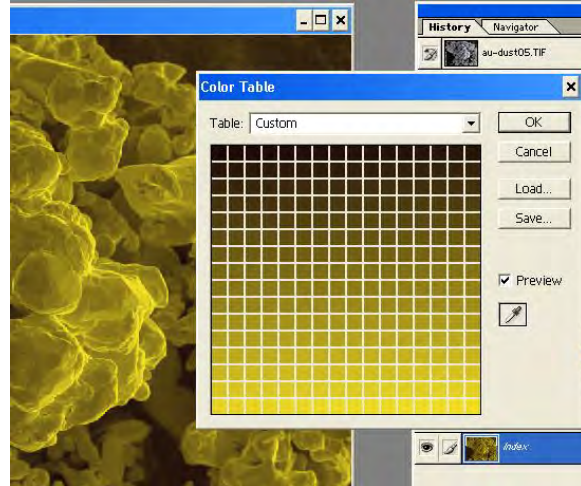

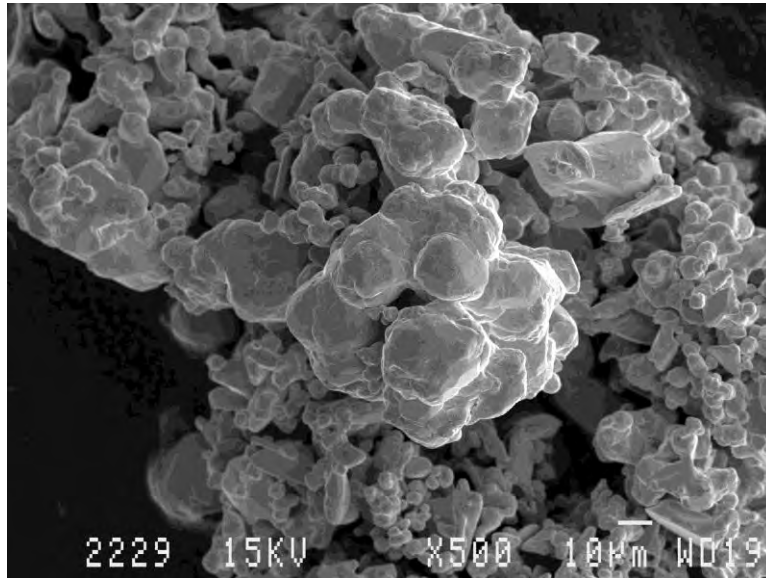

Here is the perfectly acceptable grayscale image.

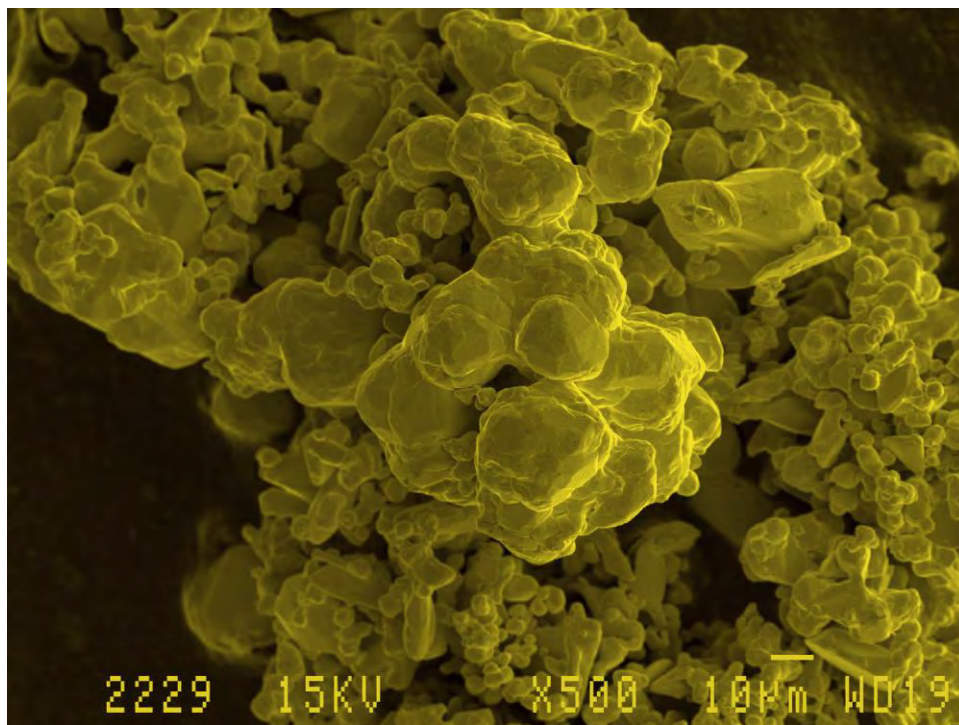

Here is the image once the pseudo color has been applied. There are many different ways to use the Color Table, but this is the simplest and quickest.

## IV. Simple Pseudo-Color – Gradient Map


The color table in Section II is great for very quick, single-color processing. This method we're going to discuss is a great "quick and not terribly dirty" method which can involve multiple colors.

**1.** We'll start by opening our image and making a copy of it by clicking:

### Image > Duplicate


Close your original image.

We don't really need to this, but it will prevent anything accidentally happening to our original image.

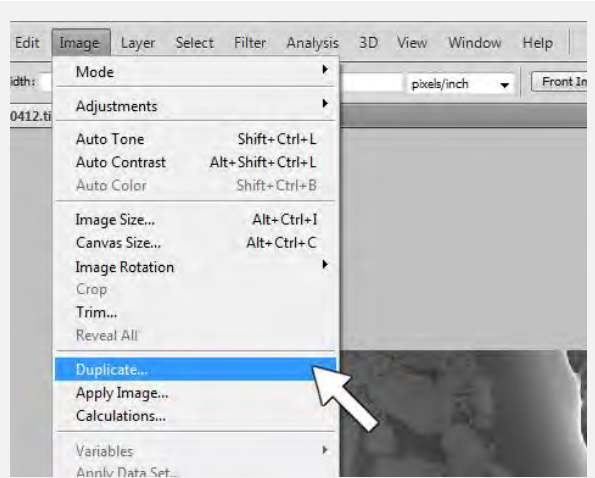

**2.** We're now going to switch the color mode of our image from **Grayscale** to **RGB Color** by clicking


### Image > Mode > RGB Color

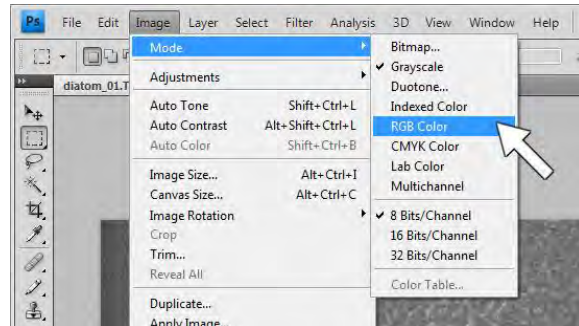

**3.** Perform the **Levels** operation on the image.

### Image > Adjustments > Levels


or **Ctrl + L**

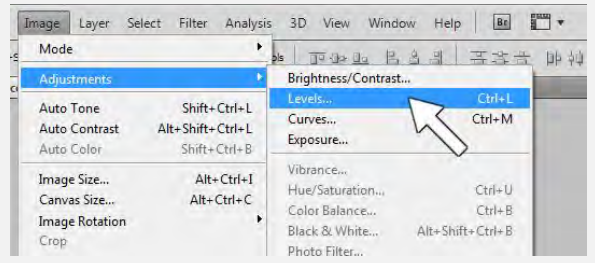




**4.** Now we're going to open the **Gradient Map** by clicking

**Image > Adjustments > Gradient Map**

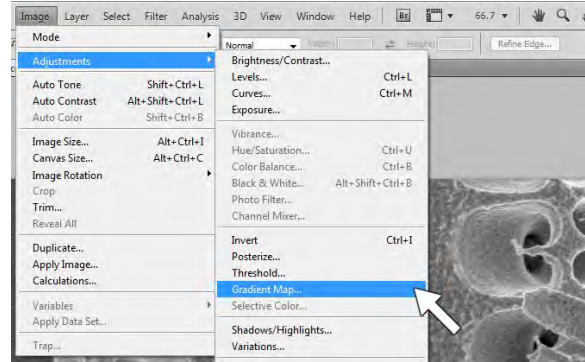


**5.** The **Gradient Map** box should pop up and our image is suddenly in whatever awful color the default settings are set to. In my case, it's this dull red. We're going to fix that by clicking that color gradient box thing in the middle there.

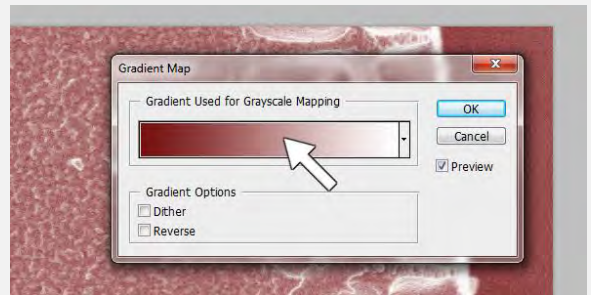


**6.** Up pops our Gradient Editor now. At the top we have a number of preset colors we can try out. I can guarantee you that using any of the presets will make your image look *awful*.


Fortunately that's not the point. The preset will give you a good starting point to work from. I want to start off with something simple, like 2 colors. So I'll click on a 2 color preset, like that green/red one, and we'll see what happens.

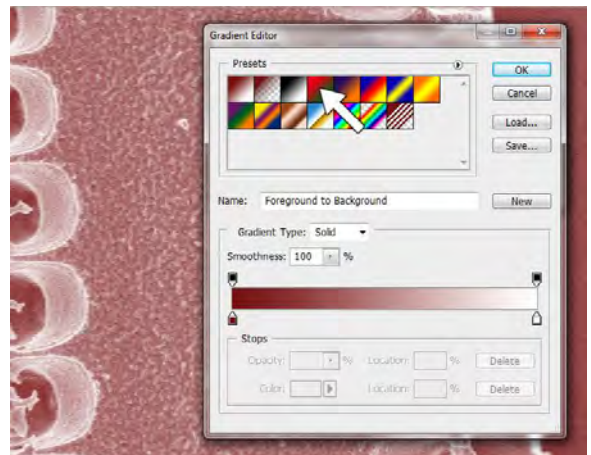


**7.** As I predicted, my image looks like some sort of Christmas decoration nightmare. But like I also said, that's ok.

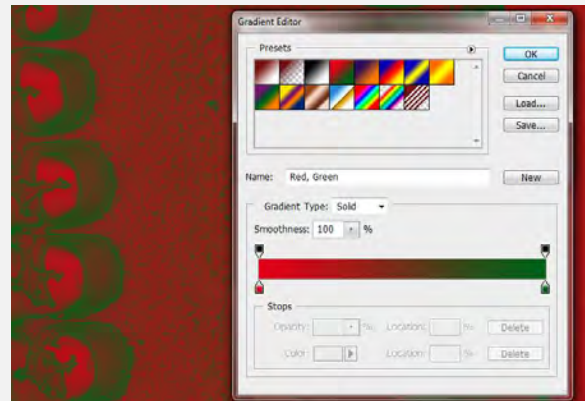


**8.** If we take a look at the bottom edge of the now Red and Green gradient bar, we'll notice a little red box and a little green box. They're called **Color Stops**. If we click on either of these boxes, we can select our color. We'll start off by clicking the red Color Stop on the left. We'll then click the box below it that says **Color** (and our current color is of course red, so the box will be red).

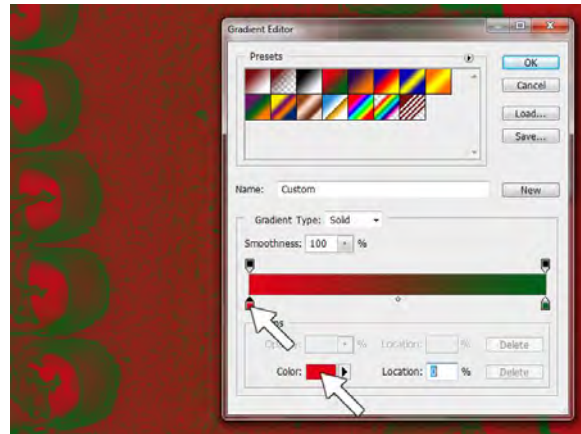


**9.** Once we click the Color box, the Color Picker should open up. Instead of red, I'm going to hunt around until I find a nice shade of dark blue that I like. Then click **OK**.

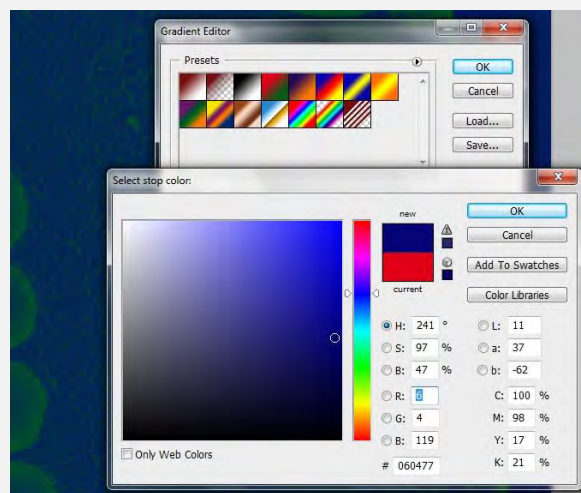

**10.** Now we'll go back and click the little green **Color Stop** on the lower right corner of the gradient bar and then click the (now green) **Color** box in the lower left of the window.

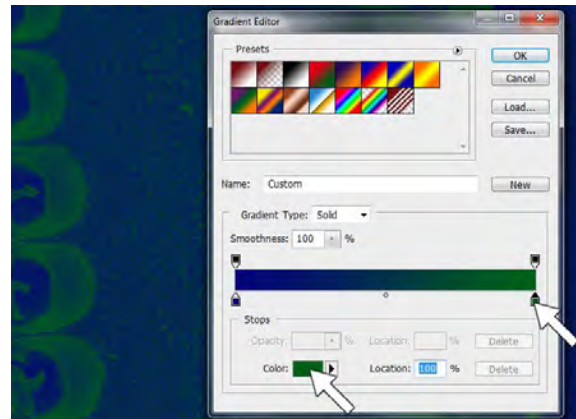


**11.** The color picker comes back and we'll poke around until we have a nice, lighter color. I just went with a lighter shade of green. Click **OK**.

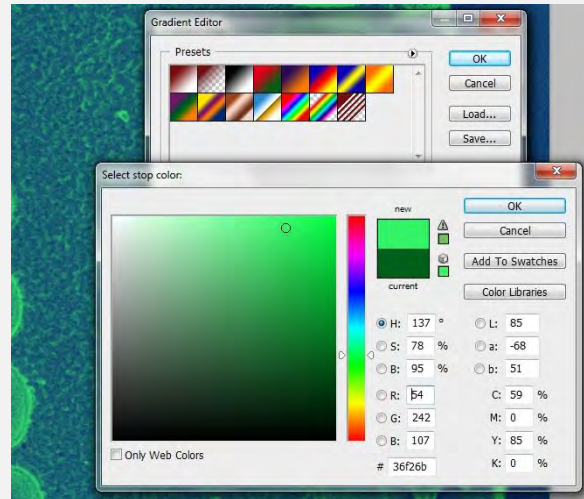


**12.** Now of course if we don't like what we've got, we can always go back and keep changing our colors until we have something we DO like.

*Another* thing we could do is move those little sliders around. If we bring the blue or the green box toward the center of the gradient bar, it will intensify that color. The little white diamond in the middle sets the midpoint between the 2 colors. (If the grays were well adjusted when you started off this method, you probably won't need to do this.) Anyway, if we like what we have, we'll click **OK**.

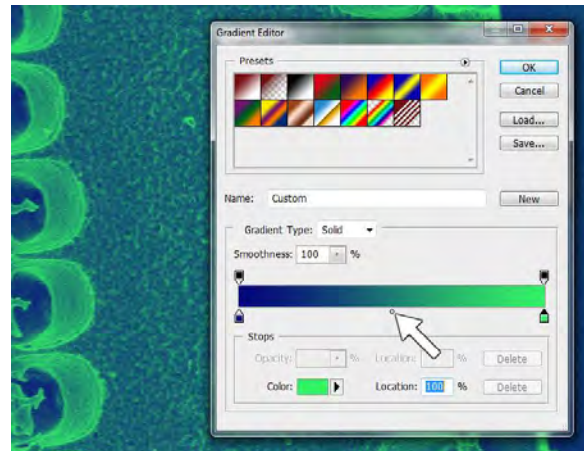


**13.** And click **OK** again.

You may want to do **Levels** again at this point. Maybe not.

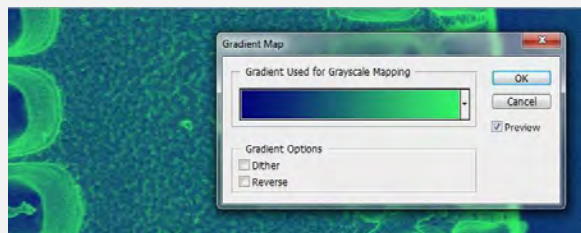

**NOTE:** There are a few other ways we can play with our image. One of them is by selecting one of the more complex presets.


As we can see here, we now have 4 colors to work with instead of 2. As per our previous instructions, we can change each of these colors and drag the little **Color Stop** boxes around to intensify or de-intensify each of the colors. You can easily spend all day picking colors and playing around with what part of your image they show up in.

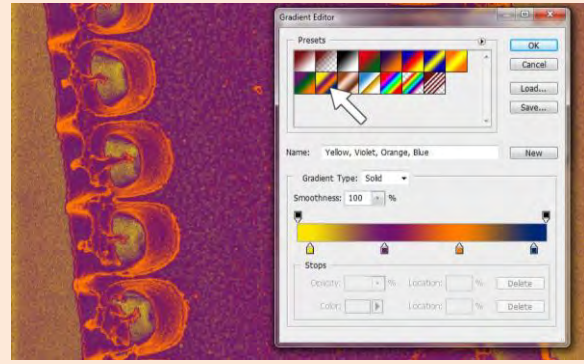


**14.** Whenever you have something you like, go ahead and **Save** it!


**File > Save**

Or **Ctrl + S**

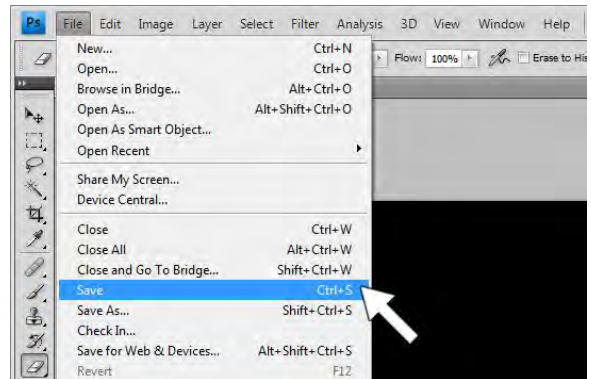

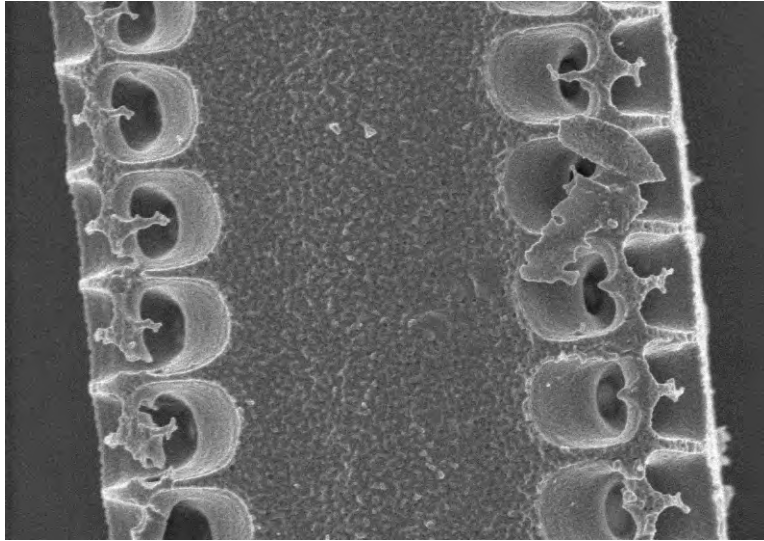

Our first image was this perfectly nice B/W offering.

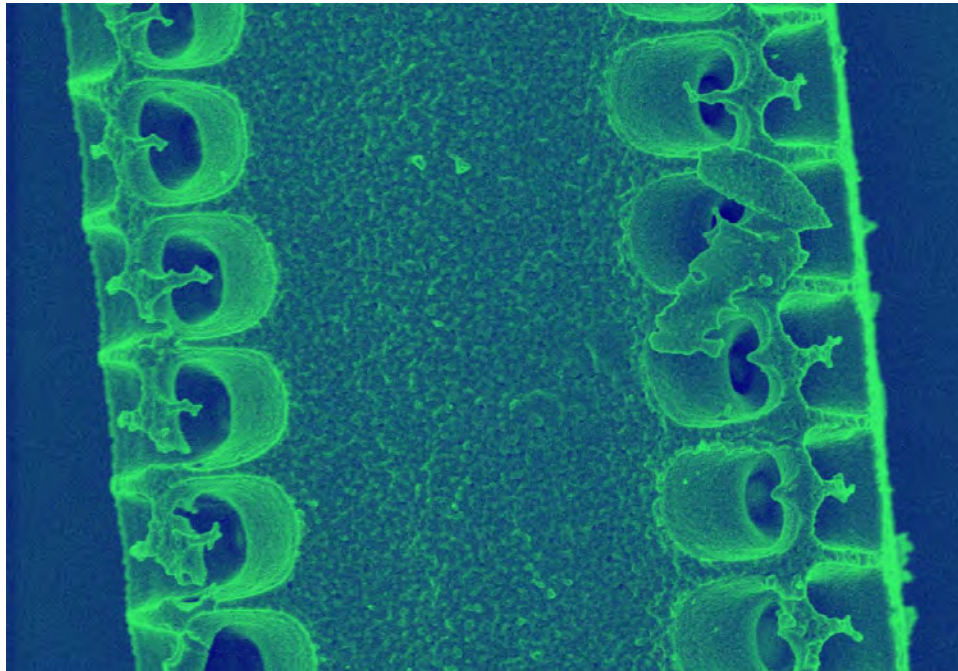

What we got in the end is a reasonable 2 color image. Is it great? Probably not! But the important part is that it was *FAST!*

## V. Complex Pseudo-Color – Color Table


Let's say we wanted to colorize our image with multiple colors, but with more discernment than the Gradient Map will give us. There are several different ways we can do this in Photoshop, but this technique should give us some of the best results.

**1.** We'll start by opening our image and making a copy of it by clicking:

### Image > Duplicate


Close your original image.

We don't really need to this, but it will prevent anything accidentally happening to our original image.

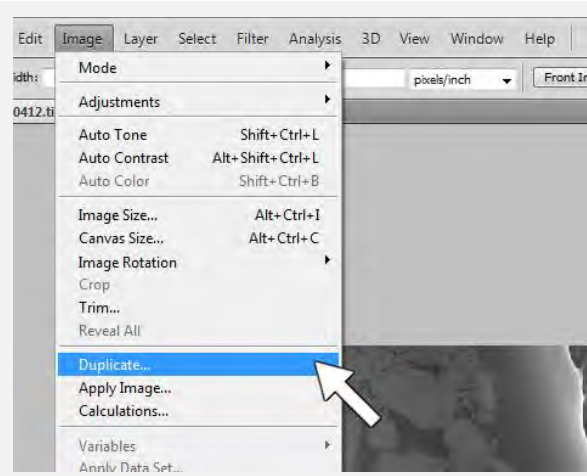


**2.** Usually we would do the **Levels** procedure here.

In this version of colorization, we **will not do Levels right now** and instead leave it until the very end. (Unless you really, *REALLY* want to do it now. There's realistically nothing I can do to stop you.)

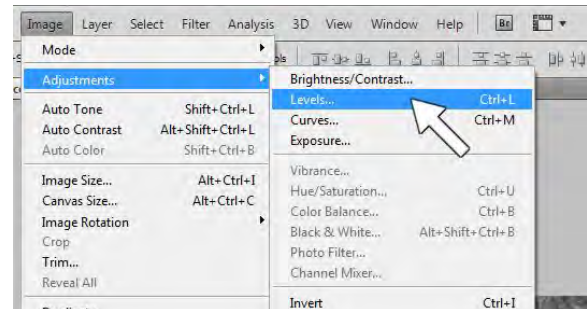


**3.** In this example, we're doing 3 colors. We'll need one copy of our image for each color. So we duplicate our image 2 more times so that we have a total of 3 images that all had the **Levels** procedure performed on them in the exact same way. (This is better than doing it separately on 3 different images.)

### Image > Duplicate

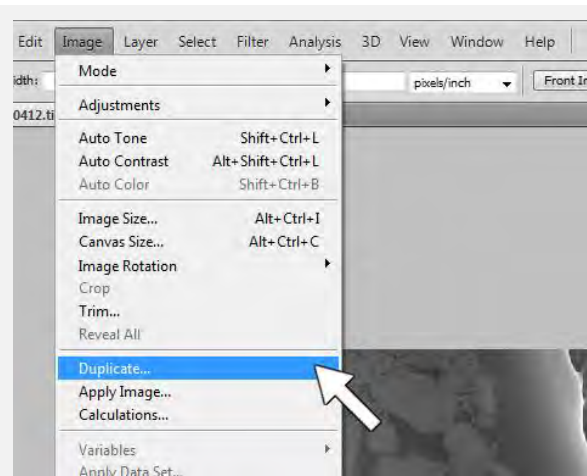

**4.** Select each image and set it to **Indexed Color Mode**.


**Image > Mode > Indexed Color**

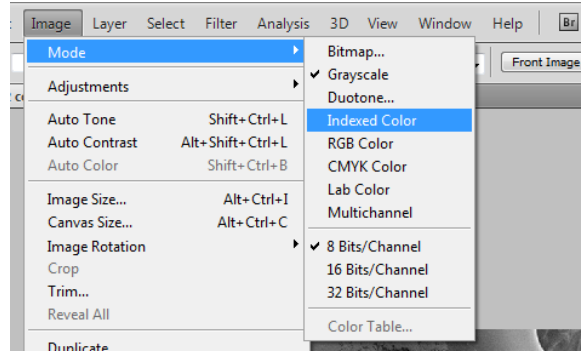

**5.** Now we will go through and colorize the 3 images in each of the 3 colors we want to end up with. (This may sound a little confusing or counter-intuitive, but trust me, it'll work.)


Select an image to start with and open up the **Color Table**

**Image > Mode > Color Table**

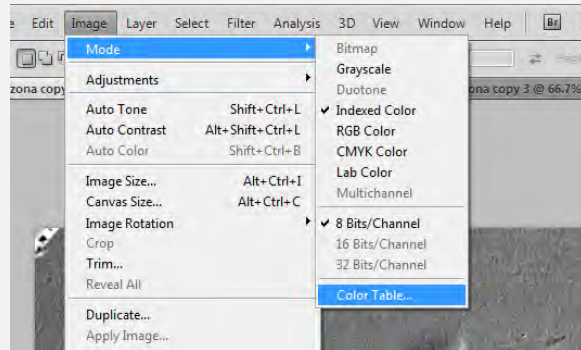

**6.** For starters, let's say we wanted to come up with the colors for the sun in the image. So we'll perform the Pseudo color procedure we discussed in **Steps 4 thru 7 of Section II**.

We'll repeat this process on the other 2 images. On one I'll choose the color for the sky and another I'll choose to color for the ground.

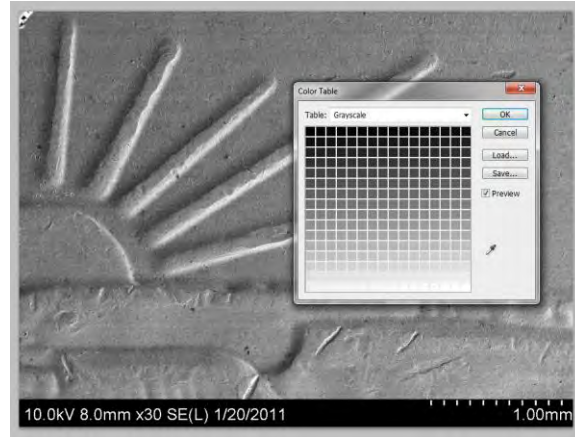

**NOTE:** If you wanted to get **REALLY** fancy, you could use the **Gradient Map** on this step instead. (See **Section III**.)

**7.** We should end up with 3 colorized images.




Now we need to decide which areas of which images we need to cut out.

The sky in the blue image seems to have the largest surface area, so I'm going to designate it as the background image. So all I'll need to do is cut the sun and the rays out of the yellow image and the ground out of the green image and then paste them on top of the blue image.

**8.** There are a couple of ways we can select these areas. It depends on the amount of time you have and the amount of effort you want to put into it.

The **best** way is by using the **Quick Mask**.


We enter **Quick Mask** Mode by clicking on the **Quick Mask** button, or **Q** on the keyboard.

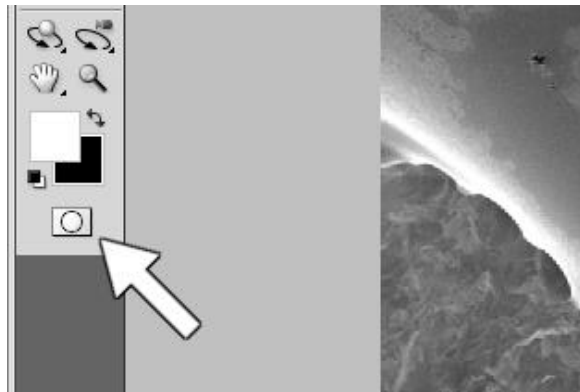

**9.** We'll do the selection process in 2 stages.


FIRST, we select the **Pencil** tool.

**Explanation of what we're doing:** By drawing over an area in the **Quick Mask** mode, the area will turn red. When we exit the **Quick Mask**, the area NOT painted over will be selected.

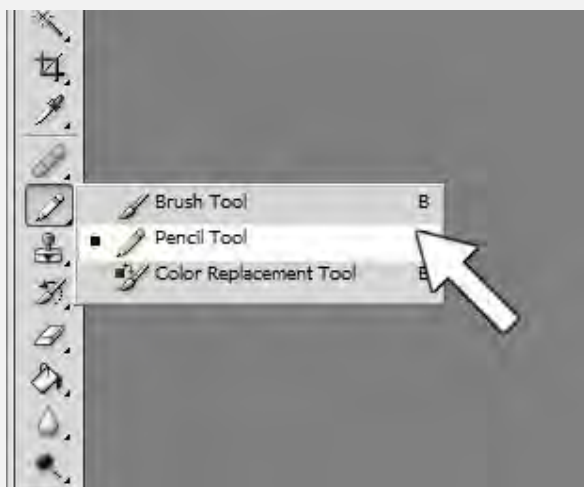



**10.** We can choose the diameter of our pencil by clicking on the little down arrow next to the **Brush** size indicator in the upper left corner of the screen.

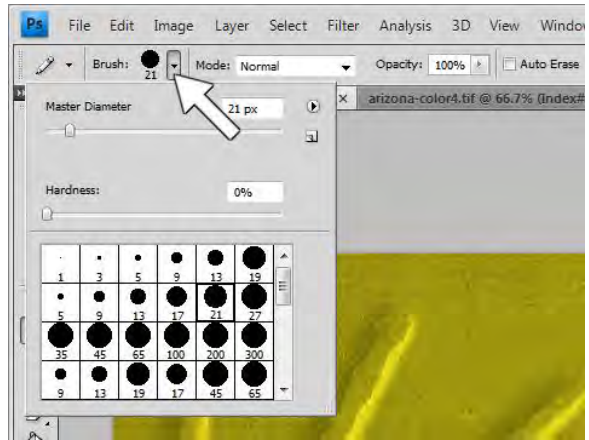

We need to select a size that will allow us to cover a lot of territory and yet still be small enough to fit inside where we want to paint.

**NOTE:** You'll probably need to change the diameter a few times as you're working.


**11.** Since we're on the yellow image, we want to take out the sun part of this image. By using the pencil tool, we're going to paint out **ONLY THE CENTER** of the sun and the rays, but **AVOID THE EDGES**.



As we can see here, areas that you select show up as red.

This can be a long, tedious process.

**NOTE:** Keep in mind that if the features you're trying to select are too small to see clearly, you can always zoom in by selecting **View > Zoom In**

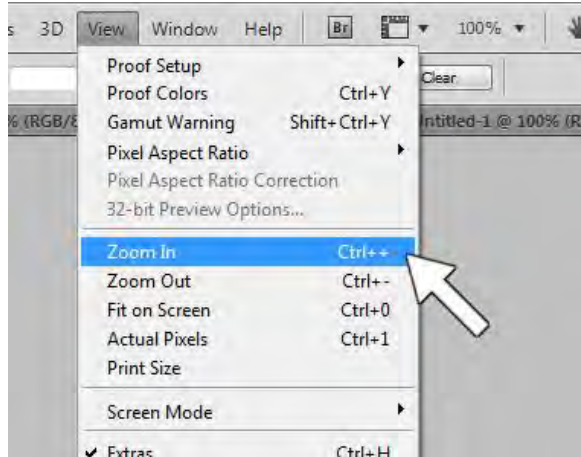

or **Ctrl + +**

You can do this multiple times.


Then, to zoom out, select

**View > Zoom Out**

or **Ctrl + -**


**12.** Once we have the center parts bulked out with the **Pencil** tool, we're going to select the **Brush** tool and use it to go over all of the edges.

And again, you may need to change the Brush size like we did with the **Pencil** in **Step 10**.

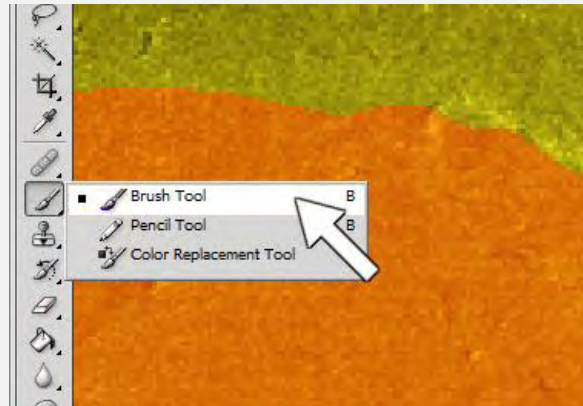

**NOTE:** There's a good reason for using the **Brush** on the edges instead of the **Pencil**.


The **Pencil** tool will make hard edges around the area we select. For some things, that's ok. But when we cut this area from this image and drop it on top of another image, the abrupt change of colors might not look very good. We can see the **Pencil** edges where things were painted over are pretty obvious and are a bit abrupt.

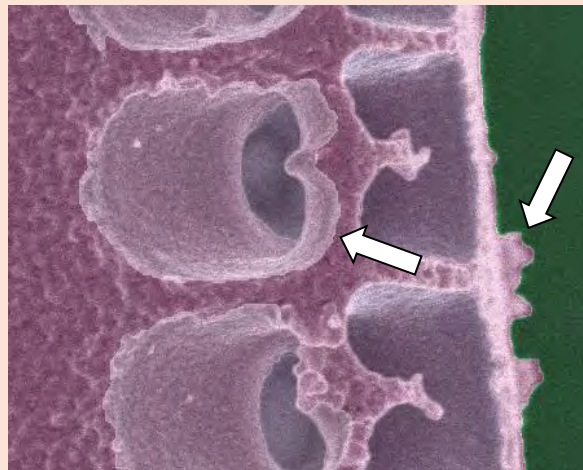


The **Brush** tool will give us a softer edge around our selected area which should blend a little bit with any adjoining color areas.

**NOTE'S 1<sup>st</sup> NOTE:** So I know your question now is why don't we just use the **Brush** to select **EVERYTHING**? Well, because of the non-solid nature of the **Brush** tool it's easy to miss things with it.


As we can see with this image, I used the brush to select this area. I went over every part of it, but you can see some parts are a little thin. Those thin areas will show up as transparent when we copy and paste this area on top of our background. We don't want that.

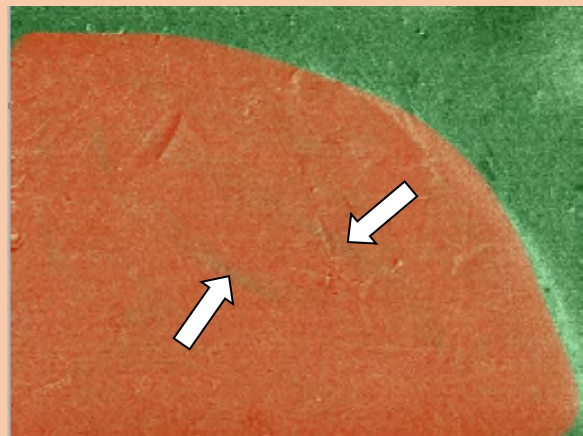

**NOTE'S 2<sup>nd</sup> NOTE:** There are certainly instances where we can use the **Pencil** tool on everything. It completely depends on how well defined your image is. This image for example, has only 2 color areas and the edges on the structure are pretty nice and sharp, so the **Pencil** tool was used for everything and it looked pretty good.

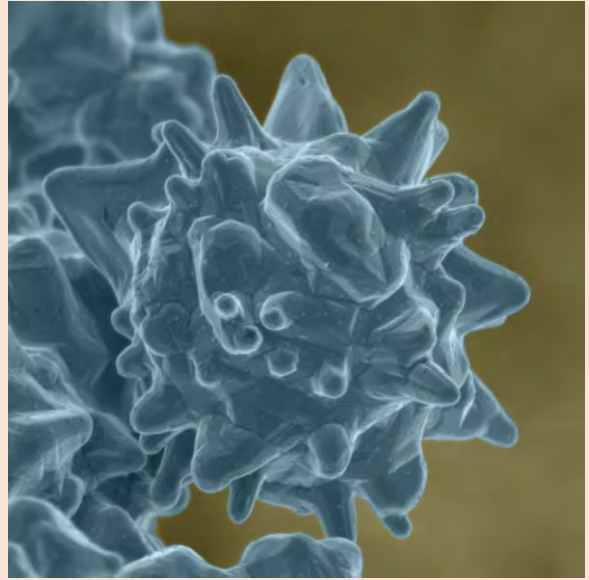

**13.** Once you're done with all the painting, you should have something that looks like this. Areas you want selected are red.

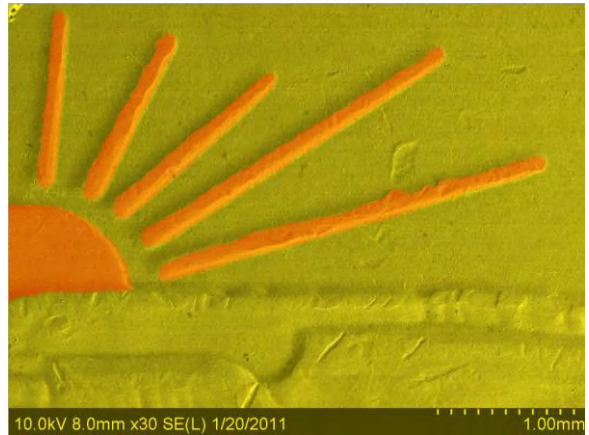


**NOTE:** If you had a major (or even minor) screw up while you were painting out your area, you'll want to fix it. We can do that easily enough by exiting the **Quick Mask (Q)** and inverting the selected area.

**Select > Inverse**

or **Shift + Ctrl + I**

Turn the **Quick Mask** back on and you can fix any areas you might have strayed outside the edges of your structure.

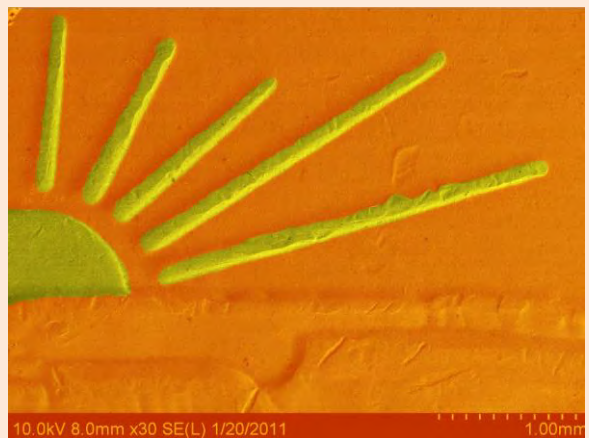

**14.** When you're done painting, exit the Quick Mask by clicking on the **Quick Mask** button or **Q** on the keyboard again. Now everything that was NOT red is now selected.

This is backwards of what we wanted, so we'll need to invert the selection.

**Select > Inverse**


or **Shift + Ctrl + I**



**15.** Now we'll copy the area we selected.

**Edit > Copy**


or **Ctrl + C**

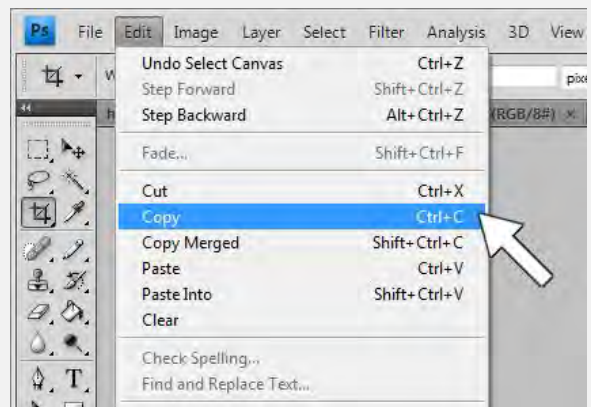

**16.** We said earlier we were going to use the blue image as our background, which means we want to paste our just copied parts from the yellow image into it.

Before we do any pasting, we need to select our blue background image and switch it from **Indexed Color** to **RGB Color**.

**Image > Mode > RGB Color**

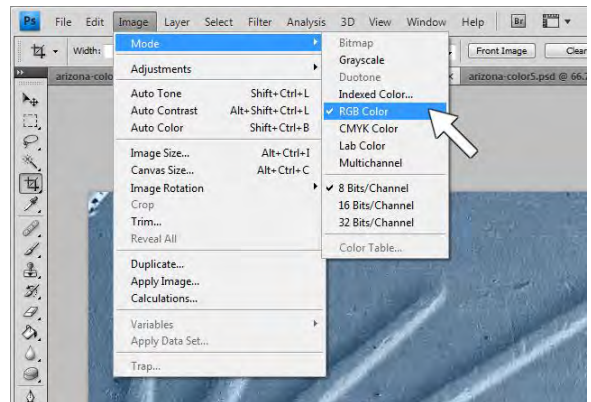


**17.** Now with the blue image selected we can go ahead and paste our yellow selection onto it.

**Edit > Paste**


or **Ctrl + V.**

It probably won't be lined up correctly, but we can leave it alone for now.

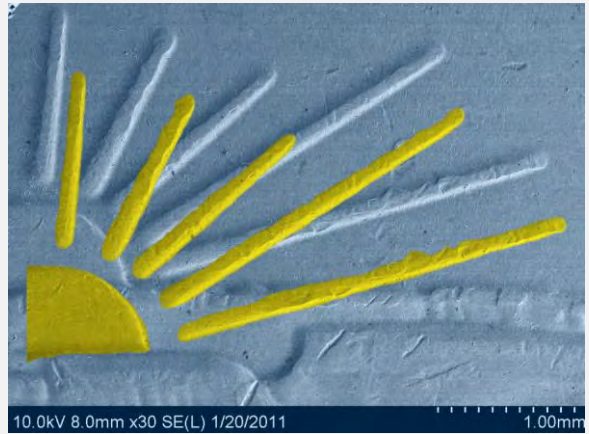


**18.** For now we need to go back to our green image and paint out the parts we want from it.

As before (**Steps 9** thru **17**) we enter the **Quick Mask**, use the **Pencil** for most of it and the **Brush** tool for the edges, exit the **Quick Mask** mode, **Invert** the selection, **Copy**, then **Paste** it into the blue image.

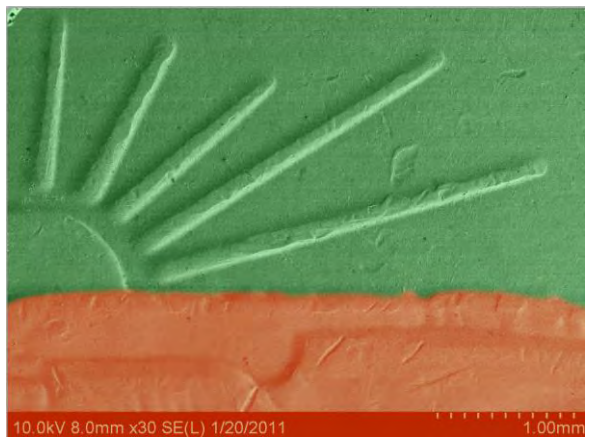


**19.** So now we have our image with all 3 colored areas in it. BUT, they're a bit misaligned.

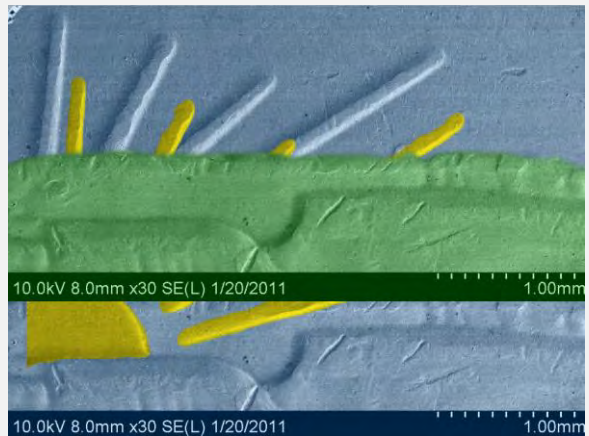

**20.** We'll start aligning them by going to the Layers tab in the lower right corner of the screen. We start off by selecting the layer we want to move around.

The most recently created layer should be selected (the one from the green image) and that's as good a place to start as any.

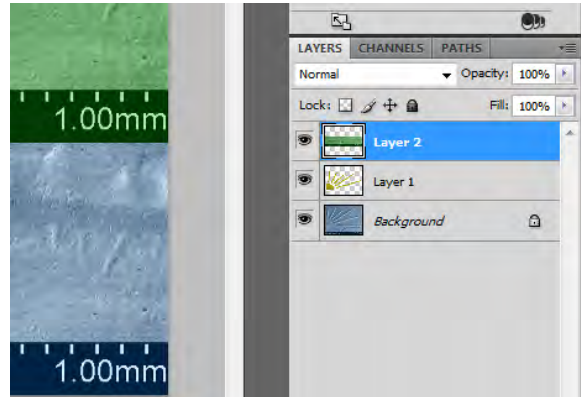


**21.** Select the **Move Tool**, or **V** on the keyboard.

Using the **Move Tool**, you can click on your layer and drag it around the image and position it where you want.




**22.** This layer was easy to position, since it shares 3 edges with the background image. (The edges of the layer should snap to the edges of the image.) The next layer will be a little trickier.

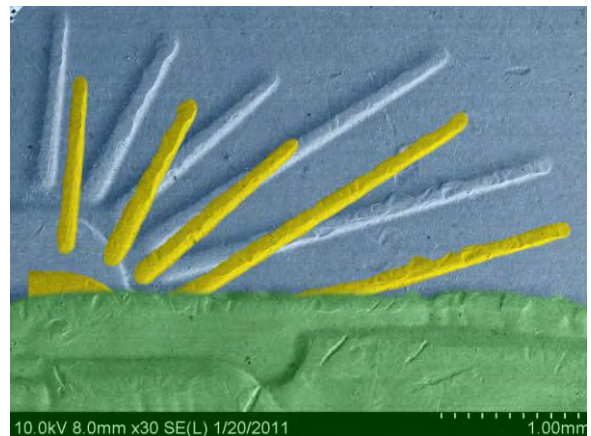


**23.** We'll select the yellow layer from the **Layers** tab as in **Step 20**. Because this layer doesn't fit perfectly against at least 2 corners of the image like the last layer did, we'll need to do some fine positioning. To do that, it's easier if we set the layer to around 50% opacity, so we can see it and the layer behind it.


You can do this by clicking on the **Opacity** box there in the **Layers** tab and either typing in 50 or dragging the slider bar around. (It doesn't need to be exactly 50%).

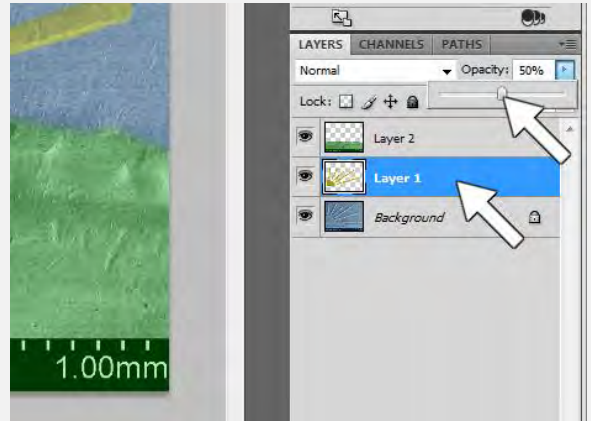


**24.** We can then use the **Move** tool to get the layer mostly lined up, then we can use the arrow keys on the keyboard to nudge the layer around by single pixels until it's all lined up. Then we'll bring the opacity of the layer back to 100%.

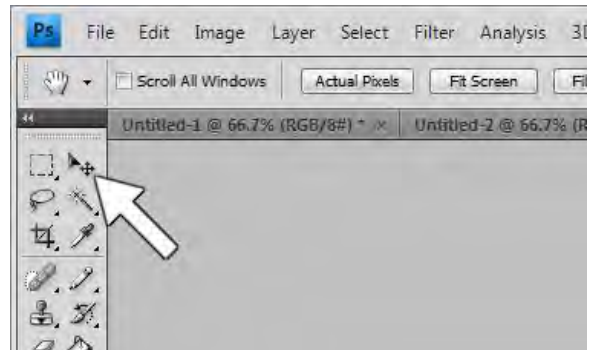



**25.** Before we save our image, we have one last thing we'll want to check. Some of our color layers might be overlapping a bit. We might want to check to make sure that they overlap the way we want them to. Let's look closer at this area on our image where the sun and the ground meet.

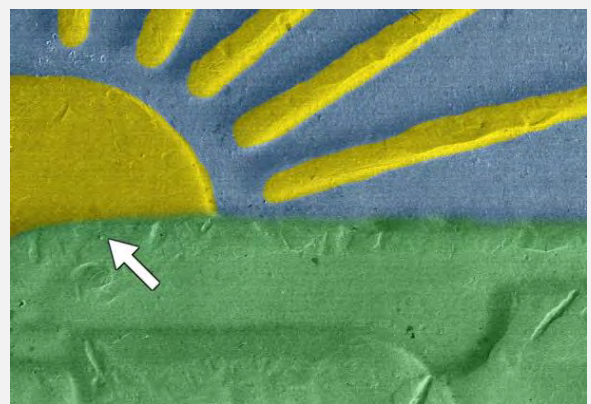


**26.** We can arrange which layer is on top of or below others by dragging them around in the **Layers** tab.

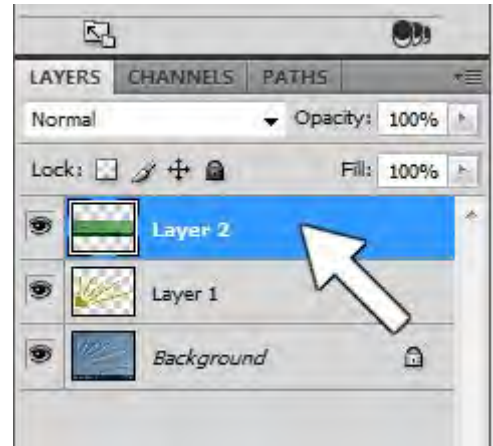

**26.** The image on the left has the **yellow** layer on **top** with the green layer behind. The image on the right has the **green** layer on **top** with the yellow layer behind. The right image looks the smoothest, so I'll stick with that arrangement.

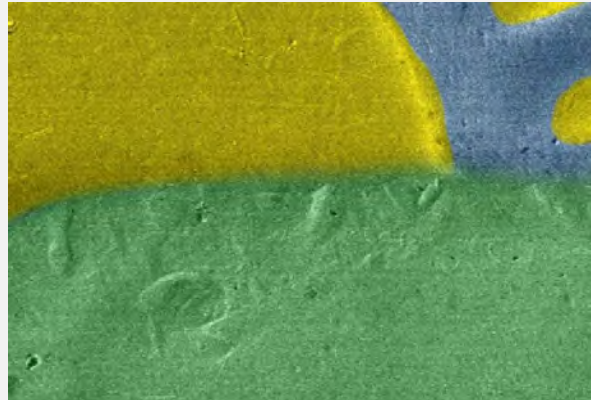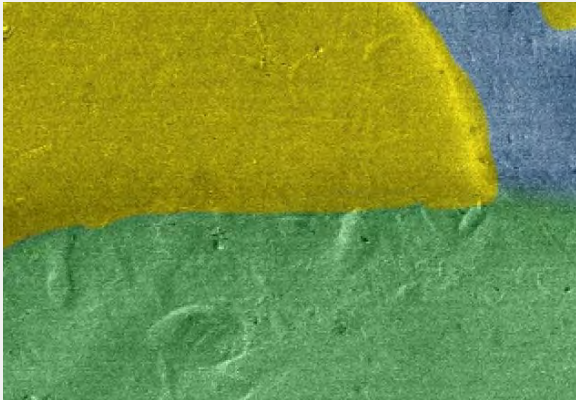

**27.** OK, we've made it to the end here and **DISASTER!** We don't like how the color turned out in one of our areas. (Like, EEW my sky is too purple!) No sweat. We have 2 ways to fix this.

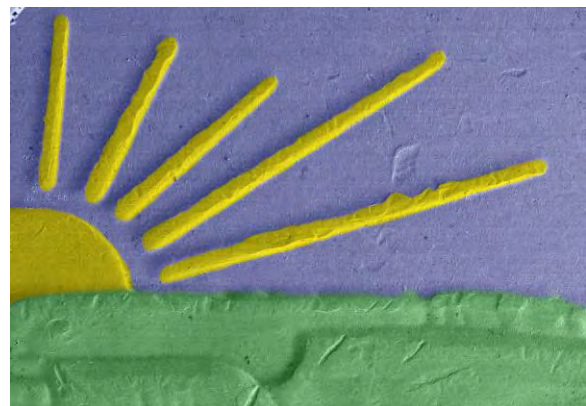




**27. Option 1** Go back to the image whose color you want to change. The area you wanted should still be selected. Leave that alone, don't deselect anything! All you need to do is go back to the color table, change the colors, copy the already selected area and paste it into the final file. You can then determine if that was better or not. If it was, you can throw away the less acceptable color layer.




**27. Option 2A** The other way we can tweak our layers is by adjusting the **Color Balance**. Start by selecting the layer you want to adjust, and click:

**Image > Adjustments > Color Balance...**

Or **Ctrl + B**

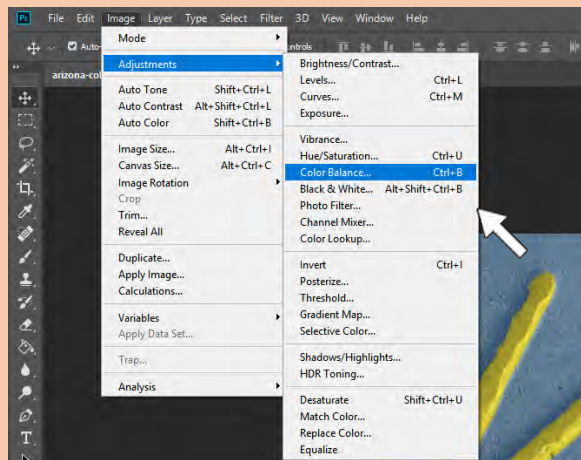


**27. Option 2B** The **Color Balance** box will come up. Make sure the **Preview** box is checked and then by moving around the 3 sliders we can adjust the colors in that layer quite a lot. (If you really know what you're doing you can select shadows, midtones or highlights and adjust those.)

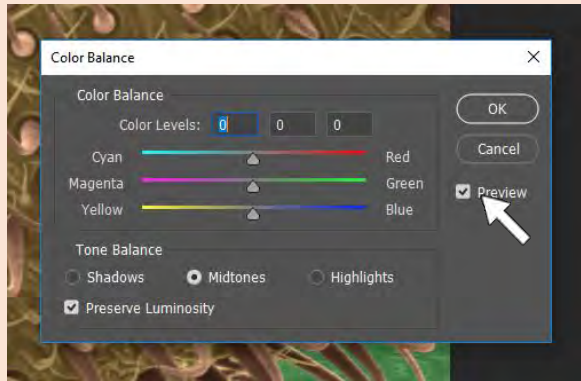

**28.** It is *FINALLY* at this point that we want to adjust the **Levels**. We will select each layer and perform **Levels** on all of them individually, only now, *AFTER* they've been colored. (The **Levels** operation is described in Section I. 3-6).


(This may not make sense now, but we'll see the example at the end to see how it turned out.)

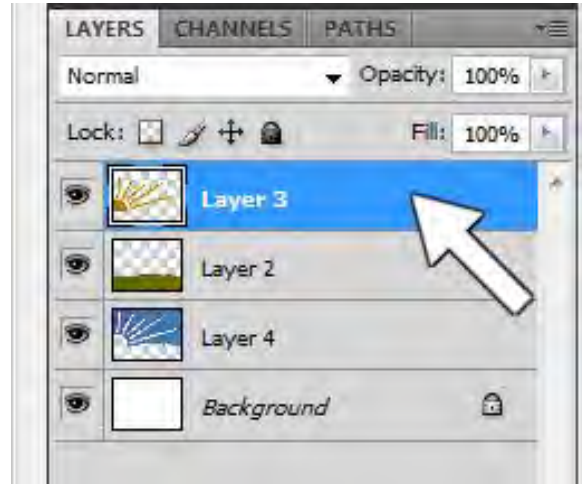


**28.** We should have a completed image now, we can go ahead and save it as whatever we want.



**File > Save**

Or **Ctrl + S**

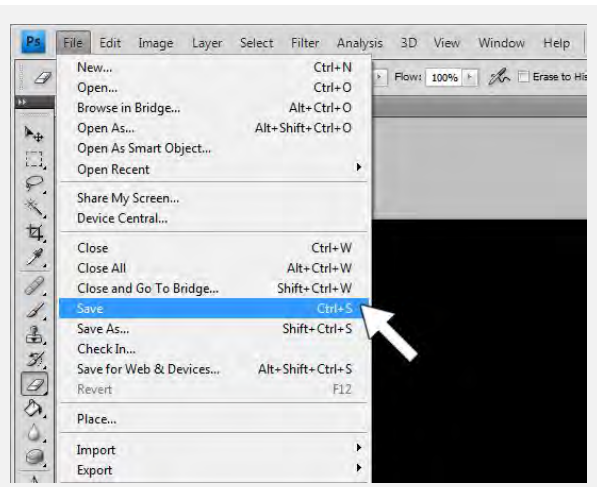

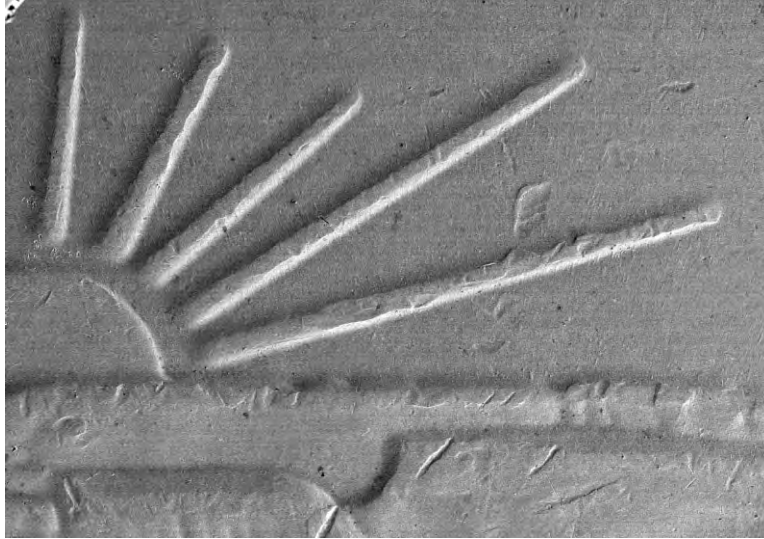

Here we have our original (slightly boring) sunrise over Arizona.

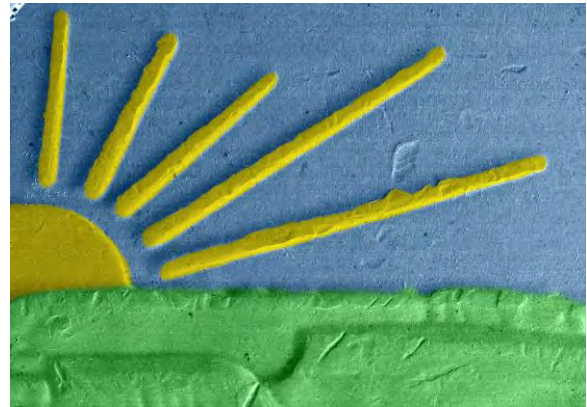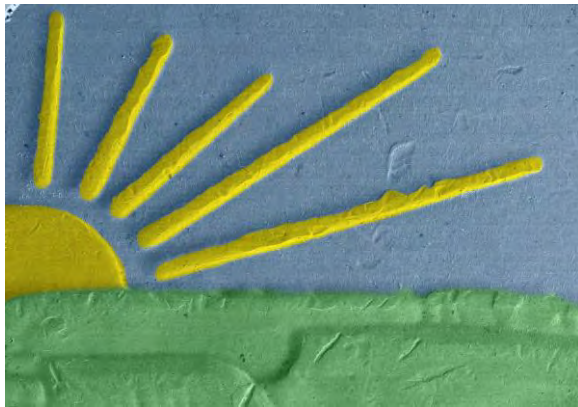

In the end, we have ourselves this nice colorful image that really kind of looks like something! (The colors don't always need to be this happy, that's just how I ended up making them for the purposes of demonstration.)

So on the left we have the image if we had done the **Levels** at the very beginning. It looks fine, but in comparison to the image on the right (which had the **Levels** done on each layer at the very end, AFTER coloring) it looks dull. Doing the **Levels** at the very end on each layer makes the colors pop a little bit more, but more importantly, it makes the shadows pop. Next to choosing the right colors, doing the **Levels** at the end on each layer will help to give your image that extra bit of reality.

## \*ADDENDUM to Section V – Complex Pseudo-Color

So I mentioned at the beginning of **Section V** that there are a number of ways that we can go about grabbing different areas of our images to throw together later into a final colorized image. Here we will discuss some of those other ways. They probably won't be better than using the **Quick Mask**, but they *may* end up being a little faster.


**TO START:** We'll follow the instructions in **Section IV** up through **Step 7**. Now instead of using the **Quick Mask**, we'll look at the use the **Quick Selection Tool** (**Section A**) and the **Lasso Tool** (**Section B**). We'll start off with the **Quick Selection Tool**.

### A. Quick Selection Tool

**A1.** In CS4 the **Quick Selection Tool** is grouped with the **Magic Wand**. If the **Quick Selection Tool** is not showing, you will need to right-click the **Magic Wand** to select it.


Once it's selected, make sure the **New Selection** option is selected on the toolbar.

Also on the toolbar, make sure the brush size is set small enough to fit into any small areas you want to select, but large enough that you're not painting with it all day.

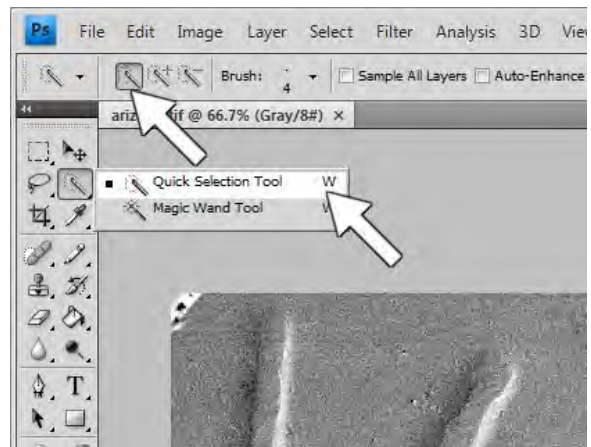

**A2.** Lastly before stating, we need to change all of our already colored images from **Indexed Color** mode to **RGB Color** mode.


**Image > Mode > RGB Color**

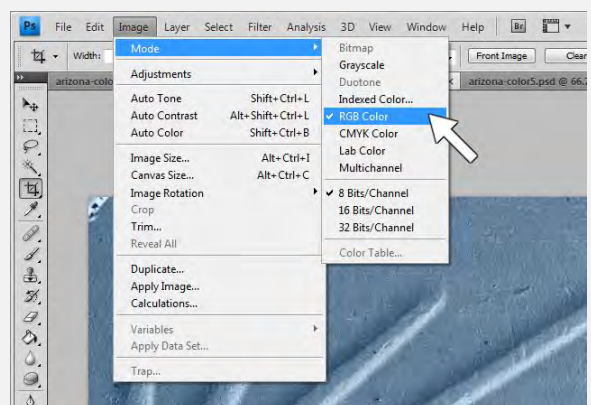

**A3.** Take the **Quick Selection Tool** and draw over the areas you want to select. As you draw, the selection will flood out across your image until it finds what it perceives as the edges of your object. This is not always perfect. When doing this on my image, I tried to just select the sunbeams and ended up with the whole top half of the image selected.



**A4.** To fix that, we just go back up to the toolbar and click on the **Subtract from selection** option. We might also select a smaller brush size by clicking the little down arrow just to the right of where it says **Brush**.

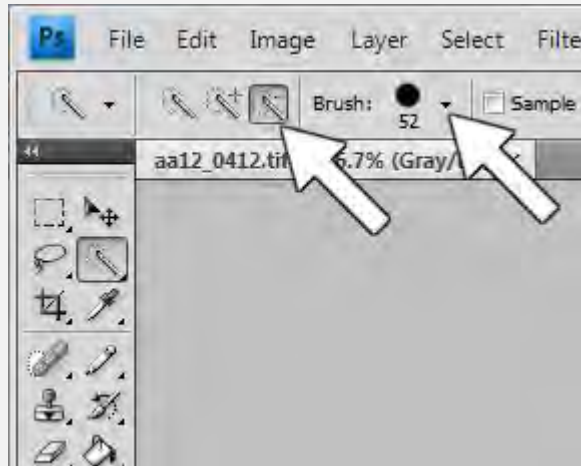


**A5.** By going back and forth between the **Add to** and **Subtract from** options we can make sure that only the areas we want selected are. The speed of this step depends greatly on how well defined the area you want to select is and how nice you want it to be.




From this point, you can continue on with the normal procedure lined out in **Section IV** from **Step 15** onward.

**NOTE:** Electron micrographs are not often known for their well-defined edges, and that will make this technique difficult. Another major issue you might have in using this tool is that in going back and forth between adding and subtracting different areas, sometimes you'll end up with stray bits left selected that you didn't want. In this image for example, we have a big green line going through the image. If I pasted this into the main image with the other color layers, I'm going to have this green line along the edge of one of my sunbeams. That's not cool. So I would need to go back through with the **Quick Selection Tool** and make sure I got rid of those adequately. But sometimes they're very difficult to see until it's too late.

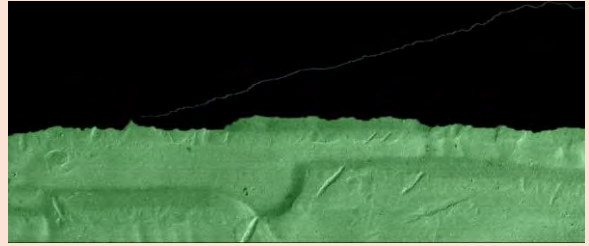


**SECOND NOTE:** Instead of the **Quick Selection Tool**, why not use the **Magic Wand** dummy? You **CAN!** But you'll probably end up with a border of some kind around the area you selected. (This sometimes also happens with the **Quick Selection Tool**.)

I took out the background of this image and you can see a dark line all around our structures. If it were colorized and recombined with a colorized background, that line is going to be **EXTREMELY** noticeable.

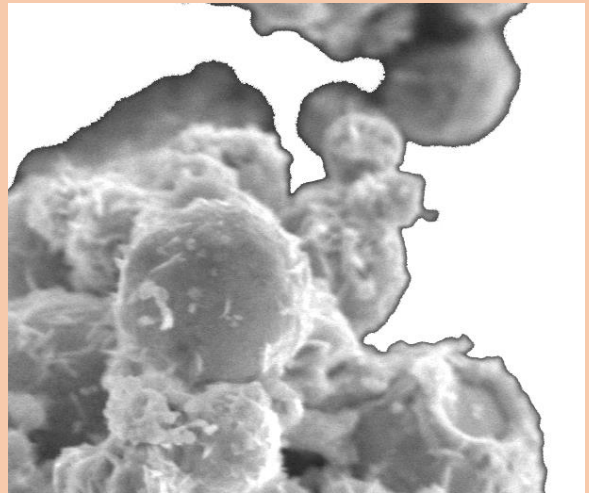

**SECOND NOTE'S NOTE:** As another aside, you **CAN** select large areas of your image with the **Quick Selection Tool** or the **Magic Wand** and then go over all the edges again with the **Pencil** or **Brush** tool. *Sometimes* that might save you time, sometimes just sticking with the **Pencil** is best.

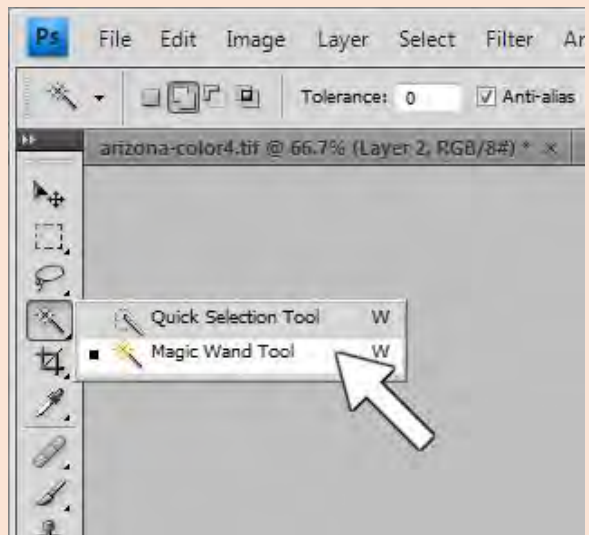


## B. Lasso Tool

**B1.** Now we'll look at the use of the **Lasso Tool**. This will probably be the fastest technique and can give pretty decent results.


Select the **Lasso Tool**.

Set the **Feather** to ~5 pixels to start. (This will make the edges of your selection area softer. You might need to play around with this setting depending on your results.)

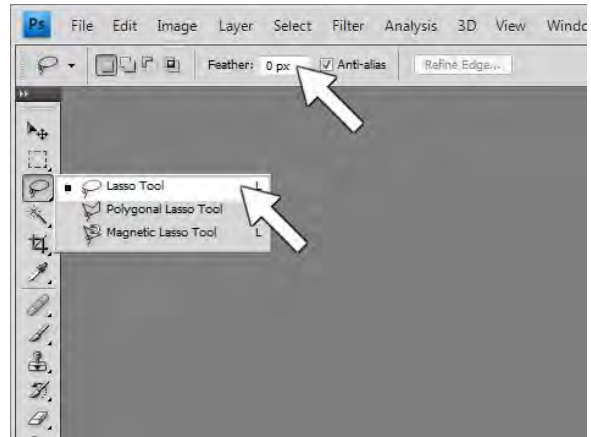

**NOTE:** Another option might be to use the **Magnetic Lasso Tool**.

With this tool, the selection area will automatically snap to what the software determines are edges.

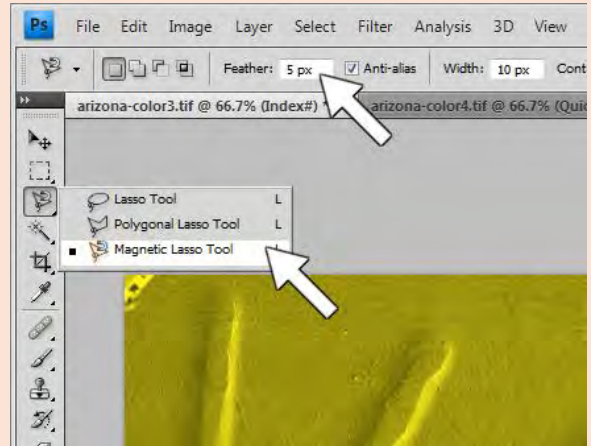


**NOTE 2:** As cool as that sounds, this might end up being more of a hindrance than a help, depending on how sharp the edges of your structure are.

We can see that this selection area will need a lot of cleanup work to make it acceptable.




**B2.** Take the **Lasso Tool**, hold down the left mouse button and **CAREFULLY** draw around the area you want to select. When you release the mouse button the selection will snap back to where you started. So be careful where you do that. We can see in this example, I released the button too soon and it cut off part of the sun that I wanted to select.

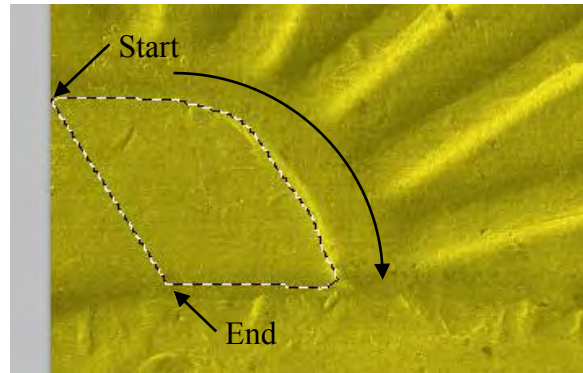


**B3.** To select multiple areas **OR** to fix missed areas like in **Step B2**, hold down the **Shift** key, then start drawing over the new area to select.

If you're doing this to fix missed areas, you'll want to overlap the existing selected area as you draw.

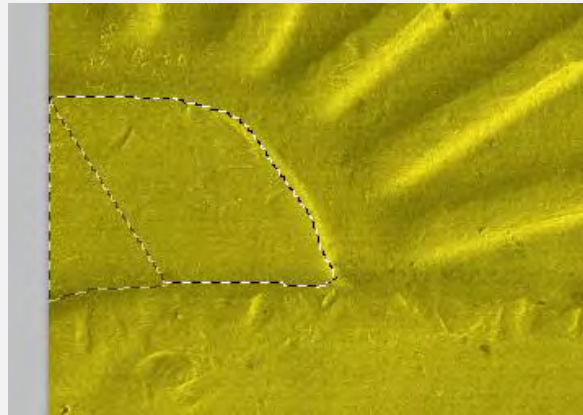

**B4.** Then when you release the mouse button, the overlapping areas combine to form a single area.

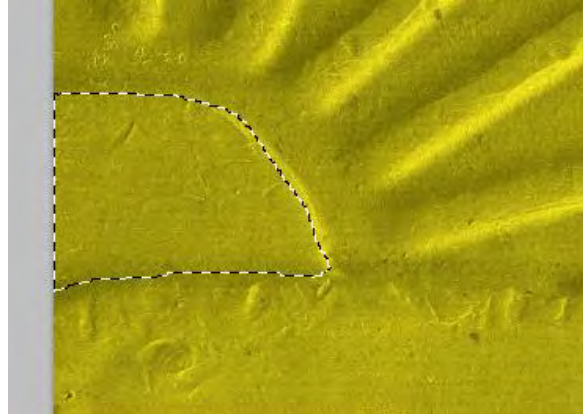



**B5.** OK, once we have everything on this image selected that we want, we might find that we selected extra areas that we didn't want.

To get rid of those areas we have a couple of options.



**B6.** If we want to stick using only the **Lasso Tool**, we'll **Invert** the selection


**Select > Inverse**

or **Ctrl + Shift + I**

Then do what we did in **Step B3** and **Lasso** out the offending areas.




**B7.** Another option might be to just turn on the **Quick Mask** and fix up the edges that way.

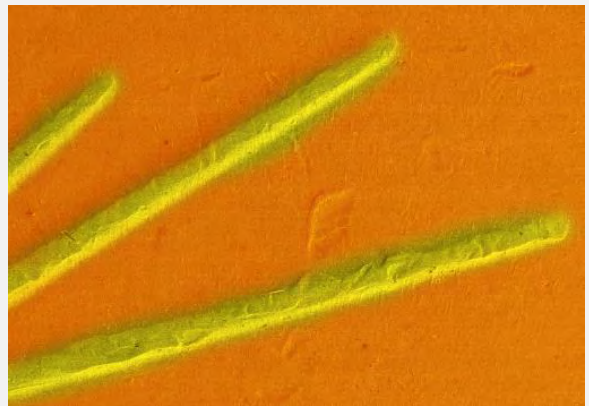


## Results


This is our first image that we combined using the **Quick Mask**. All the edges look nice and smooth and pretty.

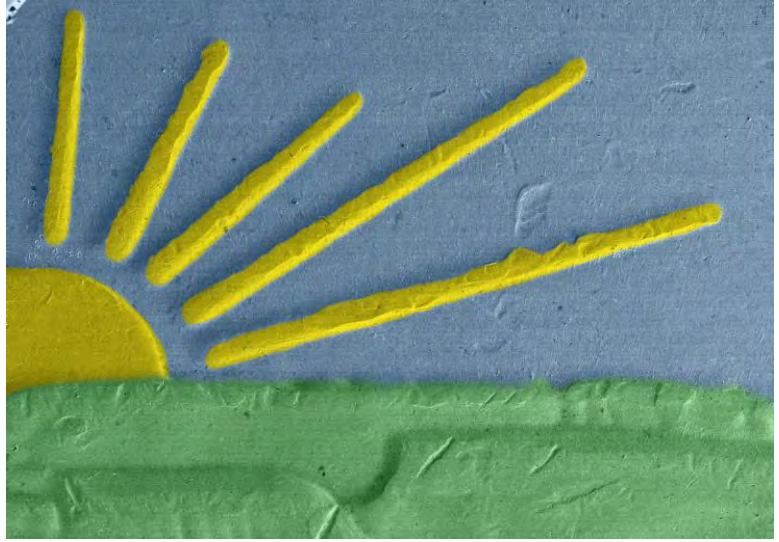

This is our image we did with the **Quick Selection Tool**. The edges are a bit rougher and more stark. With a quick glance it seems ok, but upon closer inspection it really doesn't look that great. As I said previously, it all depends on how much time you have and how much work you're willing to do.

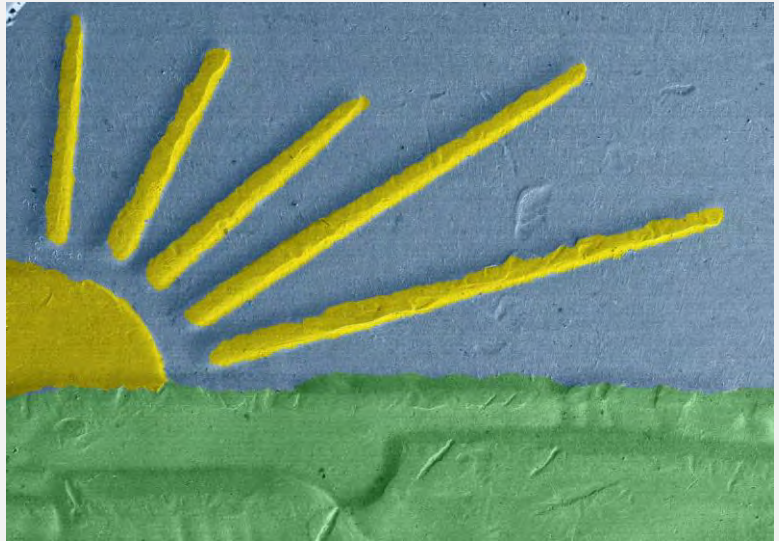

Lastly we have the image we combined using the **Lasso Tool**. It looks like the Feathering was set a little bit too high, as all the color edges look very fuzzy. There's also an area where the yellow sun and the green ground overlap that the blue sky is showing through. All these issues are fixable by backing up a few steps. This method with the **Lasso Tool** is probably the fastest.

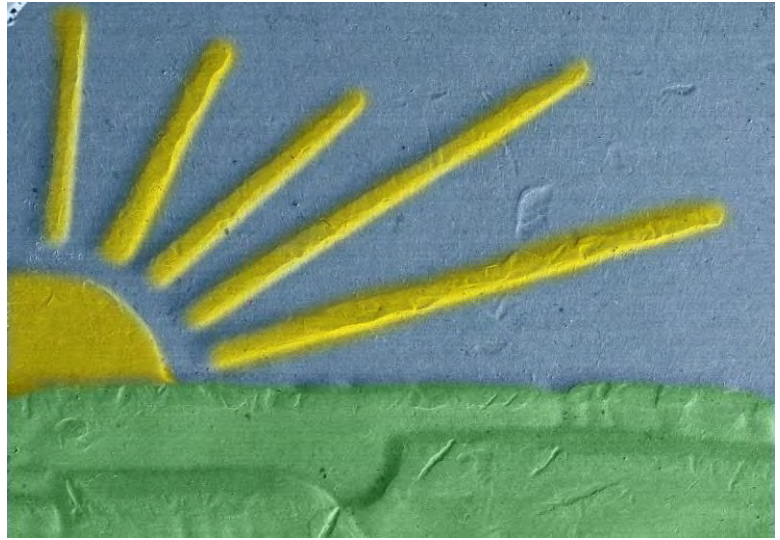

## VI. Pseudo-Color – The Easy Way


The pseudo-color technique described in the previous section I believe is the best way to colorize an image. However, it can be a little tedious, time-consuming and complicated. A large percentage of colorized electron micrographs I see are done using the technique that will be discussed in this section. I don't believe it will give you the quality of color that the previous section's instructions will. Your images look like one of those old-timey photos that someone painted to look color-y. They have color on them, but it looks a bit unnatural. However, this technique will definitely be quicker and simpler.

**1.** We'll start by opening our image and making a copy of it by clicking:

**Image > Duplicate**

Close your original image.


We don't really need to this, but it will prevent anything accidentally happening to our original image.

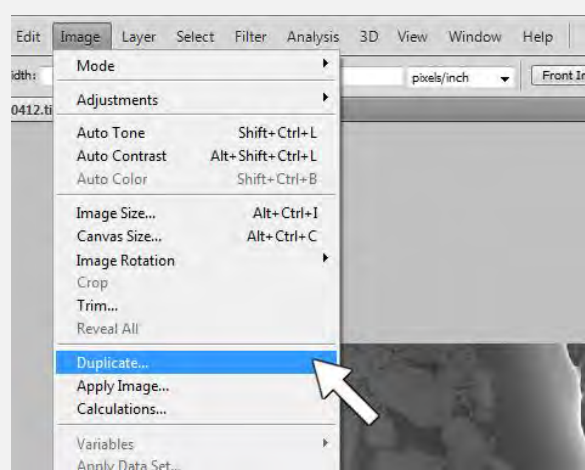

**2.** Perform the **Levels** operation on the image.


**Image > Adjustments > Levels**

or **Ctrl + L**

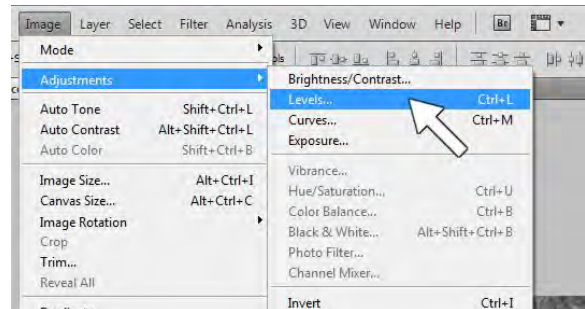

**3.** We will change our image to RGB Color.


**Image > Mode > RGB Color**

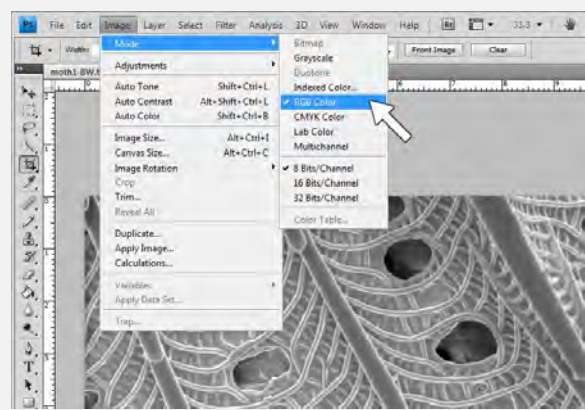

**4.** We will add a new layer to our image.


**Layer > New > New Layer**

or **Ctrl + Shift + N**

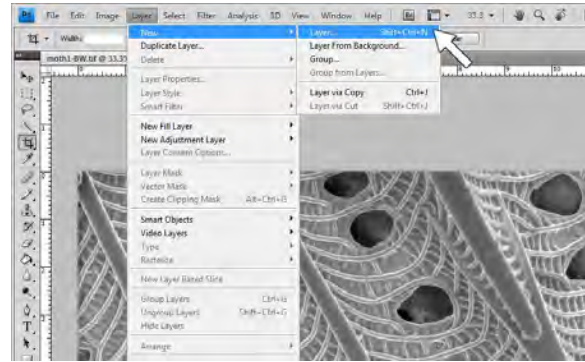


**5.** The layer name will come up as “**Layer1**” as default. We can keep that name if we want, or we could rename our layer to keep track of it better. In this case, I’m going to color some things red, so I named the layer “**Red**”.

Click **OK**.

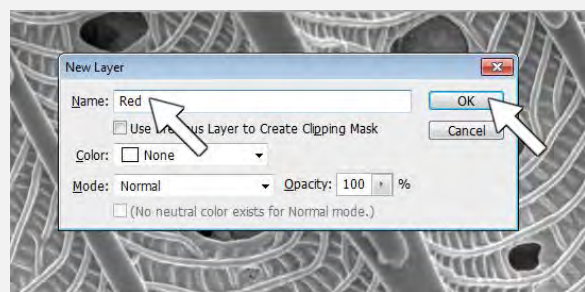


**6.** We’ll change the **Opacity** of this layer to around **50%**.

(When we made the new layer in the previous step, we could have set the **Opacity** to 50% then as well. It doesn’t matter.)

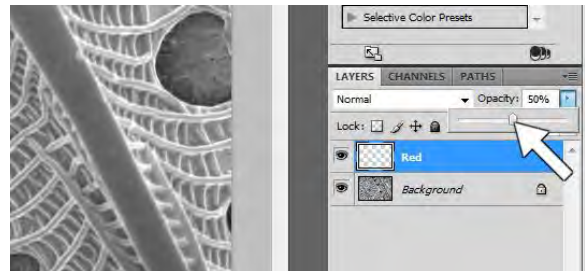


**7.** We’re going to choose our first color, which in this example will be red. So we’ll go to the **Color Picker** by clicking on the **Set Foreground Color** box (which by default is Black, with the background being White). Then we’ll choose our color by clicking on it.

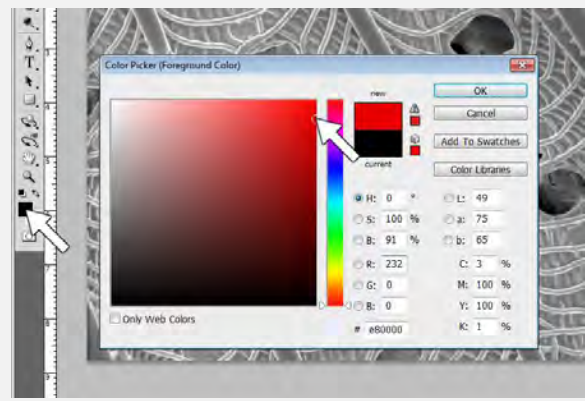

**8.** We're going to paint the areas we want red using whatever coloring tool we want. In this example, since most of the edges are going to be fairly sharp, I'm going to use the **Pencil** tool.

If the edges of my structure were less sharp, I would use the **Pencil** tool to color the center of the area, and use the **Brush** tool to go over the edges.

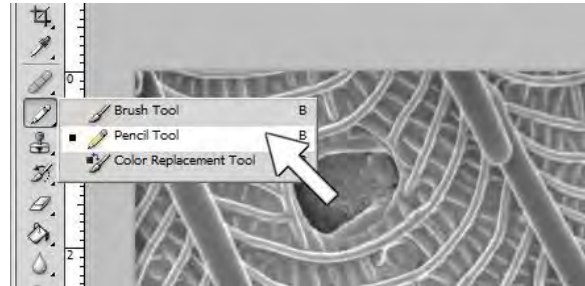

**NOTE:** If we accidentally draw outside the lines, we have 2 options to fix this.


**Option 1:** Stop drawing and undo your last action. Try:

**Edit > Undo** or **Ctrl + Z**

or


**Edit > Step Backward** or **Alt + Ctrl + Z**

**Option 2:** Select the **Eraser** Tool (**E**) to erase all your screw ups.

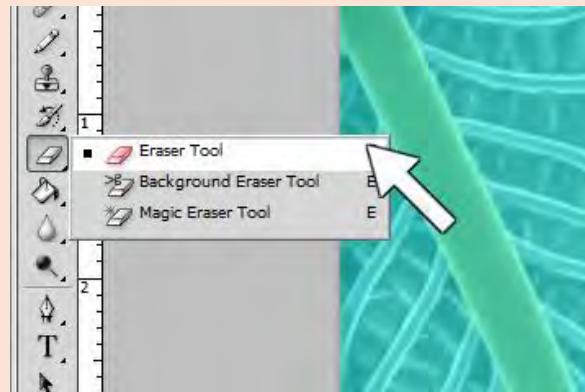


**9.** Once I have all the areas painted that I wanted, I'll move on to my next color, which I think I'll make a light orange color.


I can continue to paint my new color on this existing layer, but I like to keep things separate, so I'll make a new layer, repeating steps 4 thru 8 for the new color.

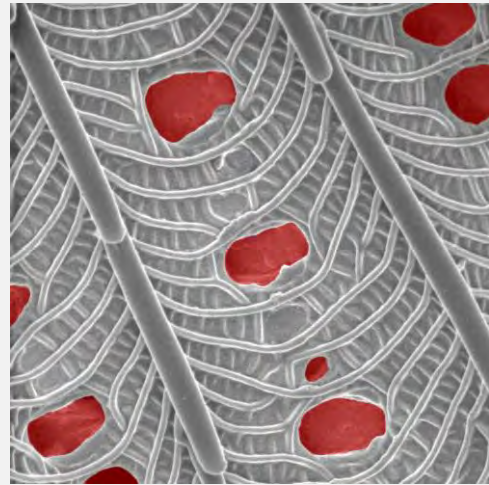


**NOTE:** Now that we have the 2 things we wanted colored done, we'll want to adjust the **Opacity** of our 2 colored layers so they look a little better / more subtle. 50% opacity is good to work with while painting, but for viewing it makes the color a little in-your-face. Usually something between 20-40% is more palatable.

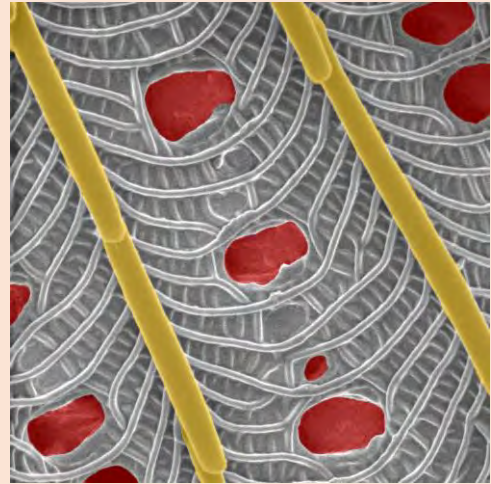

**11.** To change the **Opacity** of our color layers, make sure the correct layer is selected, then type a number in the box, or click the little arrow and drag the slider around.

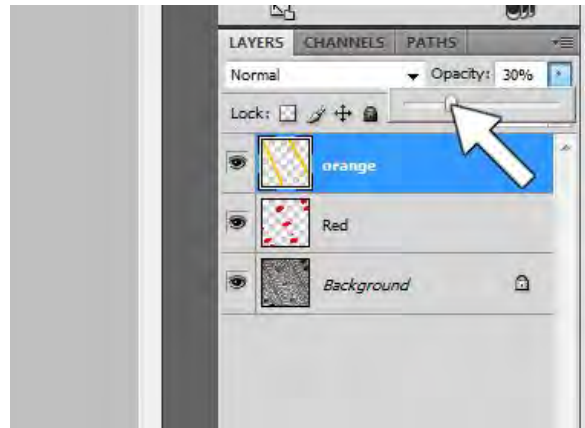


**12.** At this point, we could be done with our image, but check with the optional procedures below to see if you want to keep going.

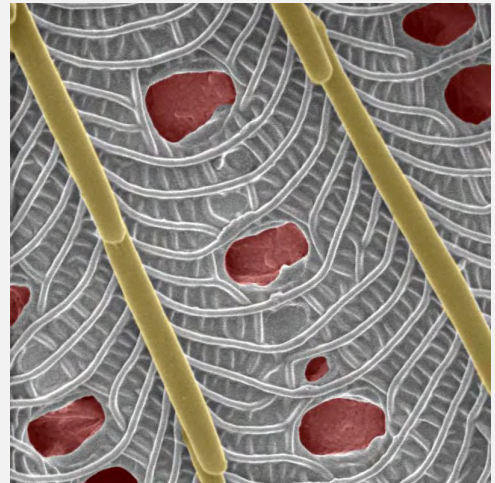


**OPTIONAL**

Most of the image is still in grayscale. This is one of the things I dislike about this technique. Some spots/areas of the image are in “sortof” color while most of the image will still be in grayscale.


We could make a new layer and try coloring all the rest of the image a new color. That presents 2 problems. 1. That’s a lot of work. 2. All the areas where any 2 colors meet might end up a bit messy. To quickly avoid messy edges, here’s something else we can try.

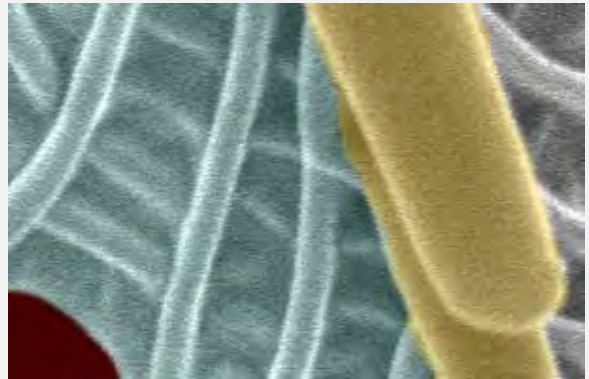

**13.** We’ll make a new layer (**Step 4**) and choose our color (**Step 7**). We’ll then select the **Paint Bucket** tool.

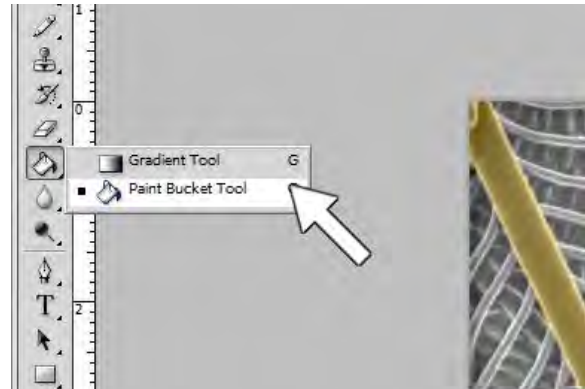


**14.** With our new layer selected, we’ll click anywhere on the image, and that should fill that layer with the color we chose. At 50% Opacity, the image looks washed-out and very blue, with our other previous colors being covered up and/or changed to some third color.

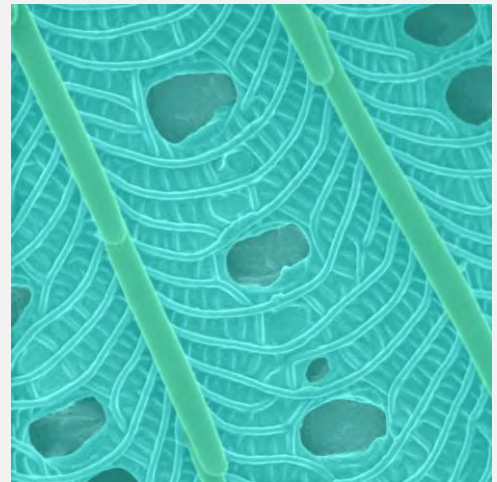





**15.** To fix this, start by reducing the **Opacity** of the new layer until the color is very subtle and will therefore cause minimal interference with our existing color layers. 10% is probably a good place to start (which is where I put the blue layer on this image).

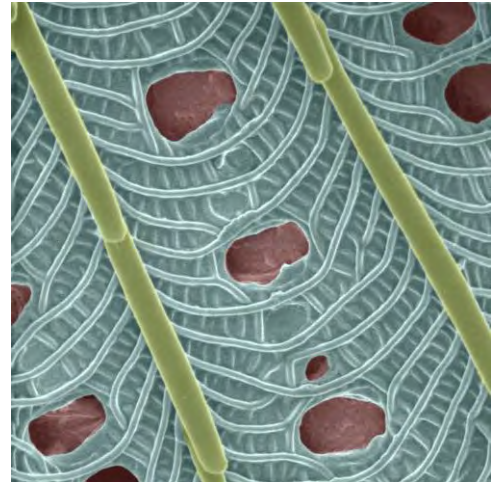


**16.** By doing some more tinkering, not only with the blue layer, but by increasing the Opacity on the orange and red layers, I'm left with a pretty ok image with some reasonable color applied to it. (Keep in mind this is just an example. I know these specific colors are goofy looking.)

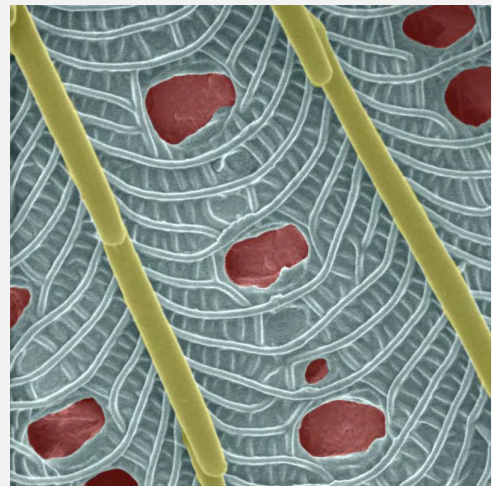


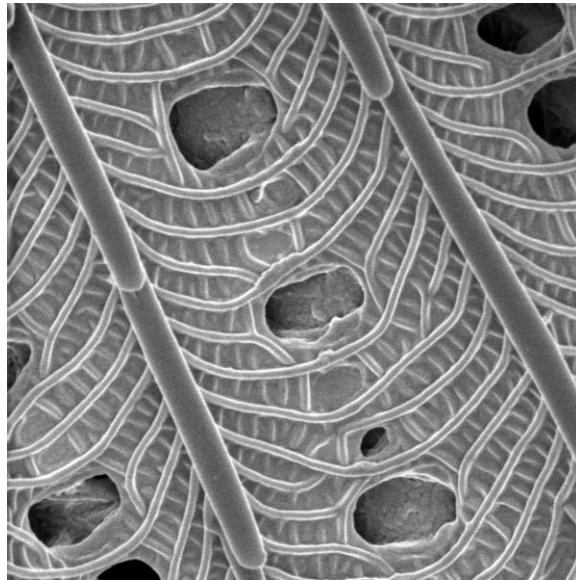

Our original B/W image with Levels adjusted.

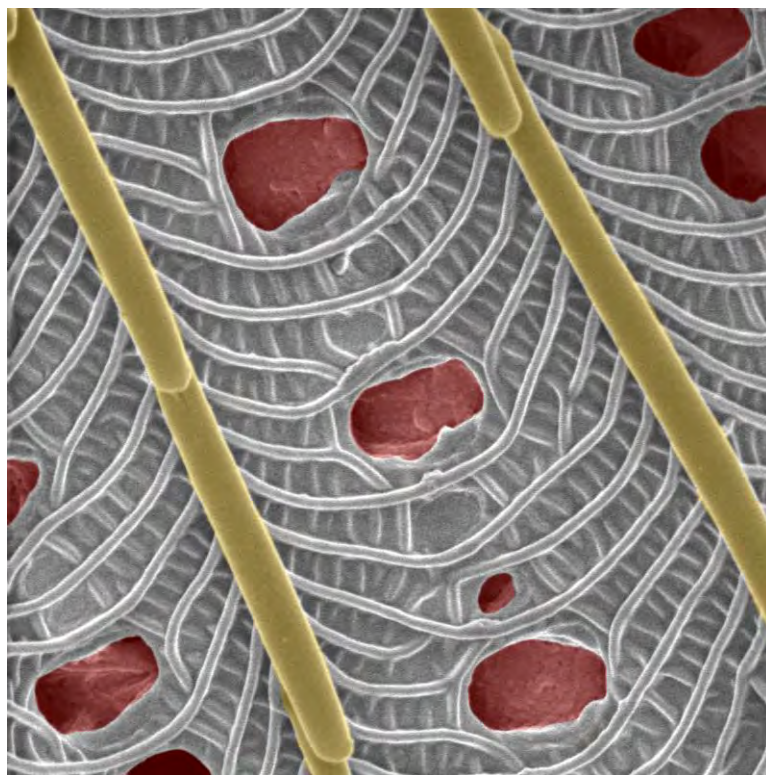

Our image with a B/W background and 2 colorized elements.

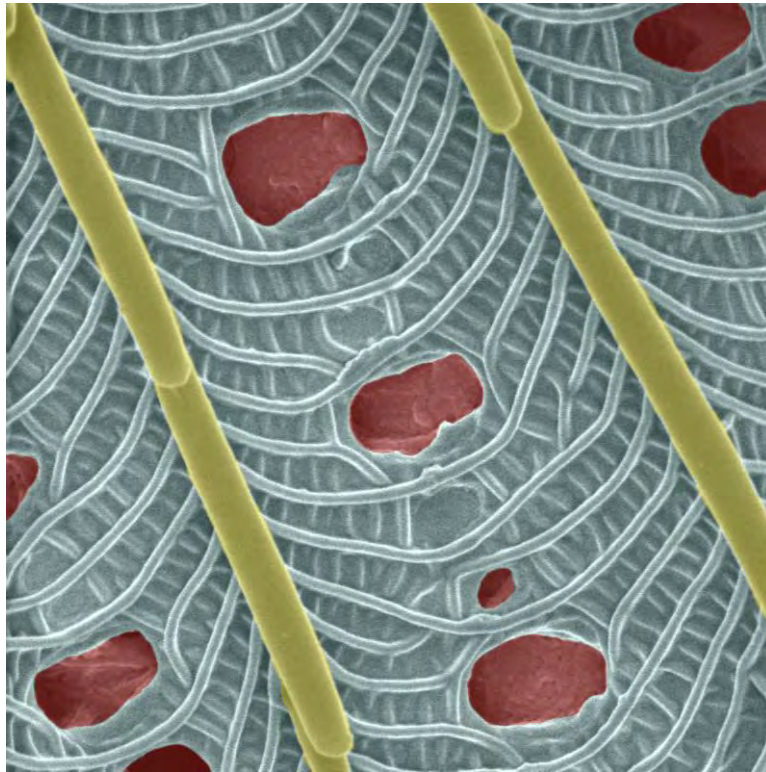


Our image with 2 colorized elements and one light color over everything.

## VII. Calculations

This process is solely for images that have uneven illumination. This is a common problem in many TEM micrographs and can also be encountered in low mag SEM images. The results of this technique can be very dramatic in repairing otherwise unsalvageable images.


**1.** Once your unprocessed image is open in Photoshop, make a duplicate of it by clicking:

**Image > Duplicate**

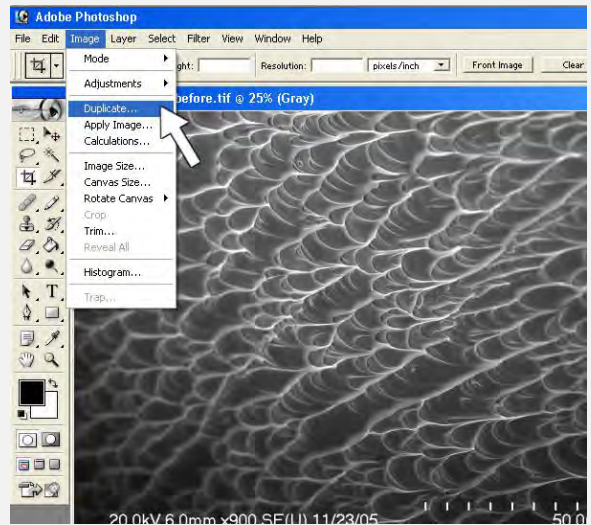


**2.** Now, make another duplicate by repeating the previous step. Close your original image. (You really only need 1 duplicate image for this, but by making 2 copies and closing your original, this ensures that nothing will accidentally happen to it.)

You should have 2 images named **filename\_copy** and **filename\_copy\_2**.

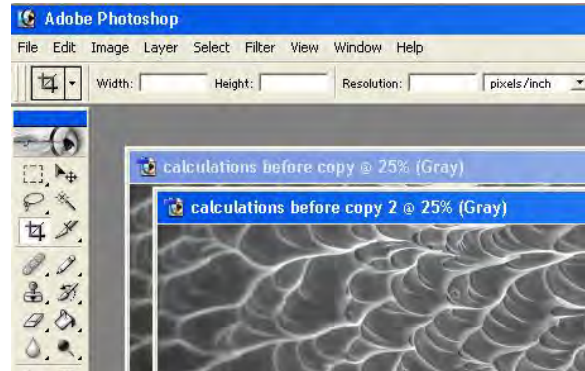


**3.** Make sure the **copy\_2** image is selected and click:

**Filter > Noise > Median**

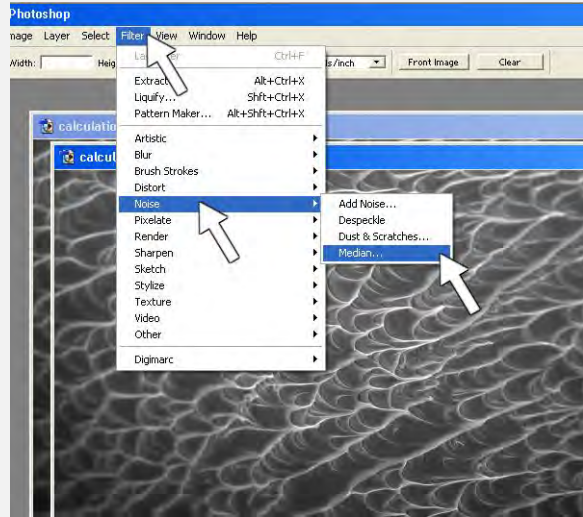


**4.** The **Median** dialog box will appear. Drag the slider all the way to the right so that the Radius is 100 pixels. Click **OK**.

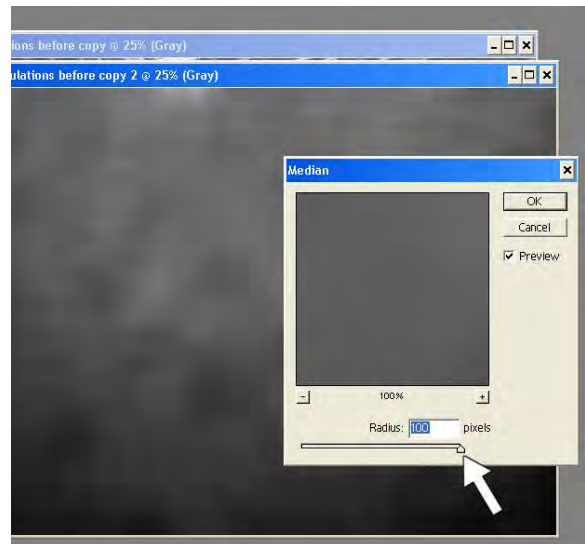


**5.** Click:

**Image > Calculations**

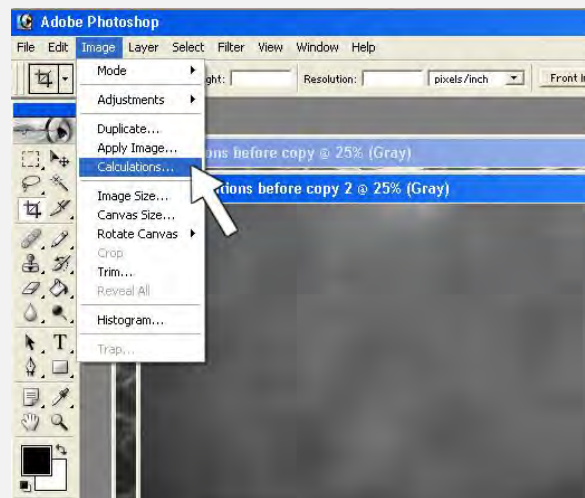



**6.** The **Calculations** dialog box will open. First, make sure the preview box is checked. Starting at the top, choose copy\_2 for your Source 1 image and check the Invert box. Next, make sure the first copy is your Source 2 image. Now, click on the **Blending** drop down menu and scroll through the different options until the illumination in your image appears the most even. In the case of the example it was the Vivid Light option. Finally under **Result** choose **New Document**. Click **OK**.

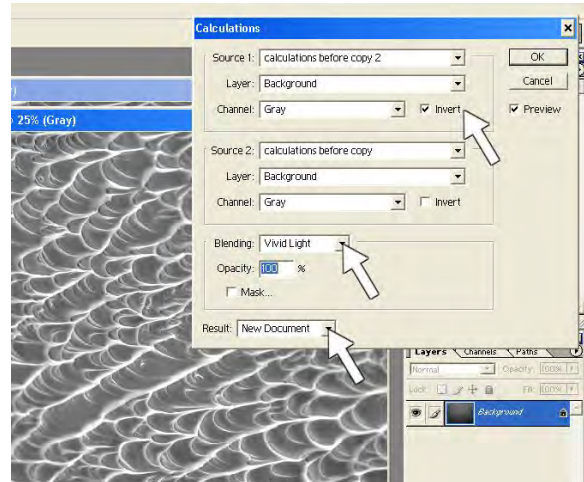


**NOTE:** This is a subjective, trial-and-error process. You may have to switch the images from Source 1 to Source 2 and go through all of the Blending options several times before you find the best result. You can even switch which image is inverted. (This rarely works though, but give it a try anyway. As a rule, the image that you performed the Median on should be inverted.)

**7.** You should now have a new third image. Now **and only now** should you finally adjust the **Levels**. Adjusting them at the beginning of the process can interfere with getting even results with Calculations.

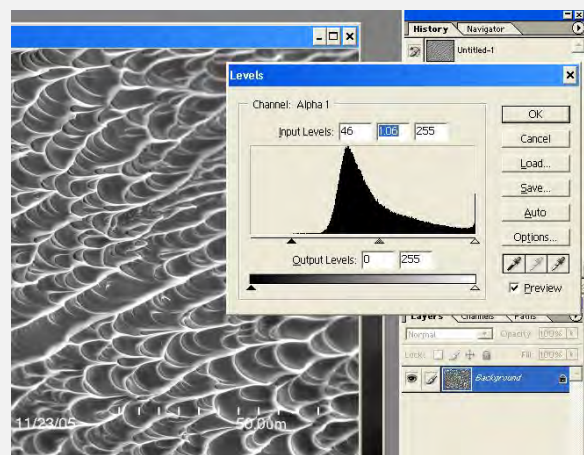

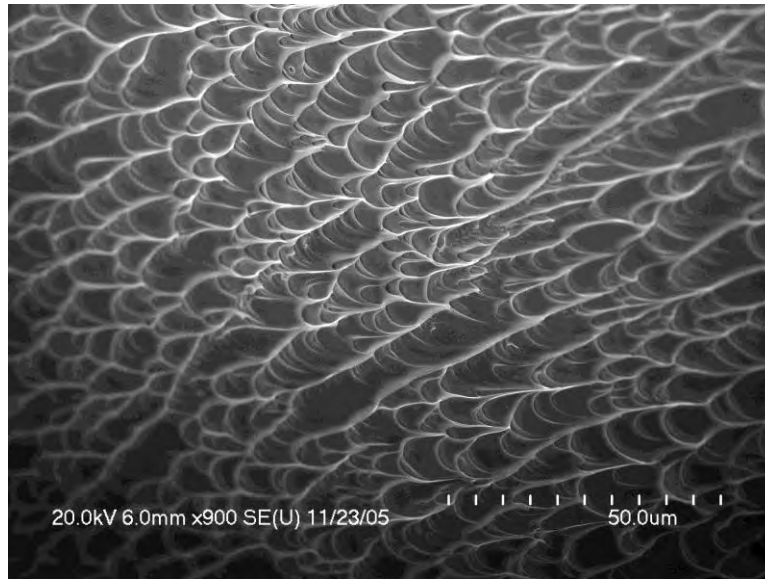

Here is the original image. There is a large amount of vignetting around the bottom of the image, nearly omitting information.

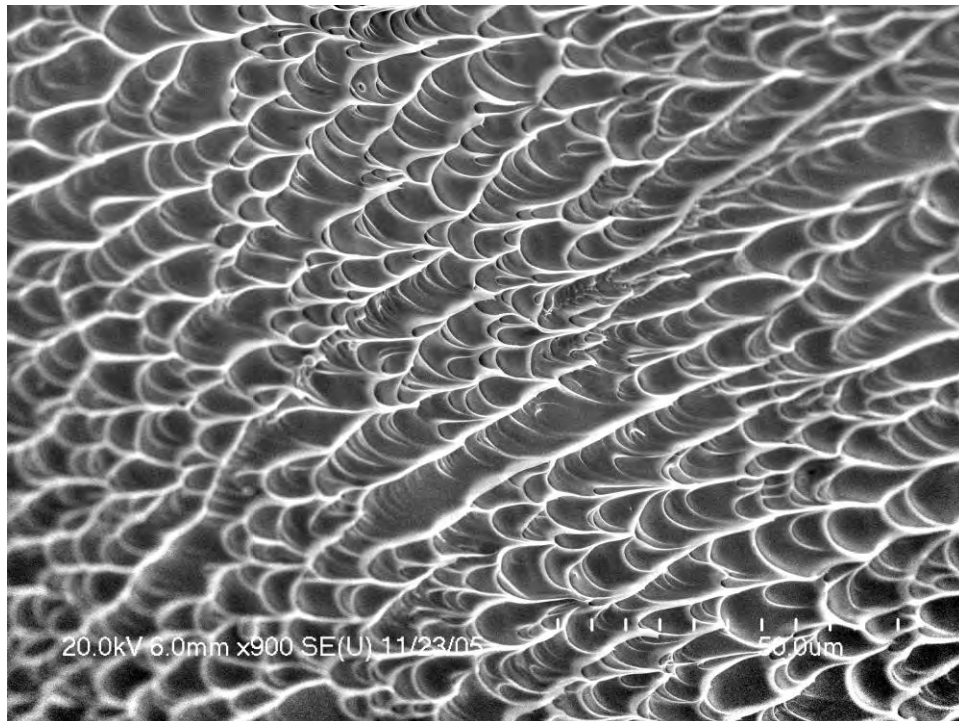

Here is the image after Calculations. The illumination is now much more even across the whole image.


## VIII. Combining EDS Color Maps - (Cut and Paste)

If you have a bunch of different EDS color maps and you want to combine them into a single image, this is one of the ways you can do that. This technique works best if your images do not have a lot of overlapping areas of color. (Most EDS software will be able to do this for you a lot better, but sometimes that's not an option.)

**1.** Open all your EDS color map images in Photoshop by clicking:

**File > Open**

or **Ctrl + O**

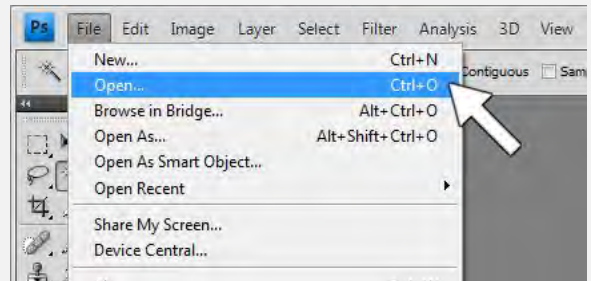


**2.** We're going to create a new image, but we want it to be the same size as the other images we just opened. By following this procedure, we ensure that they will be.

Select one of your opened color map images and click:

**Select > All**

or **Ctrl + A**

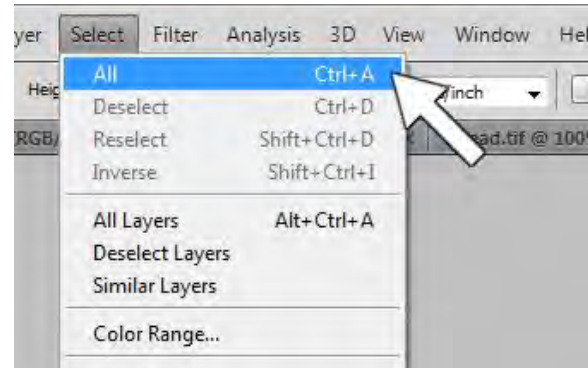


A moving dotted line should appear all the way around your image.

**3.** We now need to copy the image by clicking:

**Edit > Copy**

or **Ctrl + C**

We're not going to paste this image anywhere right now, but by copying it we've placed the image information on the clipboard. When we go to make a new image now, Photoshop will automatically make the new image the same size as the image on the clipboard.

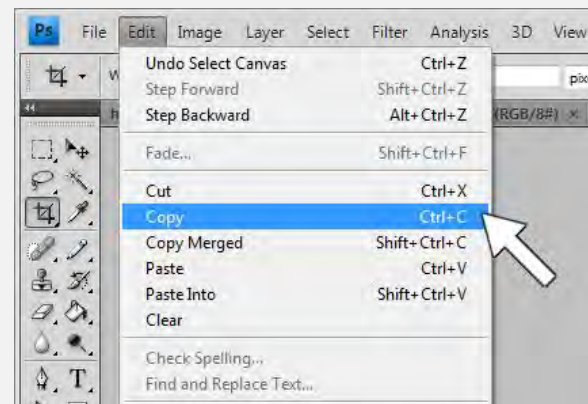




**4.** Now we can create a **New Image** by clicking:

**File > New**

or **Ctrl + N**

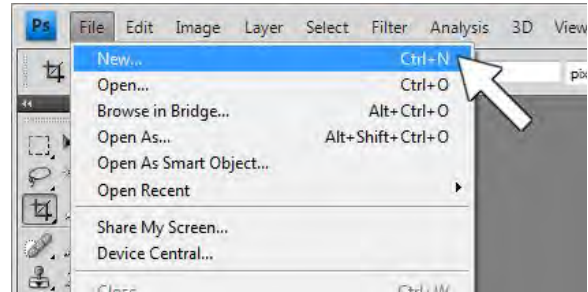


**5.** Your new image is probably white. Black makes a far better background for our color maps, so we will need to do a fill.


First we will need to make sure the **Foreground Color** is set to black by clicking the **Default Foreground and Background Colors** button, which is the small black/white box below and to the left of the actual foreground/background color display. Or hit the **D** key.




**6.** Select the **Paint Bucket** tool. (The **Paint Bucket** is often grouped with the **Gradient Tool**. If the **Paint Bucket** is not showing, you will need to right-click the **Gradient Tool** and select the **Paint Bucket**.)

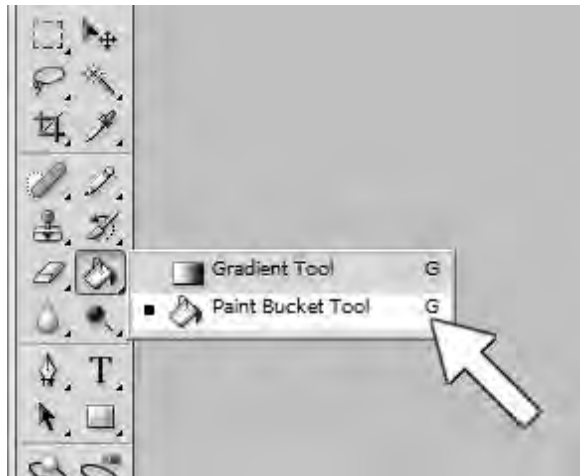

**7.** With the **Paint Bucket** cursor, click on your new image to make it black.

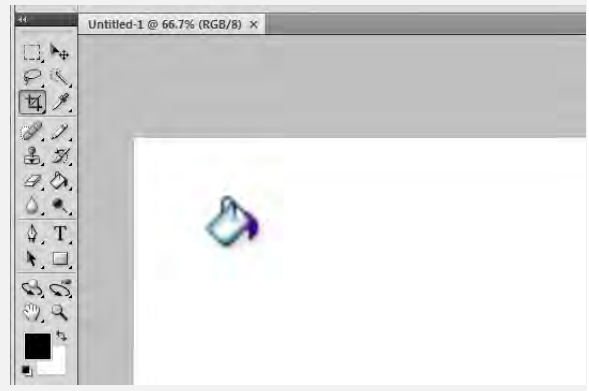


**8.** Select one of your color map images.

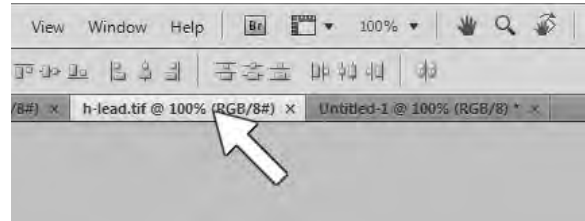


**9.** On the left side of the screen, select the **Magic Wand** tool. (In CS4 the **Magic Wand** is grouped with the **Quick Selection Tool**. If the **Magic Wand** is not showing, you will need to right-click the **Quick Selection Tool** and select the **Magic Wand**.)

With the **Magic Wand** selected, set the **Tolerance** to 0 by selecting the text box with the cursor and typing it in. Setting the Tolerance to 0 means that the Magic Wand will **only** select the exact color you click on. Sometimes that is too stringent. If that is the case, simply increase the Tolerance level until you get what you need.

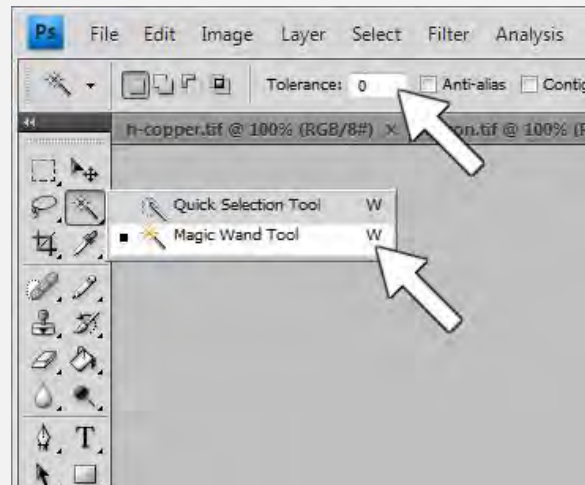

**10.** With the **Magic Wand** tool, click on an area in your image that is black.

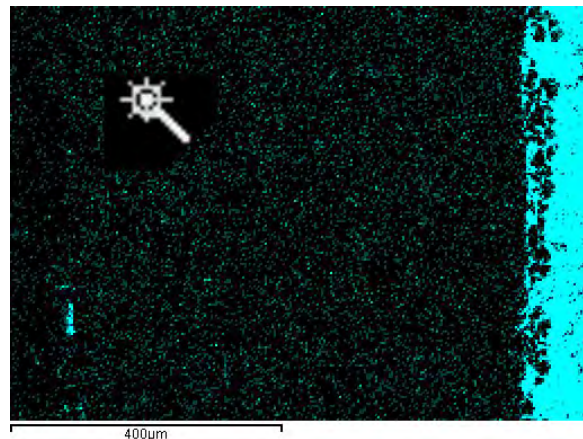

**11.** If your image is too small to tell if you're clicking on a black area or a color area, zoom in by selecting


**View > Zoom In**

or **Control (+) +**

You can do this multiple times. To zoom out, select


**View > Zoom Out**

or **Control (+) -**

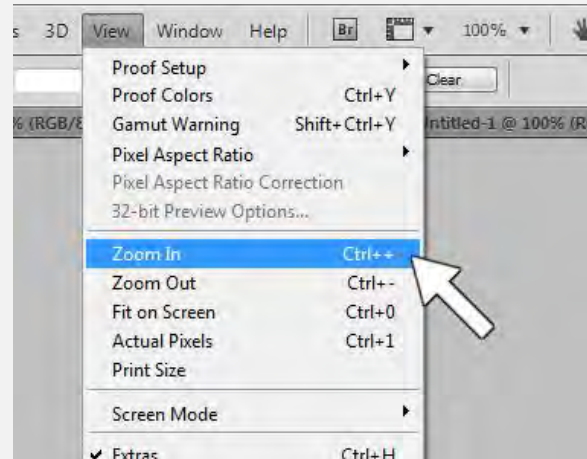

**12.** An animated dotted line should appear around all the black areas in your image now.


If it does not select **ALL** the black areas, hold down the Shift key and select any areas that were missed.

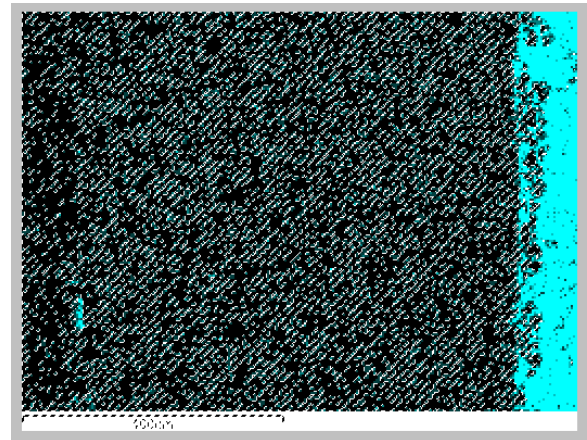

**13.** In reality, we don't care about the black part of the image, we want the color selected. So, we will need to invert the selection by selecting

**Select > Inverse**


or **Shift + Ctrl + I**

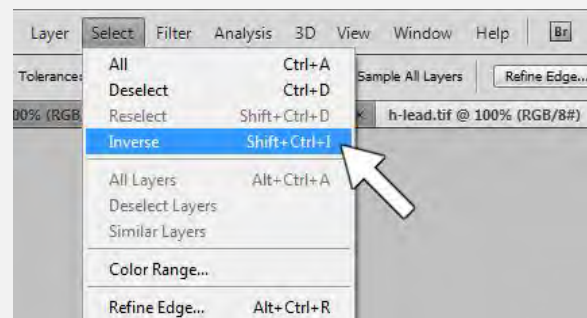

**14.** We now want to copy the color portion of the image by clicking on

**Edit > Copy**

or **Ctrl + C**

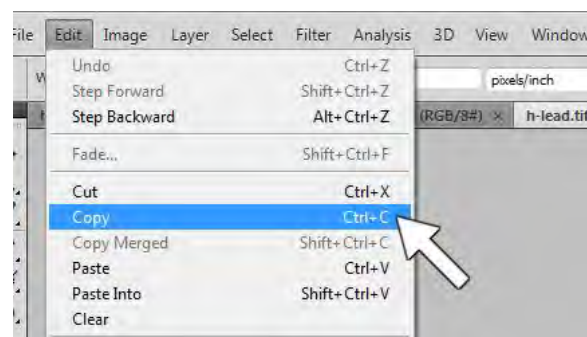


**15.** Select your new image and **Paste** the copied image into it by clicking


**Edit > Paste**

or **Ctrl + V**

Repeat steps 8 through 15 on your remaining images. I have 3 in my example.

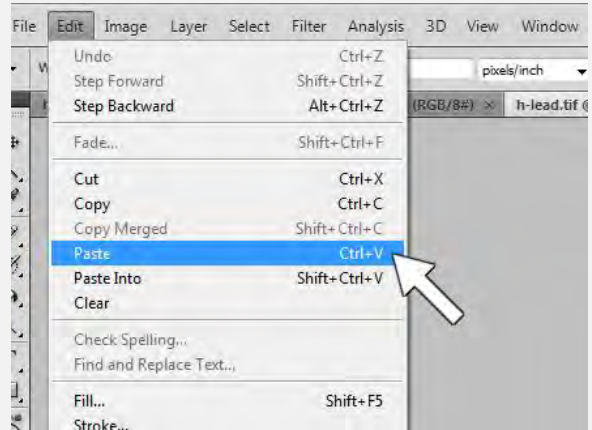


**16.** You will now need to arrange your layers to find the best order. You can order your layers by dragging them up and down in the Layers menu in the lower-right part of the screen.

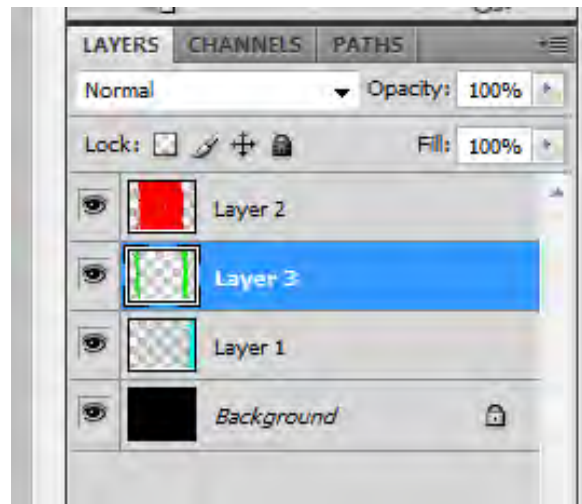


**17.** What does this do? For example, I know that my blue layer has a lot of extra noise across the rest of the image. Here is what it looks like with it in front.


Surprise! It looks noisy.

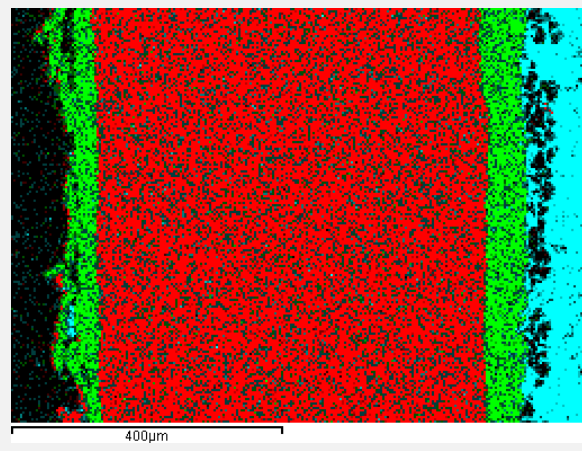


**18.** By putting my red layer on top, we don't have to see all that blue noise. (You might even want to paint the noisy areas black so you're left with just the large blocks of color. [Do this with the **Pencil** tool before **Step 10.**] The noise is not important information in my image, but it might be in yours. The point here is, do what you need to do to get the information across in the best way possible.)

You're done, **Save** it!

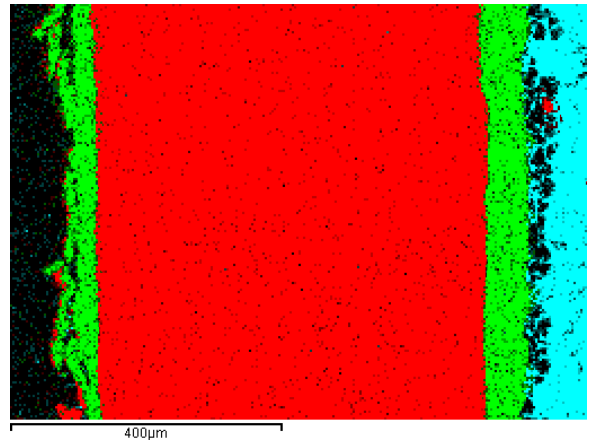

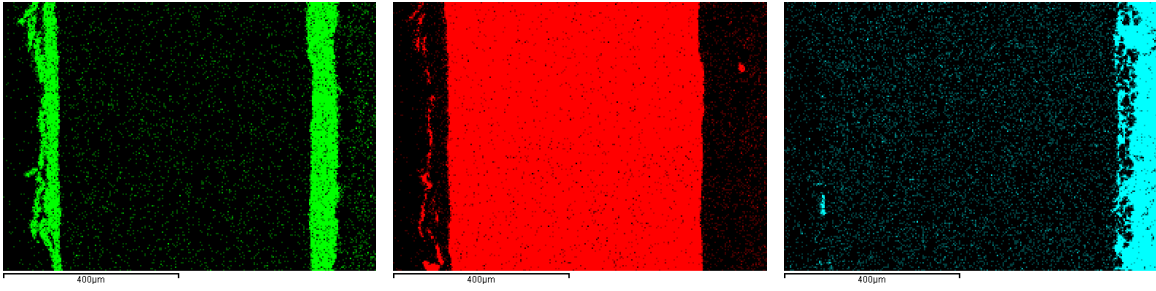

Here are my original 3 images.

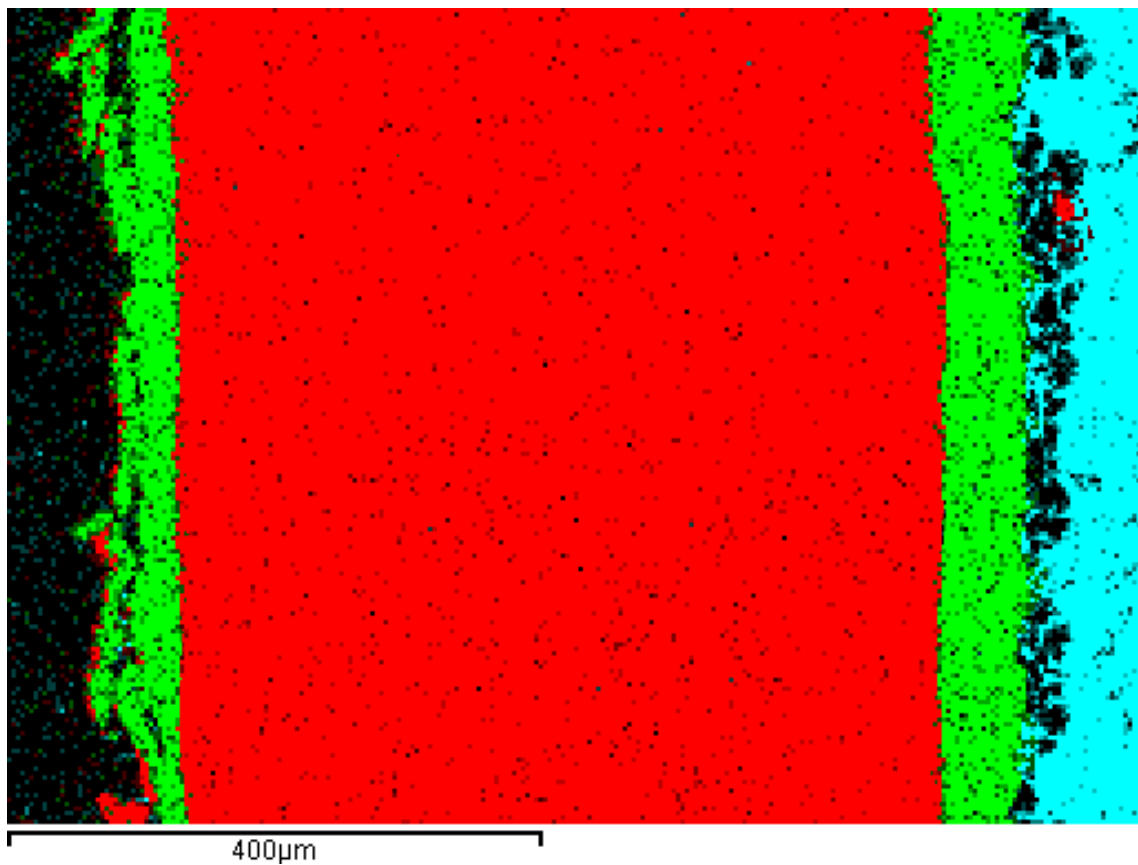

Here is the final combined image. Is it great? Not really. The best way to go is with the overlaid images generated by your EDS software. However, some of us are using some crazy old software, and sometimes we only have the EDS images to work with, and not the EDS software.

## IX. Multi-Detector Color

Let's say we took several images with different detectors (or different segments from a backscatter detector) *of the exact same area*. We can take those different images and combine them with color to get a very special kind of image. Basically, what we're doing here is making each detector a "light source" and assigning a particular color to it. It will make sense by the time we're done.


If you only have 1 SE detector, this procedure **IS STILL POSSIBLE!** Check out the ***Tips, Tricks and Suggestions*** section at the end of this section.

**ATTENTION:** This may look like a long and complicated procedure at first, but it's really quite simple and once you're familiar with it, can be done very quickly. The time-consuming part is really taking the images.

**1.** Open all your images in Photoshop by clicking:


**File > Open**

or **Control + O**

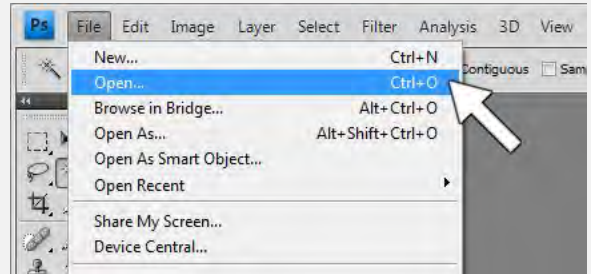

**2:** In this example, I've got 1 Secondary Electron Image and 3 Backscatter Electron Images (taken from different segments of a 4 segment detector).

Sometimes images from different detectors won't line up everywhere across the image exactly. If you have a backscatter detector with 2 or more segments, you can take all your images from it, which this makes this process a bit nicer (because unless the sample is drifting, they should line up pretty well).

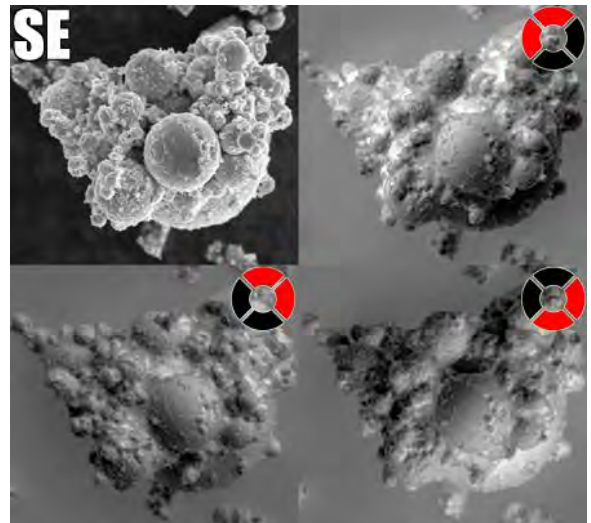

**NOTE:** Theoretically we only need 2 images to make this technique work. But I always like to have the most options available later on for processing, which means taking more images. (See the next section for tips on taking images for this technique.)


**3.** We're going to create a new image, but we want it to be the same size as the other images we just opened. By following this procedure, we ensure that they will be.

Select one of your opened images and click:

**Select > All**

or **Ctrl + A**


A moving dotted line should appear all the way around your image.

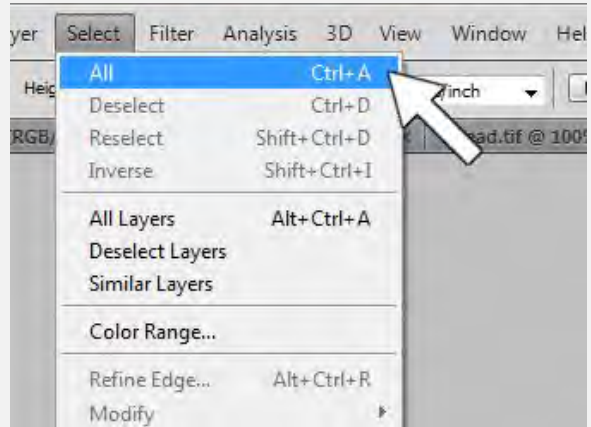

**4.** We now need to copy the image by clicking:

**Edit > Copy**


or **Ctrl + C**

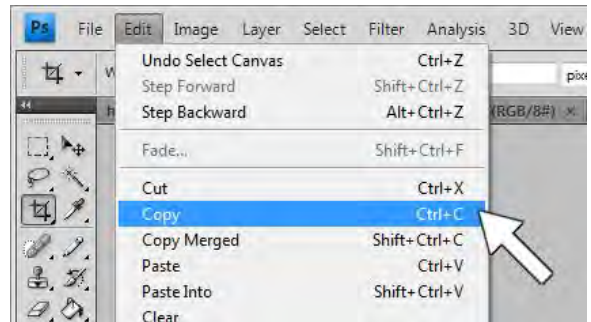

**5.** Now we can create a **New Image** by clicking:

**File > New**

or **Ctrl + N**

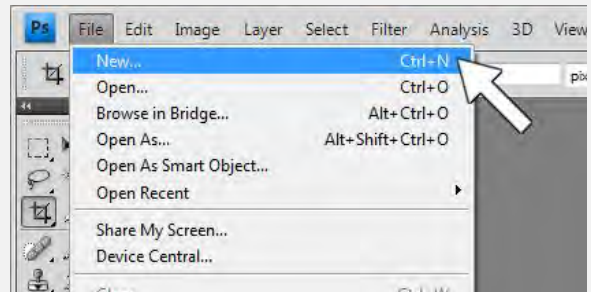




**6.** Your new image is probably white. Though it might not be necessary, black makes for a better background for what we're going to do, so we will do a fill (just in case).

First we will need to make sure the **Foreground Color** is set to black by clicking the **Default Foreground and Background Colors** button, which is the small black/white box below and to the left of the actual foreground/background color display. Or hit the **D** key.




**7.** Select the **Paint Bucket** tool. (The **Paint Bucket** is often grouped with the **Gradient Tool**. If the **Paint Bucket** is not showing, you will need to right-click the **Gradient Tool** and select the **Paint Bucket**.)

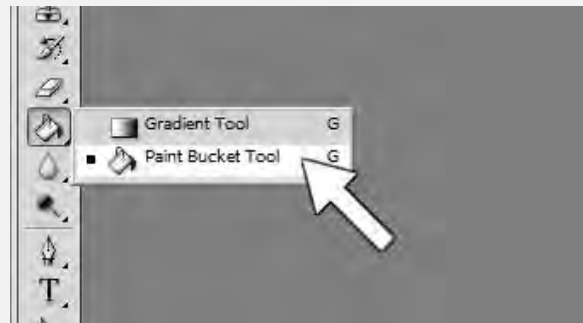

**8.** Make sure the **Opacity** for the **Paint Bucket** is set to 100%.


With the **Paint Bucket** cursor, click on your new image to make it black.

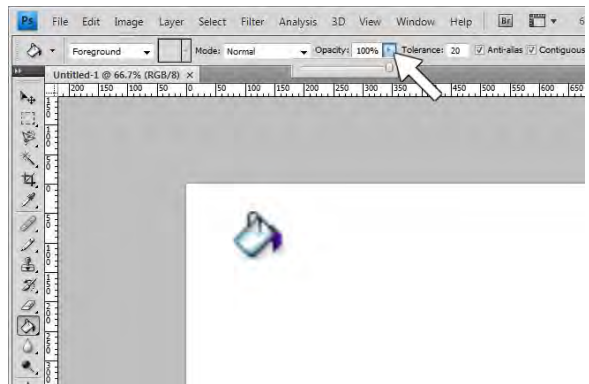

**9.** Once the new image is created and has a black background, we'll then paste the image we copied back in **Step 4** into it by clicking:

**Edit > Paste**

or **Ctrl + V**

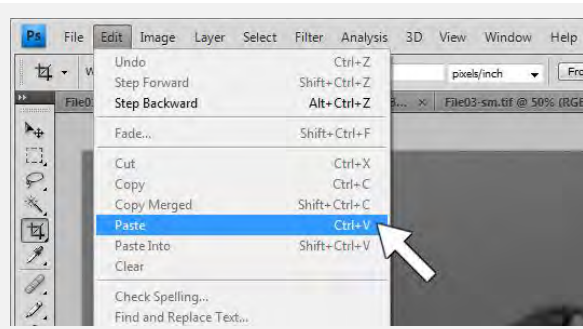

**10.** We then want to copy and paste the rest of our images into this new file, just as we did the first one.

I'm using 4 images in this example, so I have 4 image layers and a background layer.


When you're done copying all your images, you can close them, so that just your new image is open.



### **A1. Special Short Cut**

Instead of opening all our images and copying them and pasting them into a new image, we can run a script (which is already loaded in PS) to do all that for us.

**File > Scripts > Load Files into Stack...**

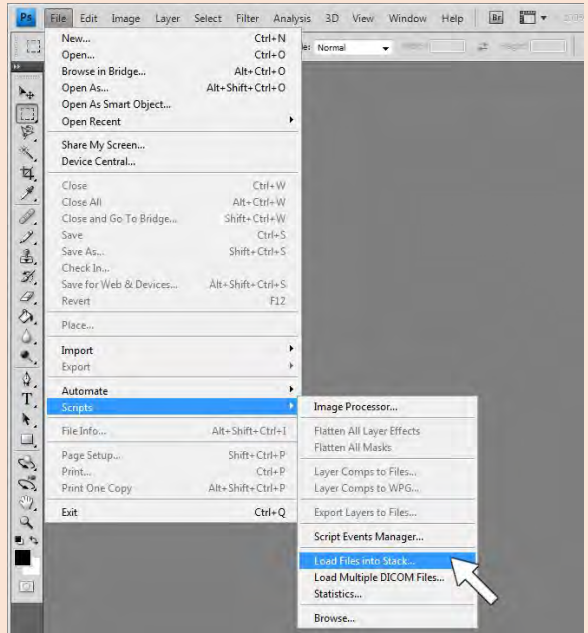


**A2.** The **Load Layers** window will open.

Click **Browse** to find and select the images we want to combine.

\*As an extra bonus, we can check the **Attempt to Automatically Align Source Images** box. This will hopefully save us from doing this manually in steps 11 thru 16.

Click **OK**.

(Your computer may take some time to think about this.)

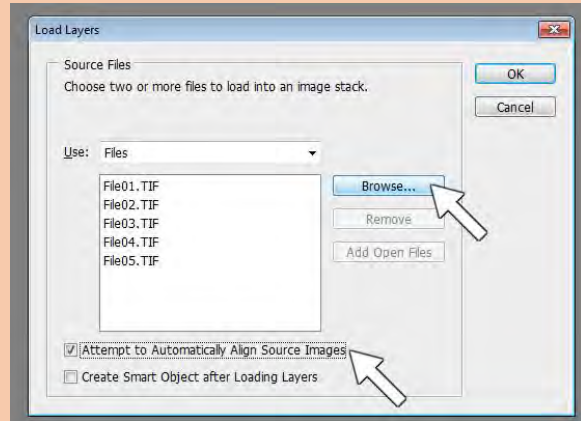

**A3.** We now have a new image with all our images in separate layers (and in the order that they were selected in the **Load Layers** window in the last step).


We have no Background layer however, and we will want to make one by clicking:

**Layers > New > Layer...**

Or **Shift + Ctrl + N**


Make sure this new blank layer is at the very bottom of all the other layers and perform Steps 7 and 8 on your new background layer.

**11.** Before we go on, we want to make sure all these images are aligned with one another. First we need to choose our central image that we will align all the other ones to. In this case, I'm going to use the Secondary image labeled **Layer 1**.

Start by deselecting the eye icon  next to layer 3 and 4.


Now we don't see those layers anymore and the only one we can see is **Layer 2**, since it's the topmost layer that's visible.

**NOTE:** If you don't need to align these images, skip to **Step 15**.

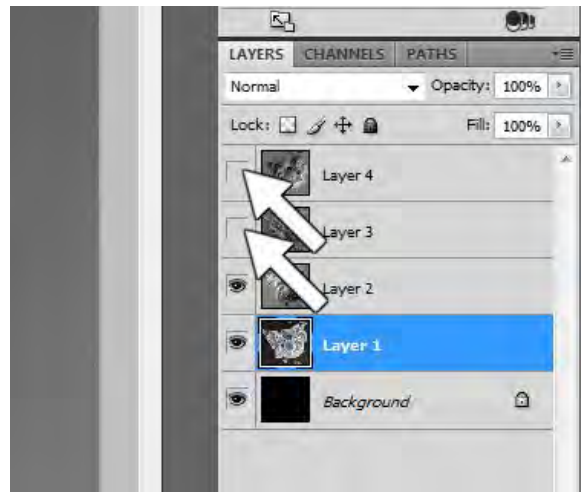


**12.** Select Layer 2 by clicking on it and then click on the **Opacity** box and change it to somewhere around 50%. (You can drag the slider bar around that appears when you click the **Opacity** box **OR** just type in a number.)

We should now be able to see Layer 1 through the halfway transparent Layer 2.

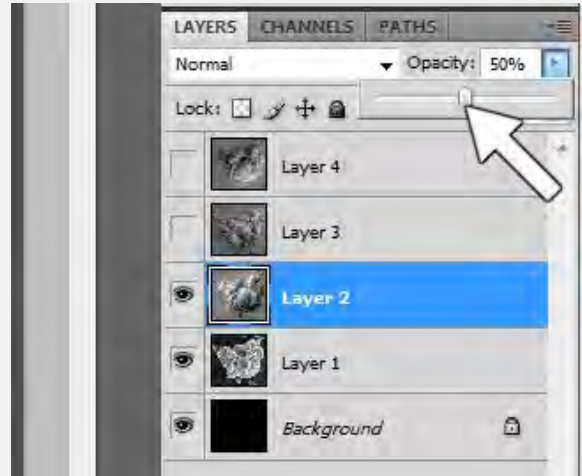

**13.** If you take a look at your image and you can see areas that look fuzzy/out of focus, or otherwise obviously mis-aligned, we'll need to line them up.

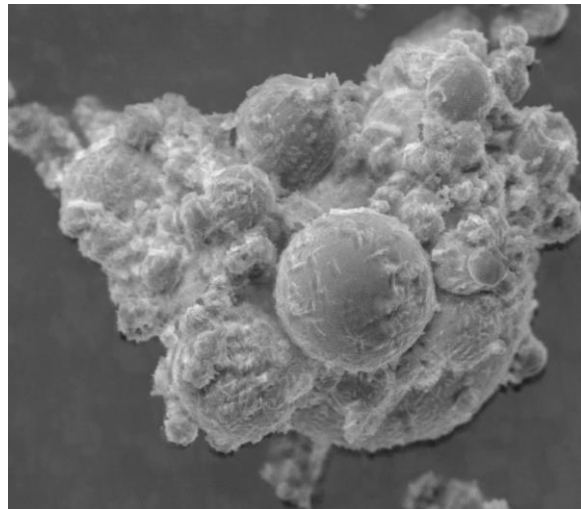

**14.** With Layer 2 still selected, select the **Move Tool** on the toolbar (or **V** on the keyboard). This will allow you to click on the layer and drag it around. (You can also move the image around by 1 pixel at a time by using the arrow keys.)


This is where you may notice that some areas of your image will line up while others won't. If that is the case, you'll need to choose a main area of interest in the image that you'll align everything to.




**15.** So once we have that layer aligned, we will **leave the Opacity at 50%**. Now we'll align the next layer by deselecting the eye in Layer 2, and reselecting it in Layer 3 (by clicking in the box where the eye should be).

**NOTE:** Make sure you have the layer you want to move selected before you start dragging things around.

You'll just repeat this process on each of the layers until they're all aligned.

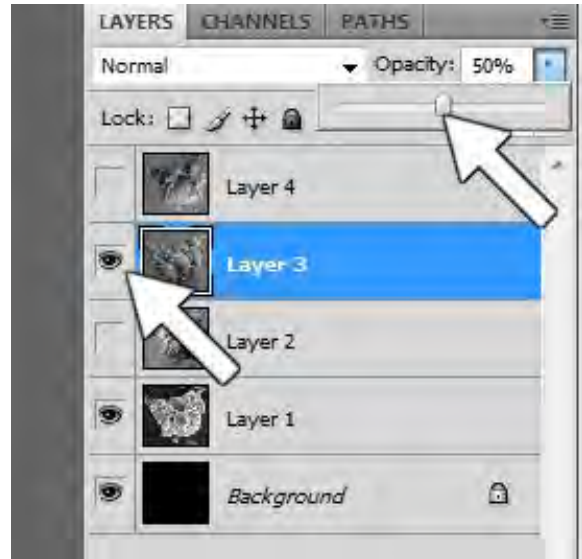


**16.** Once we're done with all of our aligning, we'll make sure that all the layers are visible.


(Meaning the eye icon  is showing next to each layer. If it's not, just click on the empty box so it appears.)

**17.** Make sure your image is RGB Color and not Grayscale.

**Image > Mode > RGB Color**

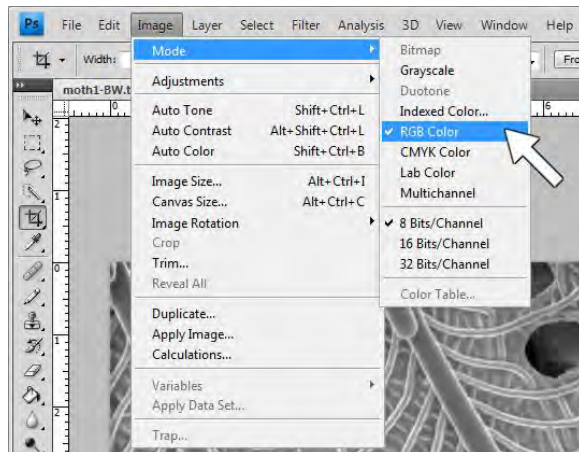


**18.** Now we move on to the real part. We'll start with Layer 4. If we right-click on it we'll get a huge pop up menu and we'll select **Blending Options**.

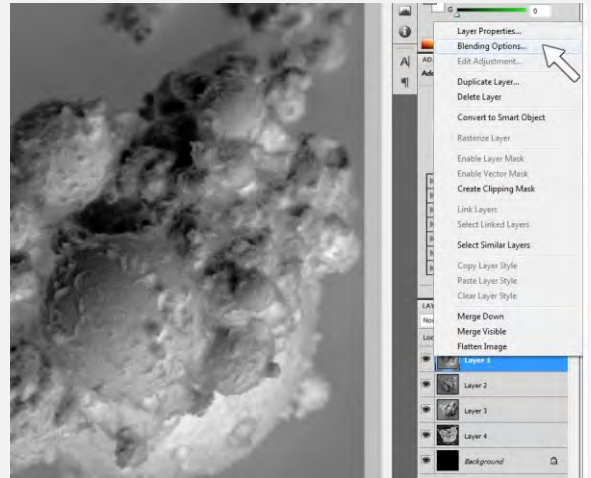

**19.** **Blend Mode** should be **Normal**. (You can change this if you want, but for this type of operation, **Normal** is probably best.)

If you didn't need to align the images earlier, you'll need to change the layer's **Opacity** to **50%** now. (Higher **Opacity** gives that layer more intense colors, Lower **Opacity** gives more muted colors.)


What we're interested in now are the 3 **Channels** checkboxes in the middle of the window. We're going to set Layers 2, 3 and 4 to be different color channels.

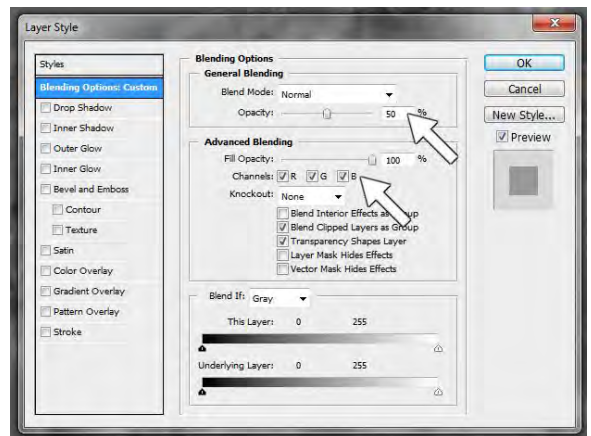

**20.** What we will do is uncheck 2 of these boxes. For this layer, I'm only going to leave the **R** box checked.

You should notice at this point that your image went from B/W to 2 crazy colors.

Click **OK**.

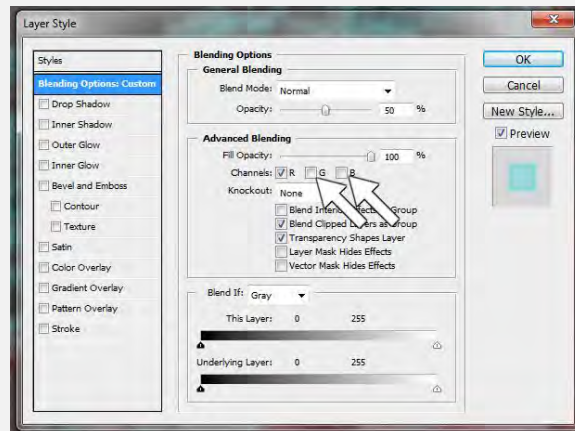


**21.** We'll then repeat this process with Layer 2 and Layer 3. Only with Layer 3 we will uncheck **R** and **B**, and with Layer 2 we will uncheck **R** and **G**.

Now, does it matter which layer has which boxes unchecked? **No**, it will be up to your personal preference.


What if I wanted to only uncheck 1 box? **Do it!**

All that is important is that these 3 layers have a different combination of Channels selected.

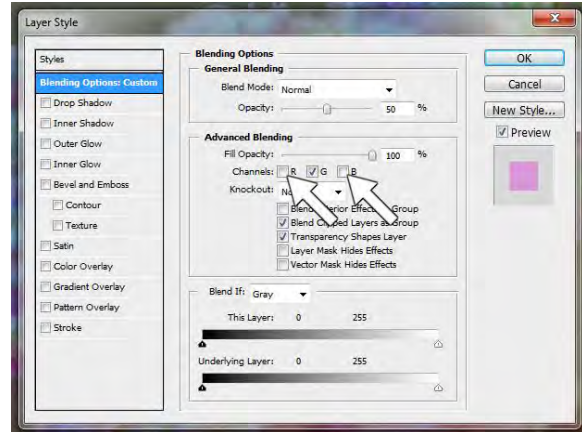


**22.** If we take a look at our image now, it looks a little different than how it started. It appears as though differently colored lights are shining at the sample from different directions. Neat!

For some people (and depending on how you took your original B/W images), this image is perfectly acceptable, and we could stop now if we wanted to. **But** there are a couple other things we *might* want to do to get the image even better.

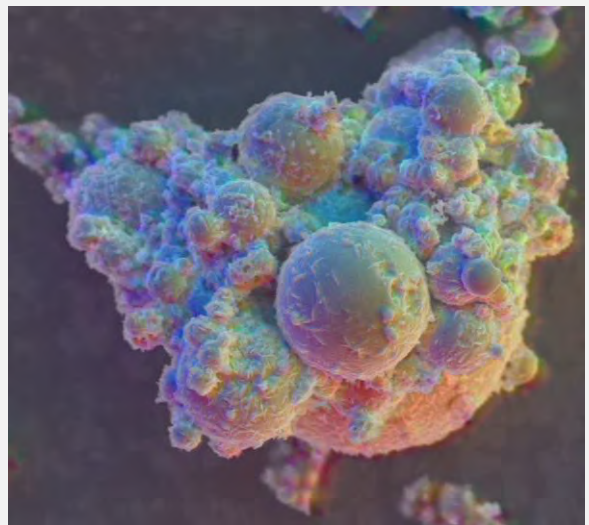

**23.** Sometimes the colors on the image come out a little... extreme. We're going to tone it down a bit. We'll start by dragging our still B/W SE image (Layer 1) from the bottom of the pile of layers to the top, and reducing its **Opacity** to **50%**.

50% Opacity is not a hard and fast number. If you want your colors MORE dramatic, use a LOWER number, like 20%. If you want your colors more subdued, use a higher number like 80%. **It's up to you!**

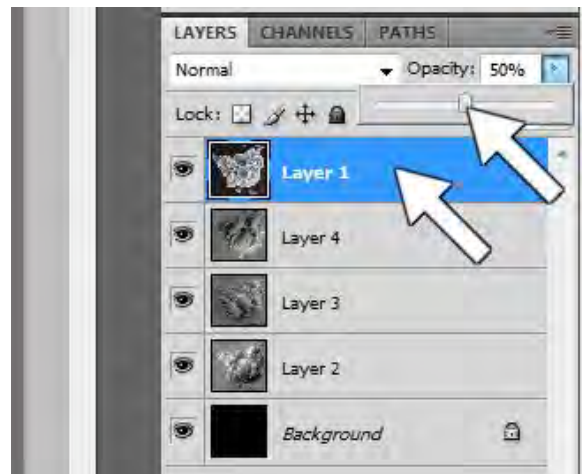


**24.** We should now have a really sad looking, washed-out image. We'll now save this image as a new file.

**File > Save As...**

or **Shift + Ctrl + S**

We want a flattened regular image (.JPG, .TIF, .BMP, .GIF, etc.) without any layers.

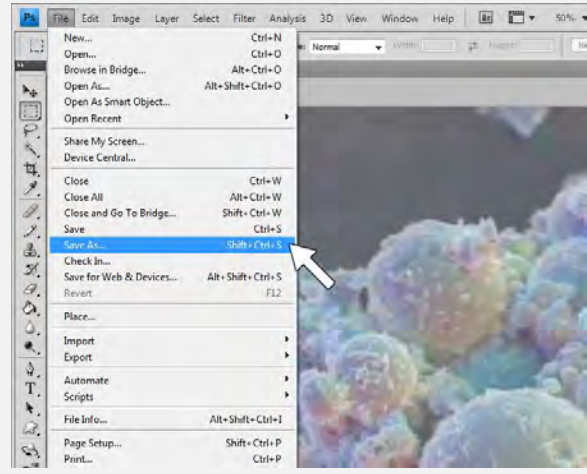

**NOTE:** If you save the file as a .TIF, PS will ask you if you want to keep the layers. You can tell it **NO** 2 different ways. First, in the **Save As** dialog box, uncheck the option that says **Layers**. **OR** after you hit **Save**, select the **Discard** radio button under **Layer Compression**.


If you **Save As** a .GIF, PS will ask you if you want to flatten the image. Say **YES**.



**25.** We'll now leave our image with all the layers alone for now and we'll open the new flattened image we just made.

**File > Open**

Or **Ctrl + O**

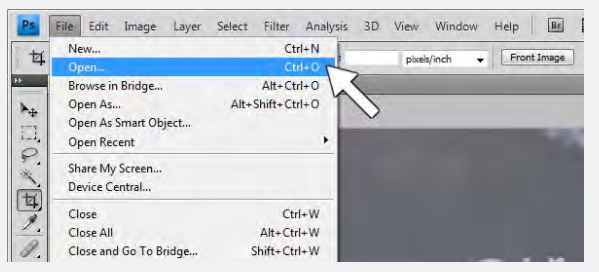



**26.** Now, once we have a flattened, 1 layer image, we will do the Levels operation as described in the first section.


**Image > Adjustments > Levels**

or **Ctrl + L**

Now because this is a color image you can also try the auto functions, they usually work great (they're not always so hot with B/W images).


**Image > Auto Color**

or **Shift + Ctrl + B**

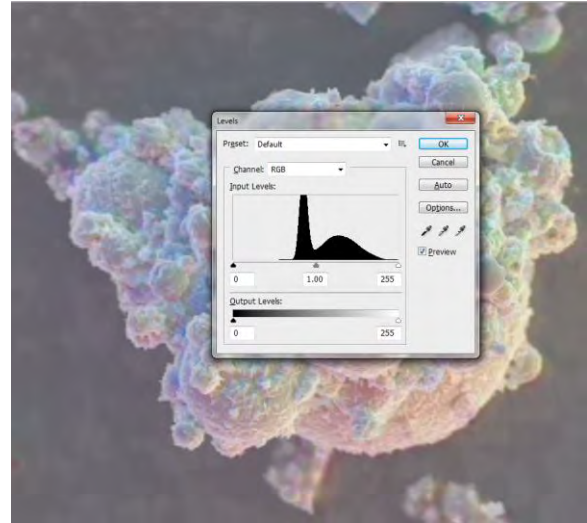


You should be done! **Save it!**

We may also want to go back to our first multi-layered image and Save it As a Photoshop file (.PSD). This will save all the layers and stuff so if we want to go back later and change things up (colors, opacity, intensity) we can.

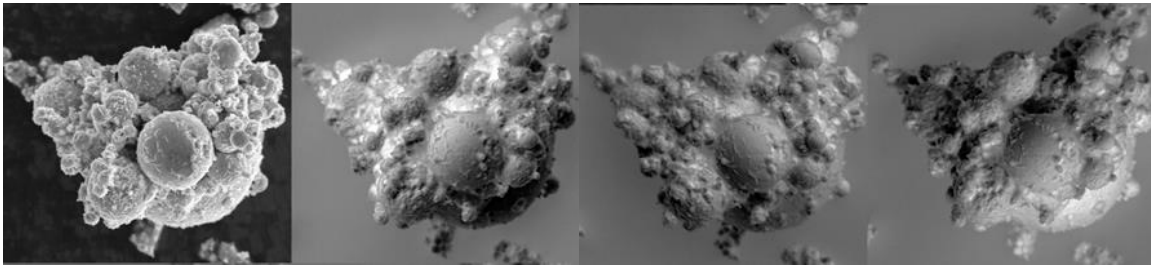

We started off with 4 regular old B/W micrographs taken of the exact same area, but with different detectors.

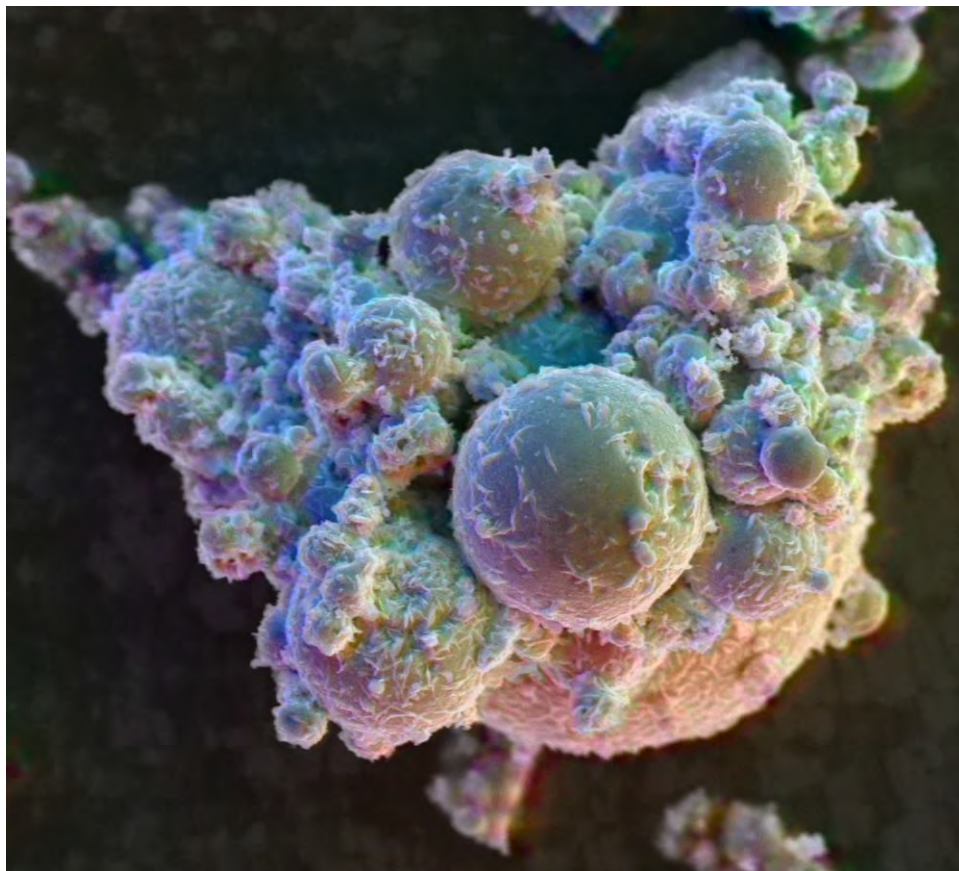

What we're left with is this rather colorful, if not slightly unusual, image with a lot of topography specific color.  
(Yes, I know the colors in this image are a bit obnoxious. It's just an example, get over it. Check out the next chapter to better understand how to get nicer images from this technique.)



## \*Tips Tricks and Suggestions For Multi-Detector Color

One of the great and horrible things about this technique is that there are a million gazillion different options of things we could do to make this final color image. It's great because we can tweak the image in a lot of different ways, but it's difficult because there are *SO* many options.


### **A. Opacity**

Here are some examples of how, by changing the opacity on the B/W layer in **Step 22**, we can really control how subtle or “in your face” the colors are.

30% Opacity will give you some SEVERE colors.

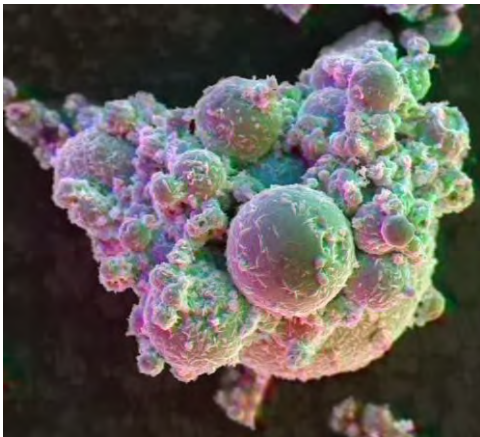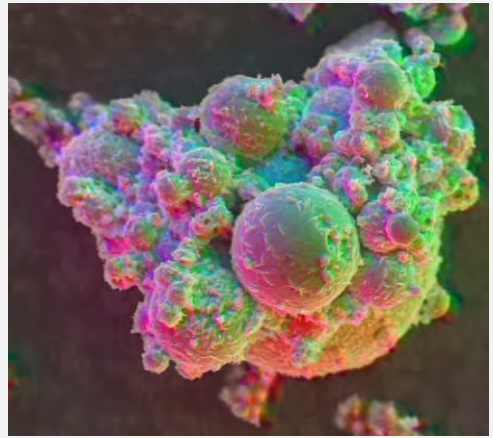

50% is usually a happy medium.

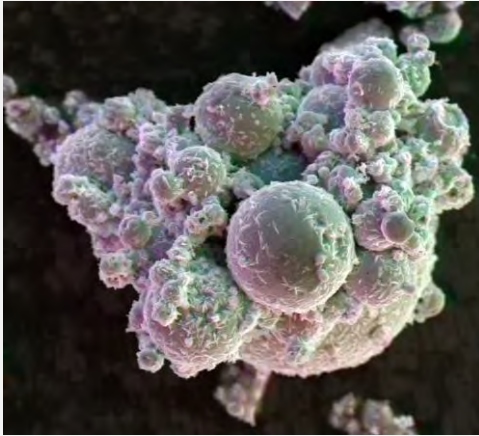

70% give us more subtle colors.


This was a heavy metal sample and the BSE signal coming off of it was pretty severe. As a result, the colors came out **severely** severe. Not every sample will respond exactly the same way to every technique. You'll need to turn down the gain on your BSE detector (if possible), or you'll need to but a big wet blanket on these colors. \*See .



## B. Taking BSE Images

One of the things we can consider is HOW we take our images. As we talked about before, I used different segments from a 4 segment BSE detector and the results were great. BUT, there are a number of different ways we could go about doing that.

First off, if we take a look at the setup for the detector I used, it actually has 5 segments. There is one extra segment (labeled E) that is off-axis and will give a little shadow to the BSE image (which, as you may know, gives a rather flat-looking image, since the detector is directly above the sample).

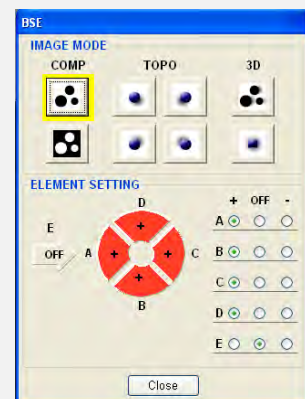


Turning on all 5 detectors is great for taking your master image that won't get colored. (On this system it's labeled as 3D.) So here we can see our sample, and despite being a BSE image, it's got a little bit of depth and definition to it, thanks to the extra segment.



Now, how should we take our other, more directional BSE images? We have a few options.

The first option, and most obvious thing to do, is to turn on one segment at a time and take a picture.

How many images you take is up to you (if your detector only has 2 segments, I guess you're going to take 2 images). I recommend taking an image from all directions it's possible to with your detector. You don't need to use them all, but it's good to have them to pick and choose from later, just in case.

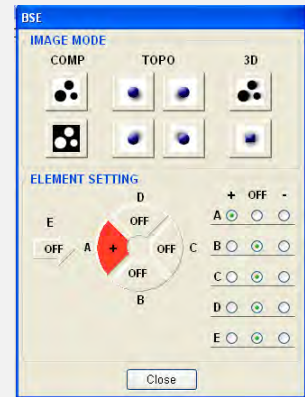


Each image will be fairly directional, so one side of the image is bright, while the other side will have a big shadow.



Our second major option (if your SEM has it) is to use the TOPO mode. This turns 2 segments positive (marked in red as +) and the other 2 are turned negative (marked in blue as -). So what does that do?

Well, positive means on, and we know what that is. Signal we get shows up as white, so the more signal we get from a certain area, the lighter colored it is.

If a segment is turned negative, it does the opposite of that. Signal then shows up as black instead of white.

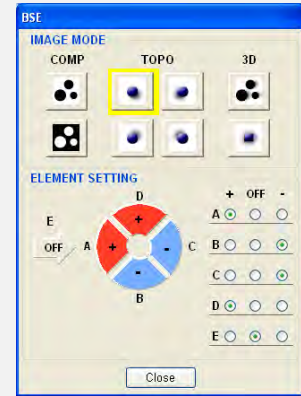

As we can see here, it creates a kind of weird, serious shadow effect.



Let's see how the 2 compare side-by-side.




The images above were each taken with the single segment POSITIVE and the other 3 OFF. The lower 3 images were taken in TOPO mode, with 2 segments POSITIVE and the other 2 NEGATIVE.



So how do these 2 different techniques affect the final result? If we work on our samples exactly the way we described in the instructions through Step 17, this is what we have:


First up we have our color image made from the single-segment images. The colors are a little bright and a little dramatic, but not overly so. The other thing we need to note about this image is that it **does not** have the B/W master image on top, like we had discussed doing in the instructions.

If you like the brighter colors, then this is great! Otherwise we'll move to the next option.

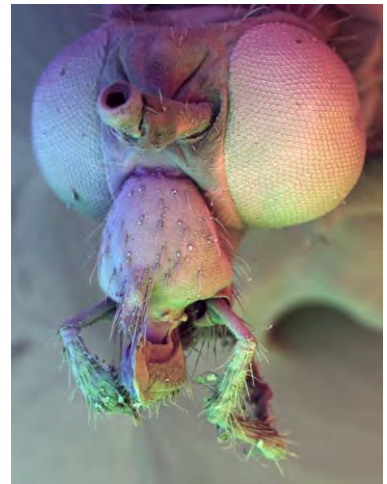


OK, here's the same image, but with the B/W master image on top at 50% opacity.

The coloring still seems dramatic, but the colors are much more muted. Again, this is something you can tinker with.

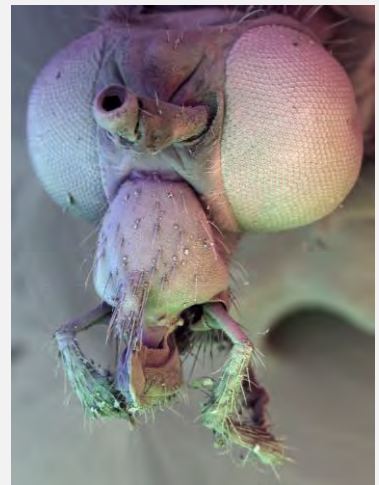


Our other option was the TOPO mode images. The differently colored areas aren't quite as sharp or serious as the image we got from combining the single segment images. (Using the TOPO mode is my preferred method, mainly because I'm using all 4 segments at the same time, rather than just 1. This should result in an image with better signal and less noise. But like I said, it's just a preference and this is totally up to you.)

Another consideration to take into account is that sometimes not every sample responds the exact same way to a particular technique. So trying things different ways *might* be beneficial.

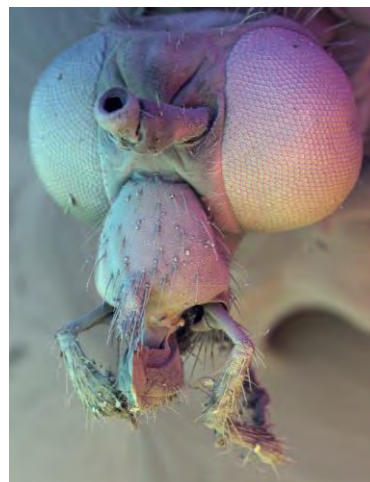

One other thing we can consider is **when** should we perform the **Levels** operation? In the instructions I said to wait until everything is finished to do them. Here are the single-segment images colorized, but I did the levels on each of the B/W images BEFORE I combined it.

The colors here came out pretty serious. (Again, this still does not have the master B/W image on top.) And since the levels of this image were pretty solid to start off with, there's not much we can do with them now (aside from maybe adjusting the color balance [**Shift + Ctrl + B**]).

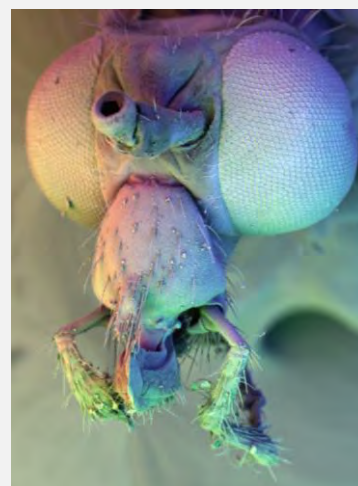

You may notice that still, given all the many little things we can do from start to finish in this process, you still will end up with a generally similar image. They have their differences, but all the examples we just looked at are not *WILDLY* different. Now I'm not saying that small details are not important (surely they can be), but as with the other techniques we've discussed, more time and effort = more better results. (Or more gooder results, whatever.)

### **C. Using ONLY 1 SE Detector**

OK, maybe you have a cheap SEM that only has 1 SE detector and no BSE detector. **OR** you took your images using the BSE detector and the image came out looking insane. How do you get images from different directions like this only using a single SE detector?

**WE ROTATE THE SAMPLE!**




If your microscope isn't *SUPER* old or *SUPER* cheap, you will probably have some kind of eucentric rotation option somewhere in the software to rotate the stage around the area you're currently looking at.

There are a number of things to keep track of when attempting this.

**First**, this won't work if you're using an "upper" (through the lens) SE detector. It's always looking straight down on your sample, so rotating your sample around isn't going to do anything. This will only work with a "lower" (in-chamber) detector.

**Second**, find an area on your sample that you think will be visible as you rotate your sample around 360°. That might not be as easy as it sounds. Don't be surprised if after you've taken 1 or 2 images, you have to start over in a new area because you lost what you were looking at. **TIP:** What I might do is center your area of interest and then take an image of it at a much lower magnification. This should make the general area easier to find after rotating the stage.

**Third**, you need a way to center your area of interest on the screen as it's rotated around. Most SEM's will have something like an area marker or crosshairs you can turn on that will help align your sample in the middle of the screen. (If your SEM doesn't have something like that, the corner of a sticky note or a piece of tape stuck on the screen will work just as well.) In the example we're going to look at I'm going to keep the middle of this pointy thing in the middle of the screen.

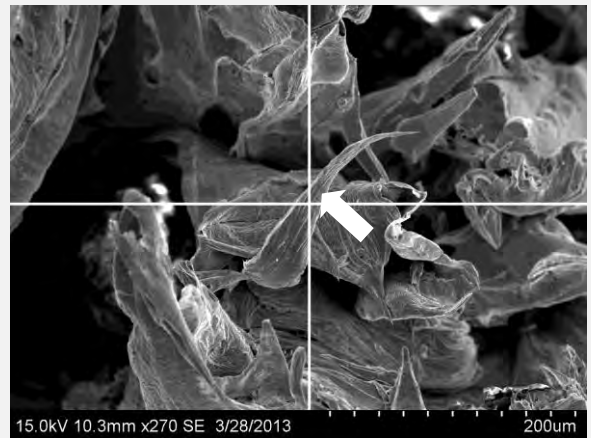

**Fourth**, we need to decide how many images we're going to take, and therefore, how much to rotate the stage by in between images. 4 images would probably be nice, so we can have 2 or 3 color images and 1 master image, like we have in the BSE images. (If you're using an SEM that has a through the lens detector, take your master image with that.) So we would just rotate the stage by 90° each time after taking our images. In the example I've cooked up, I took 3 images, so I rotated the stage by 120°.

At this juncture, we must ask ourselves two very important questions about our SEM.

1. Does my SEM have a software controlled stage?
2. Does my SEM have a Raster Rotation in the control software (and not a wobbly knob) that I can control the exact number of degrees my image is rotated?

If the answer to both questions is **YES**, then you're in luck and this will be pretty easy if you just follow the instructions under **Option A**.

If you answered **NO** to either of those questions, you're going to need to follow the steps under **Option B** for this to work properly.

### OPTION A: Software Stage Rotation

Easy.


1. Center your structure of interest, do all alignments and take an image.
2. Enter the number of degrees you want to rotate the stage around into wherever the stage rotation stuff is located in your software. (This would probably be best if it was set to Relative rotation and not Absolute.) If your stage has a eucentric rotation function (and it works the way it should), make sure it's turned on.
3. Once the stage is rotated, re-center your feature.
4. Turn on the **Raster Rotation** on the SEM and type in the same amount you rotated the stage (i.e. 120°). Your image should then be oriented in the same way your first image was. (If it's not, you probably rotated something in the wrong direction, like positive instead of negative.) Take a picture.
5. Repeat steps 2-4 again until you have all your images taken. (3 for 120°, 4 for 90°, etc.) Make sure the images are all at the same magnification, or you'll be in trouble later.
6. At that point you should have all you need to combine your images as per the regular instructions. But you may need to align the images a little bit.

### OPTION B: Manual Stage Rotation

This should work pretty good as well, but since it won't be as precise as typing a number into software, we have a few extra hoops to jump through.


1. Once you center your structure of interest, do all the necessary alignments. Because of the rotation we're doing, we're going to end up cropping out a lot of the image later, this means we'll need to take the image at a MUCH lower mag than what you think you need to. A way you can compensate for losing half your image later to cropping is by taking a higher resolution images right now. **For example:** if you normally take images at a resolution of 1280x960 pixels, just bump it up to 2560x1920.
2. Take your first image, rotate the stage around *whatever amount*, re-center your structure and take another image. Make sure all your images are taken at the same magnification, or you'll be in trouble later.

You should end up with a group of images that look kind of like this:

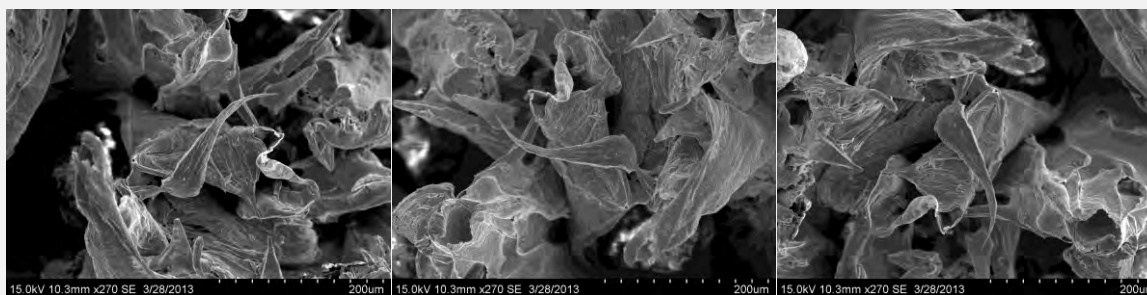

To make the color image happen, we can follow the instructions as before up through **Step 12** of **Section VIII** and then continue with the following instructions:

**12A.** In our example, we've got all 3 images pasted into a new file as separate layers. We have our correctly oriented image on bottom and only 1 of the 2 other layers is visible and set to **50% Opacity**.

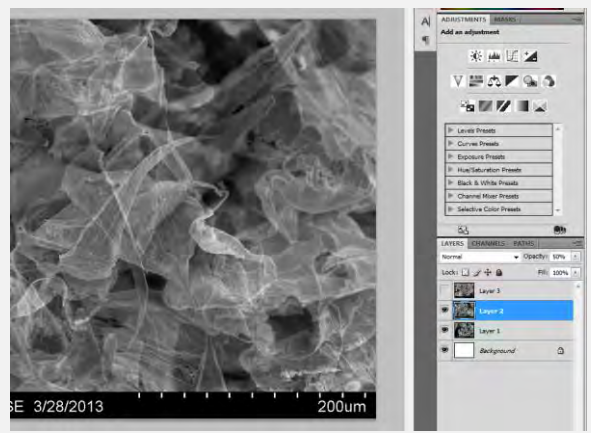


**12B.** With that second image selected, we'll need to rotate it now by going to:

**Edit > Transform > Rotate**


You can also use the **Free Transform** tool as well.

**Ctrl + T**

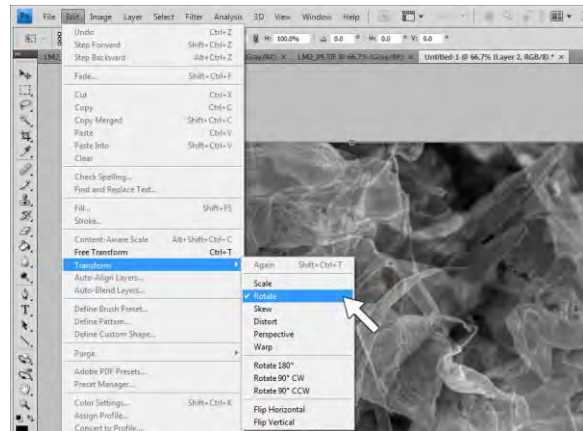

**12C.** If you move the mouse around the edges of your image, your cursor should look something like this ↻

By holding down the left mouse button and dragging the mouse around you can rotate the image. You probably won't be able to get it to line up perfectly at first. Just start off by getting it close.

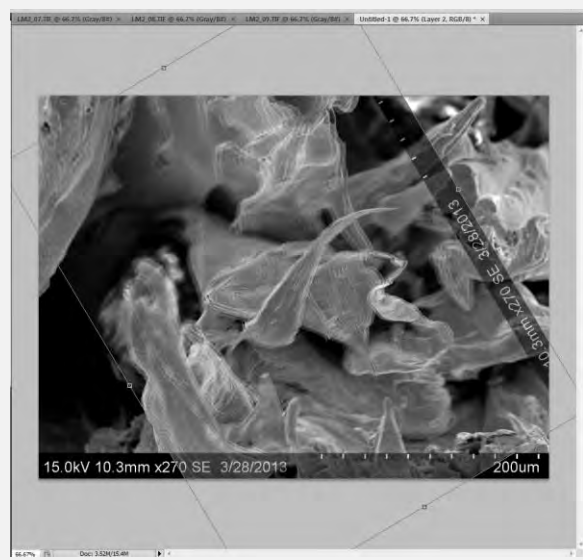

**12D.** If you release the mouse button and move the cursor towards the center of the rotated image, it will turn into the **Move Tool**, which means you can click the image and drag it around in the X and Y direction. Do this to get the 2 images better lined up, and if you need to, go back to rotating the image again.


Go back and forth between rotating and nudging the image up/down, left/right, until you get it aligned.

At some point you might need to switch from rotating it around with the mouse and change the angle in the **Rotation** box by typing in a number. You might have to rotate it in very tiny amounts, so feel free to make the fraction of an angle as small as you need to.



**12E.** Repeat **12A** thru **12C** on any other layers you have. In my example I have 3 images, and it should look something like this.

There's only a little part in the center of the image that's useable, so we'll need to **Crop** it out. We can do that **now**, or **after** we colorize the image. It doesn't seem to matter.

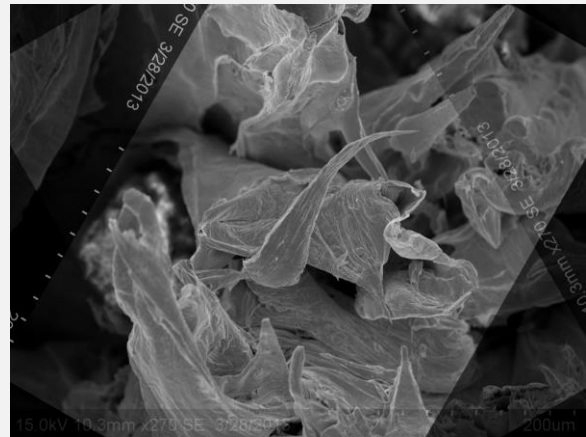

We can now continue with the standard **Blending** process as described in **Section VIII** starting at **Step 17**.

Here is what we ended up with. Nifty!

**However**, this sample, which has a lot topography action going on, makes a good example of why we might want to try **LESS** color, instead of **MORE**.



## D. Less Can Be More

**EX1.** I redid the color on this image, but I only colorized 2 of the layers instead of all 3, and we end up with something that looks slightly more normal.

It came out kind of green, and of course just by changing around which **Channels** boxes were checked, I could have made it come out purple or blue or red or whatever.

I think all those subtle gradations of color you can get with this technique are really fantastic though.



**EX2.** I liked that green image, but from that green image all I needed to do was auto-adjust the color:


**Image > Auto Color**

or **Shift + Ctrl + B**

And I end up with something that looks like a more sedate version my original image with 3 color layers.




**EX3.** Unfortunately, this technique doesn't make *everything* look cool. I had another sample that I processed. It didn't come out like I thought it would. In fact, I hated it. I messed around using 3 and 4 color layers a hundred different ways and I never got anything I liked out of it. It was just way too busy. There were *too many* colors.

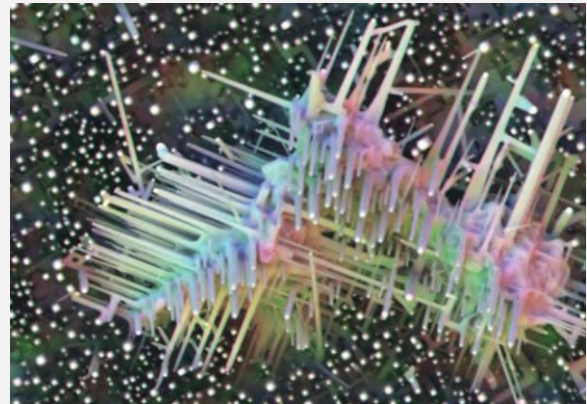

Fortunately, before I deleted the images out of disgust, I decided to use only 2 color layers instead of 3. Just by doing that one little thing, its eyesore status was revoked.

Another thing that was brought out by reducing the number of color layers was you really got a better idea about what the topography of the sample is like. With the 3 color layers, it was just a big jumble of color.


(This image also has the master B/W image on the bottom.)

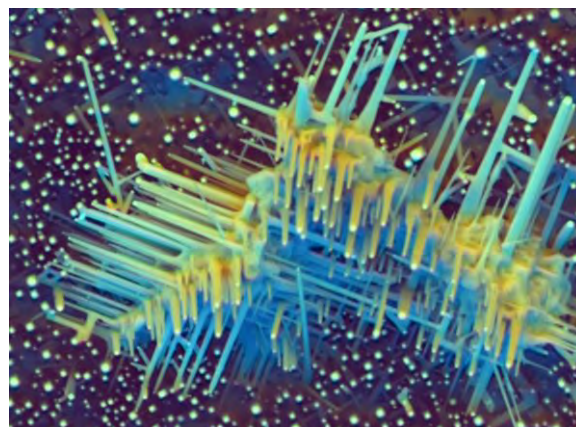

**NOTE:** When you're choosing from your different B/W images, does it matter which ones you pick when only doing 2 colors? Sort of.

What I think is the best option is to get images that were taken from opposite sides of the sample. So for example, I would take image 1 and 3 together, but not 1 and 2, or 2 and 3. (Now if your BSE detector only has 2 segments, then that's what you'll end up doing anyway.)


By doing this, both sides of the sample will be illuminated by a different color. (Preferably a lighter one and a darker one.)

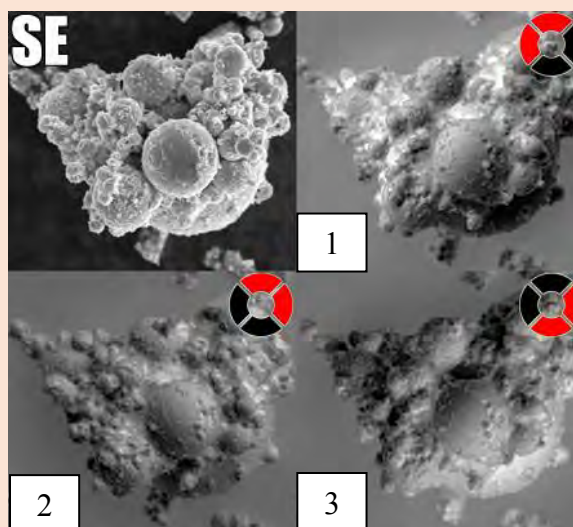

### **E. Change Up the Colors!**


The last thing we'll look at is what we get when we change the color channels around in the different layers. This will change the overall look of the image a lot more than most of these other changes we've been talking about.

This first image came out nice, with some subtle, muted colors. I really like the way these colors show off the surface structure of the sample.

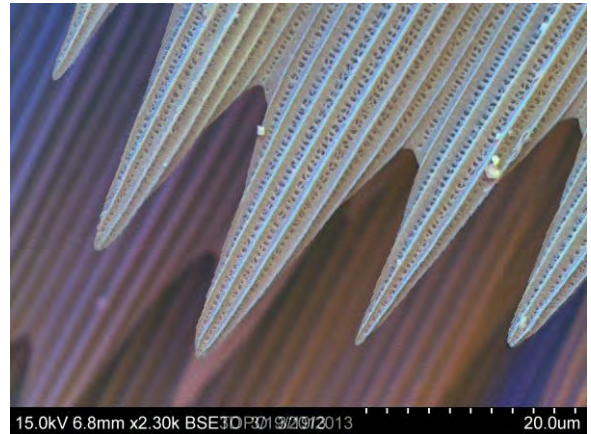

The next image looks **totally** different. We have vibrant, fun colors that are in broader range than the previous image. What did I do to make them different? Simply checking and unchecking different **Channels** boxes in **Steps 19** and **20**.

That's it.

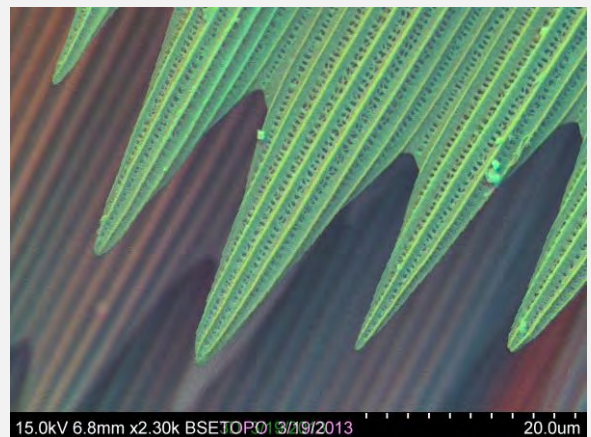

In conclusion, the fantastic thing about this technique is not only the ease and speed you can put these images together, but also the fact that all the tiny details in your image are not lost, like they can be when using something like the **Gradient Map**.

**I LOVE THIS TECHNIQUE!**

## X. Finishing Moves For Your Color Images

The main problem we're going to have with all of this fake color we're putting on our images is that it's fake. As a result of being fake, it often looks fake. Here are a few tips to de-fake-ify some of your colors.


### A. Saturation (Great for ALL fake color images)

**1.** *This is the best option* to make your colors believable.

With your final image opened, open the Hue/Saturation window.


**Image > Adjustments > Hue/Saturation**

**Shift + U**

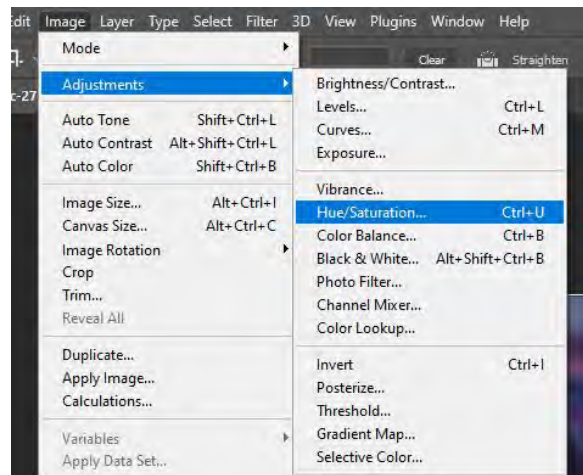


**2.** By dragging the **Saturation** arrow to the left, you can leech the life out of the image – but in a good way.

You'll want to adjust the saturation **down** until the colors are still present, but such that they look more realistic.

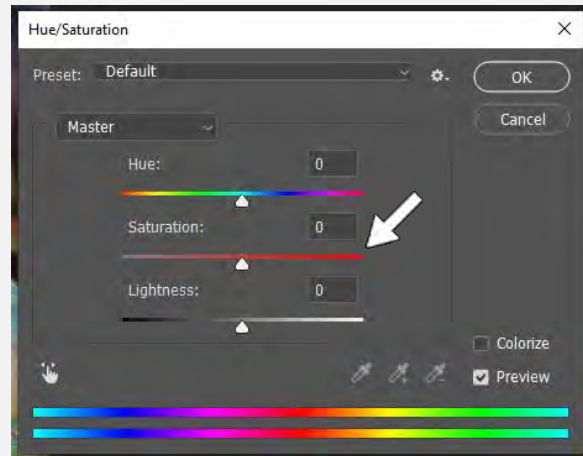




**EXAMPLE A1 Before:** Here is the original image that we colorized earlier in the book.

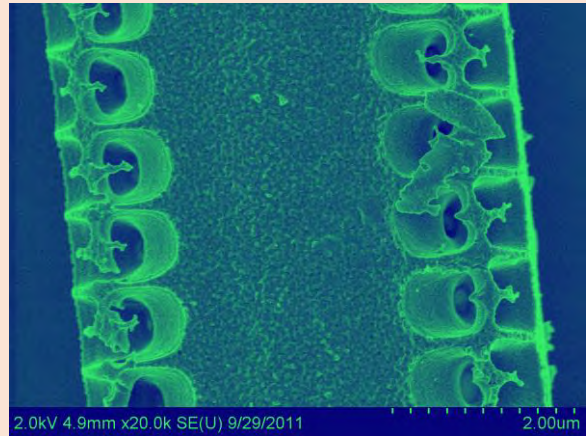


**Example A1 After:** In this image, the **Saturation** slider was moved to -80. The colors are now very dull.

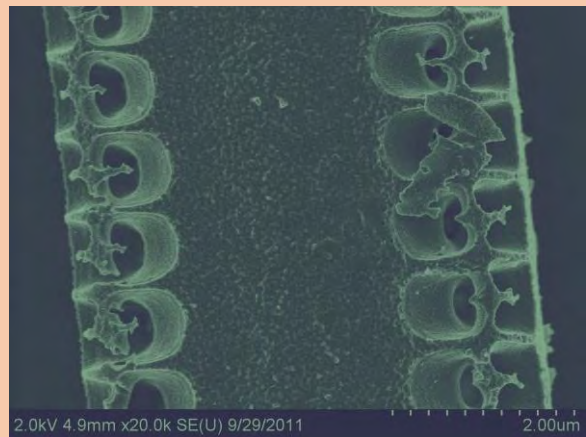

**Example A1 Final:** In this final image, we re-did the **Levels**. These colors now look much better and less offensive than they did in the original at the top of the page.

Adjusting the **Saturation** down a notch can really take your image up a notch!

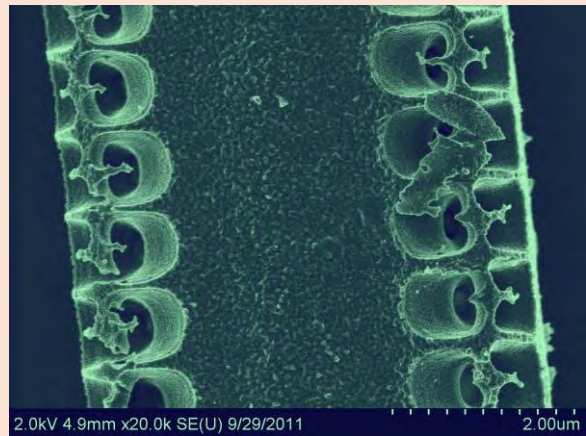

**Example A2:** The image on the left is the original “crazy color” image. The image on the right has had the **Saturation** turned down to -65 and the **Levels** we redone. We can see how this makes the colors look less cartoonish and a bit more believable.




## B. Hue (Great for Multi-Detector Color Images)


**1.** Let's say you're not quite happy with the exact colors you have in your image and want to shift all of them around – either by a little or a lot. Adjusting the Hue can help. With your final image opened, open the Hue/Saturation window.

**Image > Adjustments > Hue/Saturation**

**Shift + U**

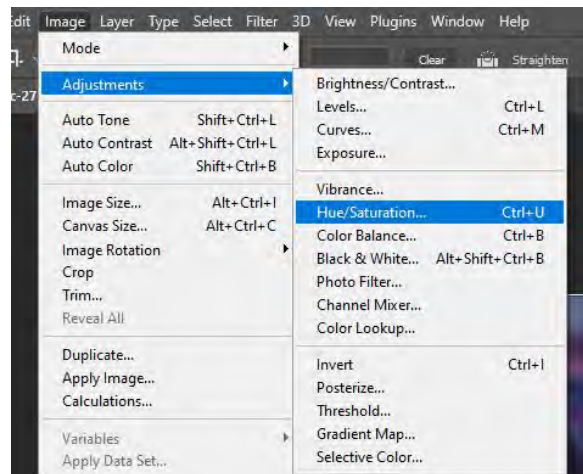

**2.** By dragging the **Hue** arrow back and forth, you can shift all the colors in the image around the color wheel. (If that doesn't make sense, you'll just need to try it and see what happens.)

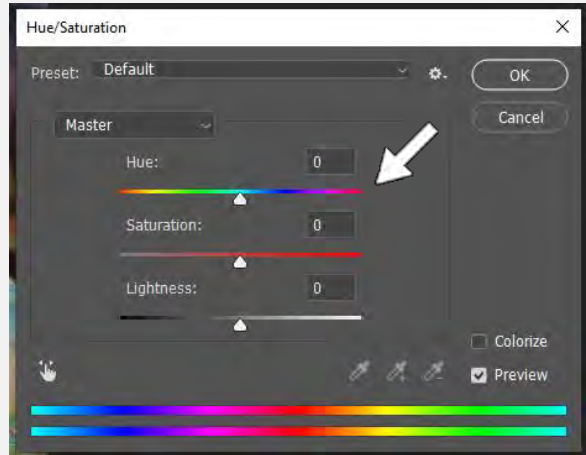

**Example B1 Before:** You'll recall this crazy image from earlier. Let's assume we don't like the amount of pink in this image, and would prefer it to be something different.




**Example B1 After:** Here we have a totally different set of colors in this image, and all we had to do was drag the **Hue** slider around. We didn't have to go back and readjust all the color layers, etc. This can be really helpful with the **Multi-Detector Color** technique, since you can adjust the colors quickly and easily.




### C. B&W On Top! (Also Great for Multi-Detector Color Images)


**1.** Another good way to tone down your colors in an image is to take the original B/W image and put it as a separate layer on top of all your other images. You can then adjust the colors by adjusting the amount of opacity on that top B/W layer.

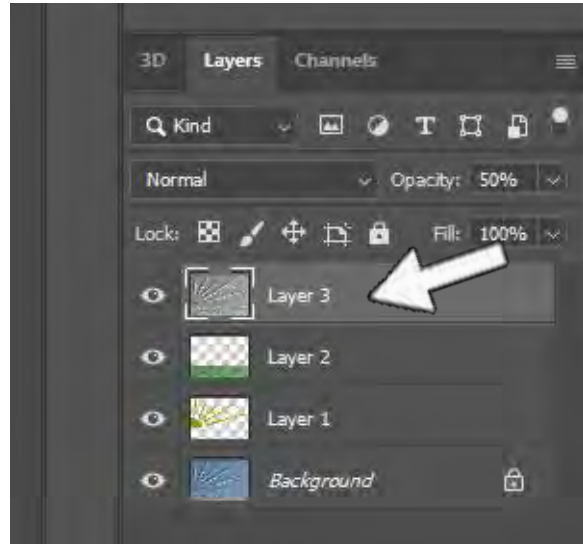


**Example C1 Before:** This is the original image we colorized earlier in the book. (It's intended to look simple and a little silly. Calm down.)

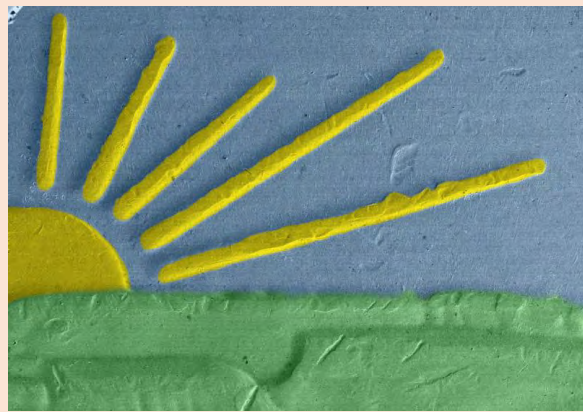

**Example C1 After:** This image had the B/W image on top at 50% opacity, and then was flattened, and then Levels were performed. These colors (though still a little silly) have a much more tangible and realistic look to them.

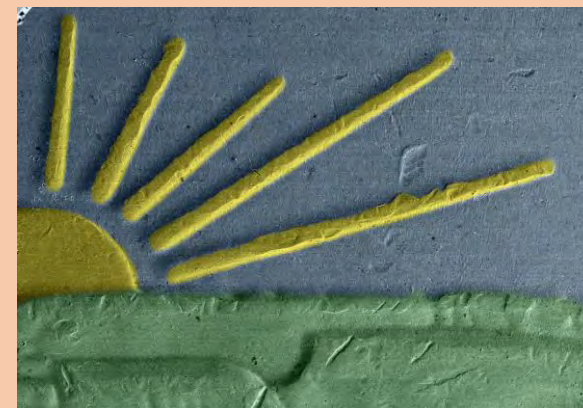

**Example C2 Before:** Here's an image I put together with 3 color layers. It's ok. But maybe I don't want to tone down the colors a bit. (Crazy colors are OK, but best used in moderation.)

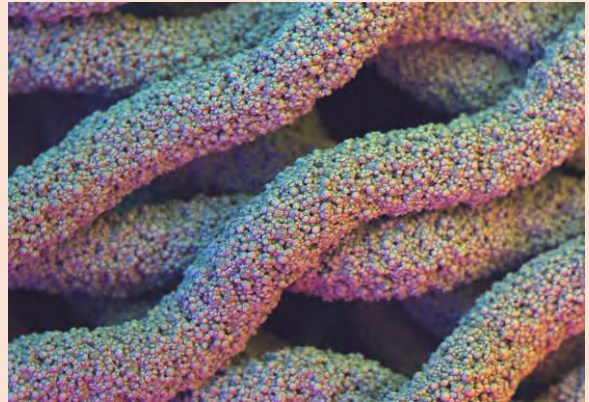

**Example C2 After:** I put the master B/W image on top of all the other layers at 50% opacity. I then flattened the image and did **Auto Color** correction (**Ctrl + Shift + B**), then **Levels** (**Ctrl + L**).




Now, compared to the more colorful image I made before, it looks pretty dull. BUT, if you just look at it, on it's own, I think it's really fantastic. Sometimes the colors come out so subtle, it's almost like you're looking at a black-and-white image at first. It may even take you a second to realize the colors are there. I think it's that aspect of discovery that makes this image exciting.

This is just another option to try in your bag of tricks.

#### D. Just Don't (Ever)

Don't put lens flares in your microscope image.

We all know there are no glass lenses in the SEM, nor can you see the sun from inside the microscope to make such a flare. This is just not good... at all.

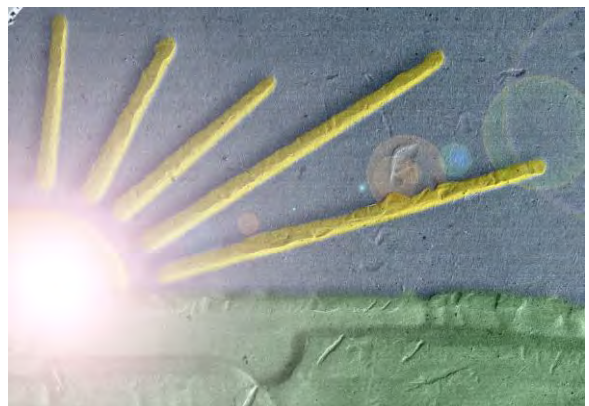

## XI. Drawing a Mag Line Marker (The Cheating Way with Pre-Existing Mag Marker)

So you want to draw a Mag Marker on your image. We have 2 options.


**Option A:** if you already have a mag marker on your image put there by the SEM/TEM software and you want to draw a different one.

**Option B:** if there is **no** mag marker and you know the magnification. Go to **Section IX**.

**WARNING:** This section will involve math!

**1.** Open your image in Photoshop and select the **Rectangle tool**.


Rectangles are nice because they always come out perfectly horizontal AND you can change how thick you want it to be as you're drawing it.

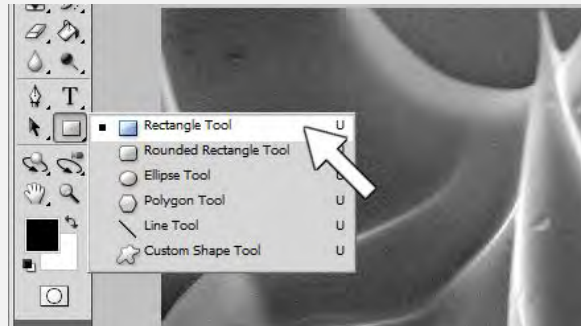


**NOTE:** You can draw a straight line using the **Line Tool** if you hold down the **Shift** key while you're drawing it. But you might not get the thickness right the first time. You can change that in the **Weight** box. But you'll need to be careful to make sure you also type in the correct unit, **px** for pixels, **in** for inches, **cm** for centimeters, etc. This will be the **WIDTH** if you line, you'll draw it as long as it needs to be. A 10 inch wide line is going to be quite a bit different than a 10 pixel wide line (and PS seems to default to inches when you type in a number and hit enter here). So you will need to pay attention to this if you're using the line tool.



**2.** The shapes toolbar should appear at the top of the screen. Check to make sure **Create new shape layer** is selected. Also, we want to check what color our rectangle will be. By looking at the box next to where it says **Color** we can see it's set to black. If we click on this box the **Color Picker** will open and we can choose our color. White is probably the best.

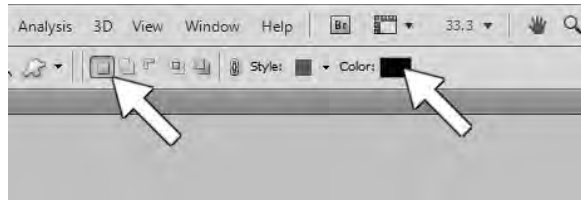


**3.** We're ready to draw our mag marker now. Since this micrograph already has one, we're just going to borrow it by drawing our new one right on top of it. That way, we know it's the right size.

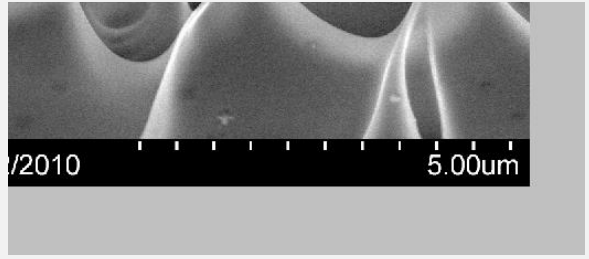

**4.** Here is the problem though: where do we start and stop drawing our rectangle in relation to the existing line? As illustrated here, we don't want to draw our new marker from the inside or outside edges of the existing marker. It needs to be from the same side of the first and last hash marks.

I know just reading that explanation might not make sense, but you should be able to figure it out eventually if you stare at the illustration long enough.

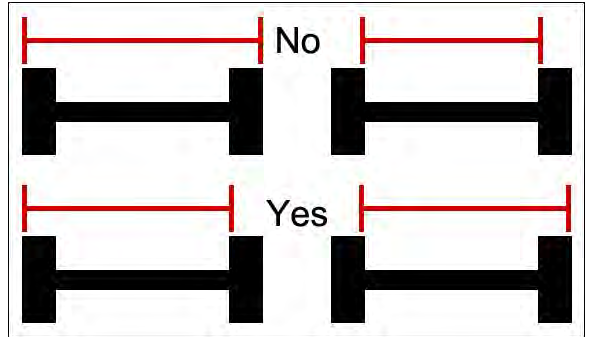

**5.** If the existing marker is too small to be able to draw your rectangle over accurately, zoom in by selecting


**View > Zoom In**

or **Control + +**

You can do this multiple times. To zoom out, select


**View > Zoom Out**

or **Control + -**

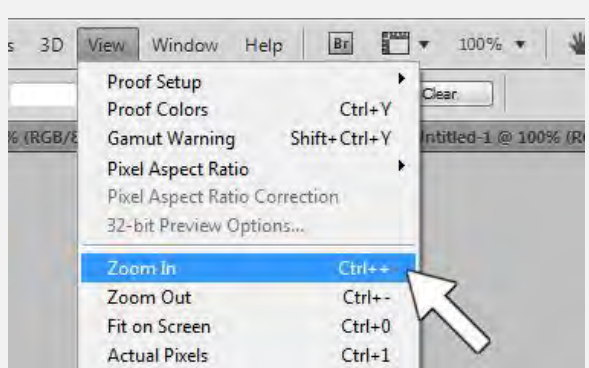

**6.** We can go ahead and draw our marker now by clicking at one end of the existing marker, dragging the cursor to the other end of the marker and releasing. The layer the new marker is in should be selected in the **Layer** tab in the lower right-hand part of the screen.


An exception to this would be if our image is a **.BMP** file which cannot accept new layers. If that is the case, Save your file in a different format, like **.TIF** before drawing the line.

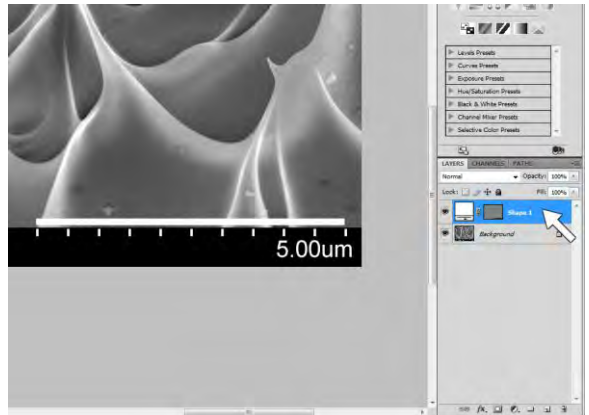

**7.** The existing line marker is 5µm. Let's assume for a moment that we wanted our new marker to be a different length, like 2µm.

With our new line's layer selected in the **Layer** tab, we click on

**Edit > Transform Path > Scale**

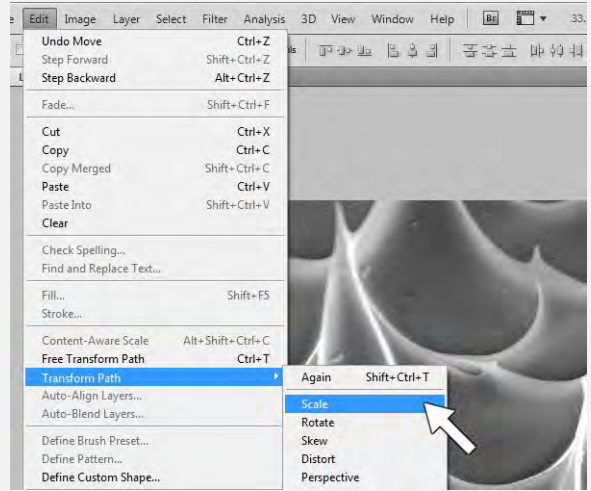

**8.** A toolbar at the top of the screen should appear that gives us a number of options. Width (**W**) and Height (**H**) are currently at 100%. To make our line shorter, we just change the **W** value. But by how much?

Original line = 5µm = 100%


$100\% / 5\mu\text{m} = 20\%$

20% of original line = 1µm.


Therefore a 2µm line should = 40%.

**YAY MATH!**

I enter 40 into the **W** box.

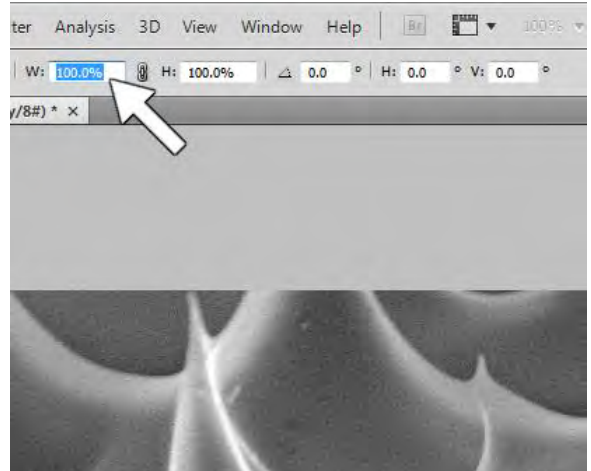


**9.** To commit the change to your line, click the checkmark on the right end of this toolbar.

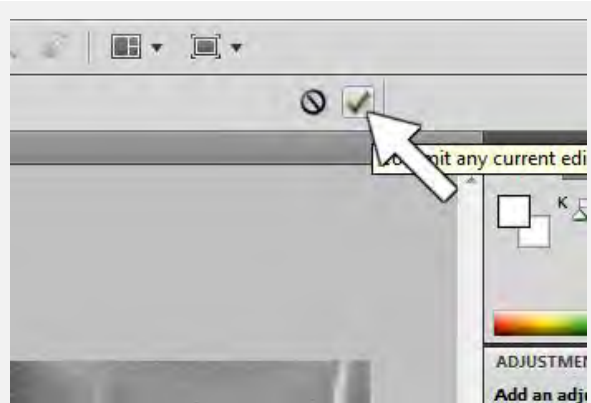





**10.** Select the **Move Tool** (V on the keyboard).  
Using the **Move Tool**, you can drag your new marker around the image and position it where you want.

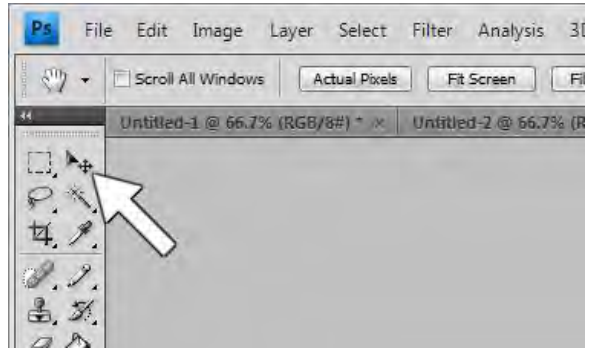

**11.** We've moved our marker and now we need to label it. Select the **Horizontal Type Tool**.

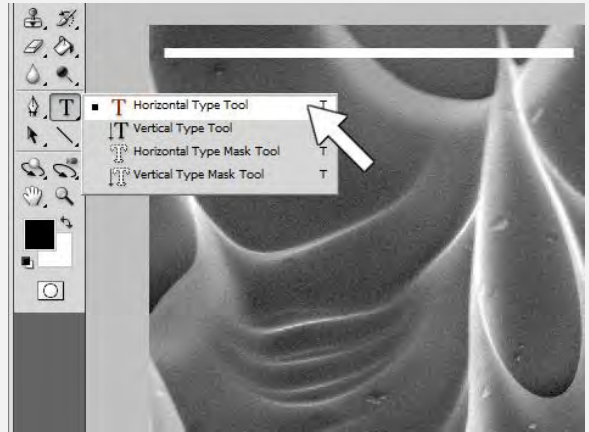

**12.** A toolbar along the top of the screen should appear where we can choose all of our text options. I can see right now I need to change the text color from black to white. The size also says 30pt, I'm not sure right now if that will be too big or too small. We can change it later if it's no good.

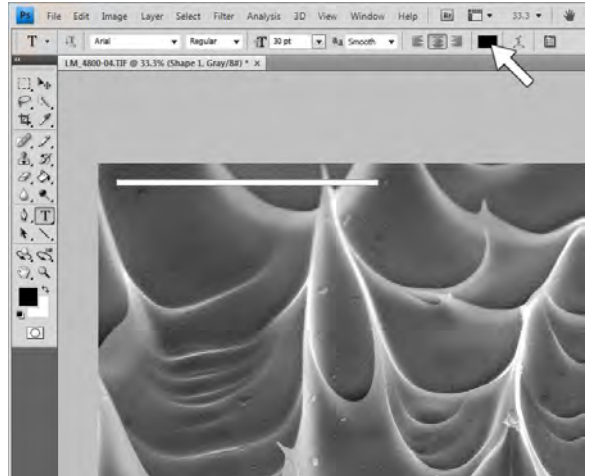


**13.** When you're done choosing all your options we should be able to click on our image generally where we want our text to go and start typing.


I can see that 30pt was way too big. So I select my text by clicking and dragging the cursor over it and then selecting a smaller size with the font size drop down menu in the toolbox.

When you're happy with how it looks, click the checkmark on the right end of this toolbar to commit the change to your line.

Now you can use the **Move** tool to position the text where you want it.

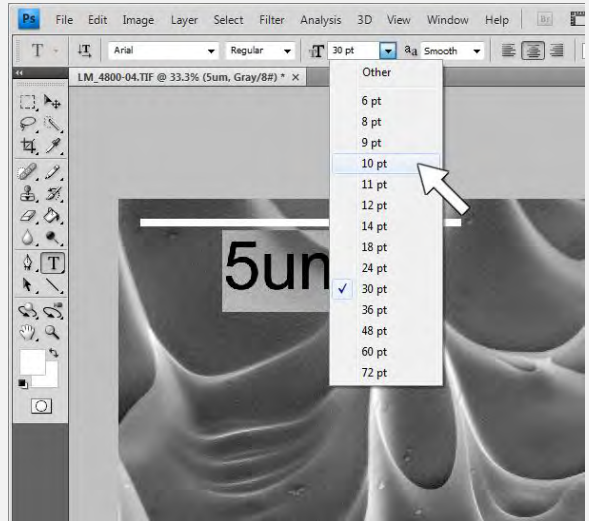


**14.** If, for whatever reason, you want to edit your text after you've hit the **Commit** checkmark, you'll need to select the layer your text is in in the **Layers** tab in the lower right corner of the screen and double-click the white box with the **T** in it. Your text should now be editable.

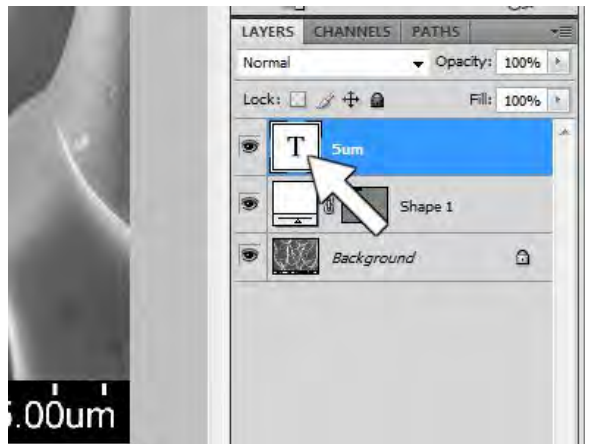

**15.** The **only** reason to put a new mag marker on your image is because you're going to crop the old one out. Select the **Crop** tool and drag the cursor across your image to select the area you want to keep. You can adjust the selected area by grabbing and dragging the small gray boxes that appear on the selection box. When you're happy with your box, click the **Commit** checkmark in the upper right of the screen.

Your image should be ready for saving.

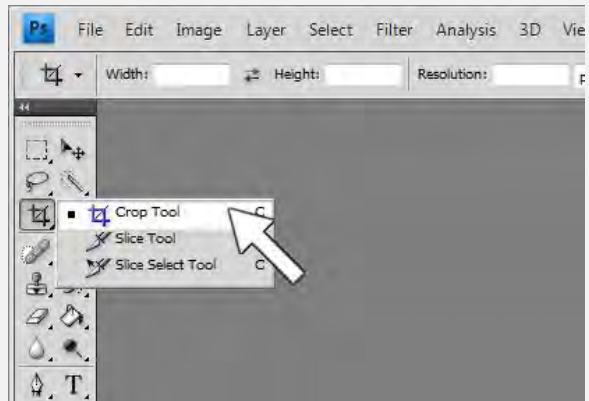


## XII. Drawing a Mag Line Marker (The SUPER Cheating Way with NO Existing Mag Marker)

**CAUTION:** This technique **ONLY** works with Adobe Photoshop. Don't try it with something else like ImageJ! It totally won't work.

**EXTRA CAUTION:** This only works if your image **HAS NOT** been resized. (Unless you know by how much. i.e. If it's exactly half the size, reduce your line 50%, etc.)

**SUPER CAUTION:** This will only work with .TIF, .JPG, .BMP and .PNG files. It **WILL NOT WORK** with a .GIF file. You have been warned.

**1.** Let's say we have an image that has **NO** mag line marker on it. What do I do?? **First**, we ask ourselves if we know the magnification. If the answer is yes, then you can proceed. If not, you may commence crying. For this example, I'm going to use an image that **DOES** have a mag marker on it. This is to prove that this technique actually works. (Because in theory it shouldn't work at all... but somehow it does.)

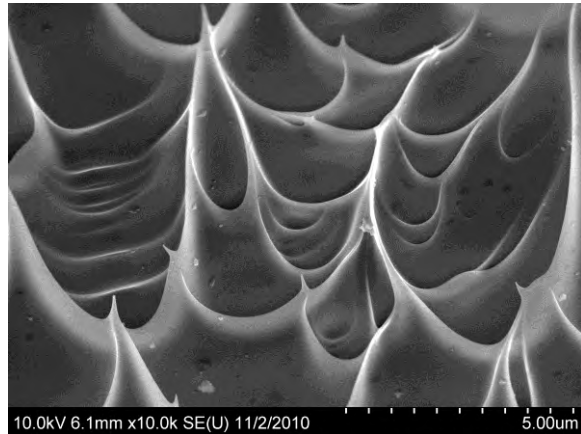

**2.** We take the known magnification of our image (in my example, it's 10,000x.) and we divide it by 10,000. The answer we get is the distance a line needs to be in cm on our image to represent 1µm. That might be confusing. So,  $10,000x / 10,000 = 1$ . This means a 1µm marker on our image will need to be 1cm long. For this example, I need a 5µm line, so I just need to make it 5cm long. Simple!

Magnification	=	The distance of a 1µm line in cm
Ex. $\frac{10,000x}{10,000}$	=	1µm per 1cm of line
Ex. $\frac{250,000x}{10,000}$	=	1µm per 25cm of line
Ex. $\frac{3,200x}{10,000}$	=	1µm per 0.32cm of line


**3.** So how do we draw a perfect 1cm long rectangle on our image?

First we open our image in Photoshop and turn on the **Rulers** (if they're not already) by clicking:


**View > Rulers**

or **Ctrl + R**

This puts a ruler along the top and left edge of your image.

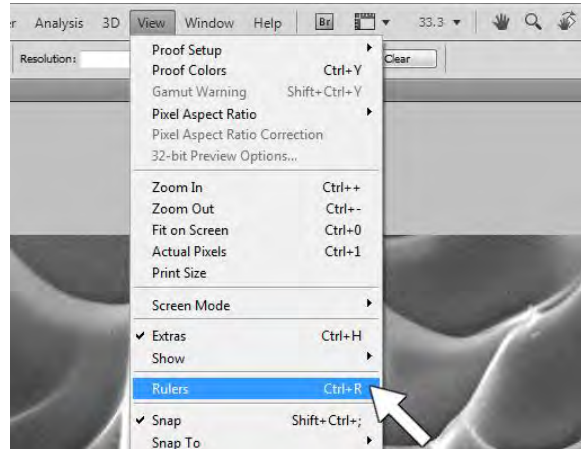

**4.** We then set the units on the ruler to Centimeters by right-clicking on the ruler and selecting **Centimeters**.

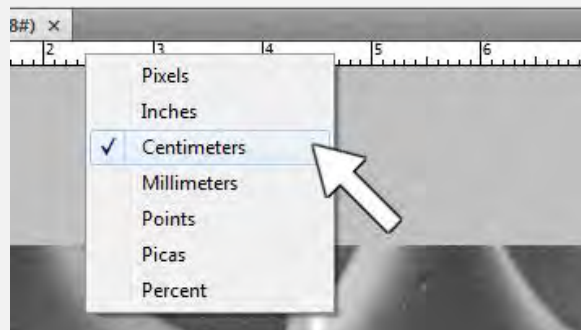

**5.** In order to draw our rectangle a very precise length we'll need to open the **Info** box by clicking

**Window > Info**


Or **F8**

The box will show the width and height of our rectangle in cm as we draw it.




**6.** Select the **Rectangle tool**.


Rectangles are nice because they always come out perfectly horizontal AND you can change how thick you want it to be as you're drawing it.

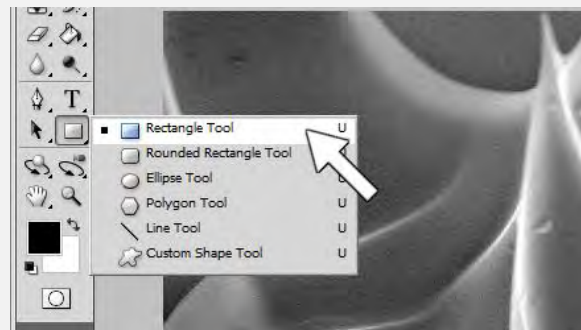


**NOTE:** You can draw a straight line using the **Line Tool** if you hold down the **Shift** key while you're drawing it. But you might not get the thickness right the first time. You can change that in the **Weight** box. But you'll need to be careful to make sure you also type in the correct unit, **px** for pixels, **in** for inches, **cm** for centimeters, etc. This will be the **WIDTH** of your line and you'll draw it as long as it needs to be. A 10 inch wide line is going to be quite a bit different than a 10 pixel wide line (and PS seems to default to inches when you type in a number and hit enter here). So you will need to pay attention to this if you're using the line tool.

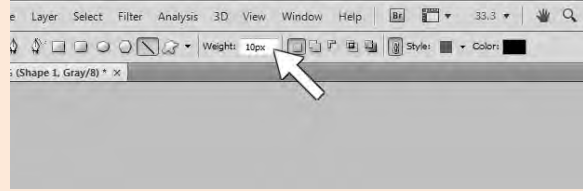


**7.** The shapes toolbar should appear at the top of the screen. Check to make sure **Create new shape layer** is selected. Also, we want to check what color our rectangle will be. By looking at the box next to where it says **Color** we can see it's set to black. If we click on this box the **Color Picker** will open and we can choose our color. White is probably the best.

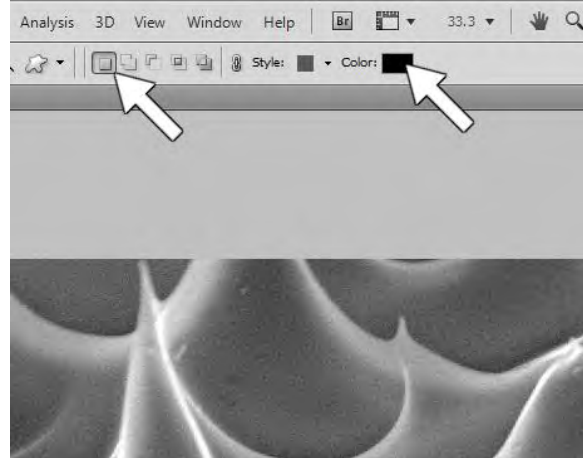

**8.** Now we can draw our marker. (It doesn't matter where we put it since we can move it around later.) Keep an eye on the info box and draw your rectangle out to the correct length. In this example, we want X to equal 5 (circled in red).

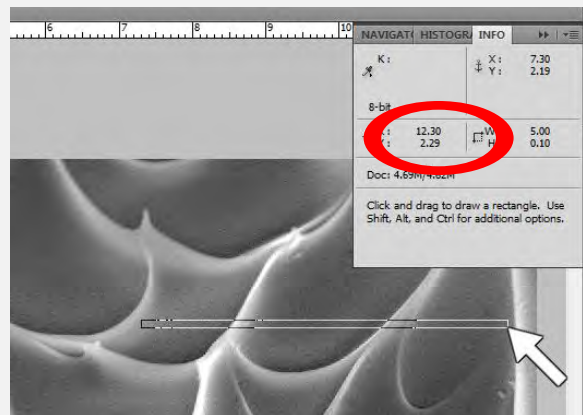


**9.** Our newly drawn layer should be selected in the **Layers** tab in the lower right corner of the screen and automatically named something exciting like "Shape 1".

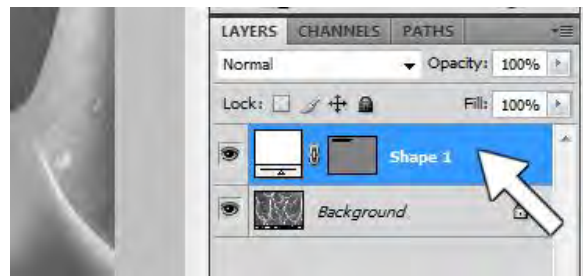


**10.** Select the **Move Tool**, or **V** on the keyboard.

Using the **Move Tool**, you can drag your new marker around the image and position it where you want.

Then you can put a label on it if you want. The text label process is discussed earlier in this section in steps A11 thru A14. Your image should be ready for saving.

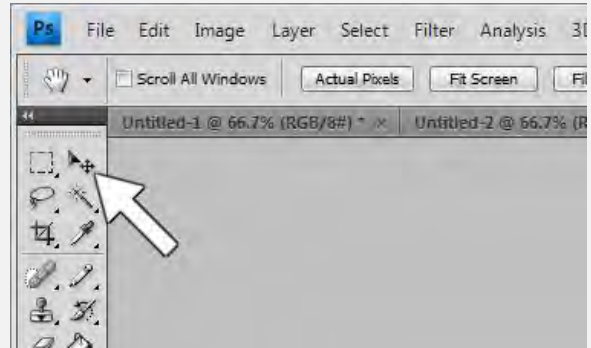

**11.** Lastly, I took the line I drew and placed it here over the existing marker to show that this actually worked. We can see that it came out perfectly.

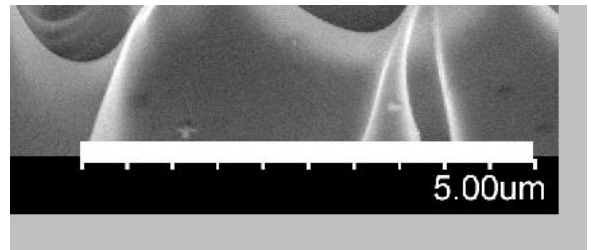

I've tried this with numerous micrographs from several different microscope manufacturers (and even on actual photographic micrographs). It didn't matter what size or resolution they were, as long as it was the **original non-resized image**, not a .GIF, and done in Photoshop, this procedure works. Just don't ask me why.

## **Appendix I** **File Formats**

**PHOTOSHOP (.PSD or .PDD)** This is the native file format for Adobe Photoshop. While you can't really use this file format for display purposes, it's great for saving files that you are currently processing. For example, if I can't finish fixing up my image in one sitting, I might save it as a **.PSD** so that I can save all the layers and changes I've made to my image so I can pick up again right where I left off at a later time. Then when I was finished fixing my image, I would **Save As** one of the actual file formats listed below. This file is a good insurance policy to have, especially if you need to go back and tinker with 1 little thing in your image, you won't have to completely redo it.

**TIFF (.TIF or .TIFF)** Uncompressed file format with a large file size. In one way you can think of the image data being very pure, in that it is not compressed, with the large file size being the only major downside. Good for scaling images up and down. This is my preferred file type when saving images from a microscope. **TIFF** files can also save layers, like a **.PSD** file.

**JPEG (.JPG or .JPEG)** Compressed file type with smaller file sizes. Image compression is good and bad. Good because it makes your file size small, bad because it degrades the quality of your image by averaging pixels together. When you save an image as **.JPG**, Photoshop will give you the option of how much you want to compress the image. Here you can see the file size as a result of the compression, and if you have the preview box checked, you can see any compression artifacts you are introducing into your image. Also, keep in mind that every time you save an open **.JPG** image, it gets re-compressed.



So this is a really bad format to edit an image in, unless you plan on only saving the image **once**. I think it's preferable to edit a **.PSD** or **.TIFF** and then save the image as a **.JPG**, which is much smaller and handier for emailing or whatever.

**BITMAP (.BMP)** Uncompressed file format with large file size. Similar to **.TIFF** in that it is not compressed, and it is also easy to open and work with since it is a Microsoft file format. Unfortunately it does not resize up or down very well.

**GIF (.GIF)** An uncompressed file format with a small file size. An 8-bit indexed color file, which is limited to 256 colors. Good for grayscale images, not so good for color images.





## **Appendix II** **The Short Versions**

### **I. Levels**

	<b><u>Instructions</u></b>	<b><u>Shortcuts</u></b>
1	<b>Open</b> image	<b>Ctrl + O</b>
2	Is image 8bit gray scale?	<b>Image &gt; Mode</b>
3	If no, change.	<b>Image &gt; Mode</b>
4	Open <b>Levels</b>	<b>Ctrl + L</b>
5	Make sure <b>Preview box</b> is checked	-
6	Adjust <b>Levels</b>	-
7	<b>Save</b>	<b>Ctrl + S</b>

## II. Resizing Images

	<u>Instructions</u>	<u>Shortcuts</u>
1	<b>Open</b> image	<b>Ctrl + O</b>
2	Is image 8bit gray scale?	<b>Image &gt; Mode</b>
3	If no, change.	<b>Image &gt; Mode</b>
4	Open <b>Levels</b>	<b>Ctrl + L</b>
5	Make sure <b>Preview box</b> is checked	-
6	Adjust <b>Levels</b>	-
7	<b>Save</b>	<b>Ctrl + S</b>

### **III. Simple Pseudo-Color – Color Table**

	<b><u>Instructions</u></b>	<b><u>Shortcuts</u></b>
1	<b>Open</b> Image	<b>Ctrl + O</b>
2	<b>Duplicate</b> Image	<b>Image &gt; Duplicate</b>
3	Close original image	<b>Ctrl + F4</b>
4	Switch Image to <b>Indexed Color</b>	<b>Image &gt; Mode &gt; Indexed Color</b>
5	Open the <b>Color Table</b>	<b>Image &gt; Mode &gt; Color Table</b>
6	Click and drag cursor from top left box to bottom right	-
7	Select black in the Color Picker. <b>OK</b>	-
8	Select the other color. <b>OK</b>	-
9	Preview Image. <b>OK</b>	-
10	<b>Save</b>	<b>Ctrl + S</b>

**IV. Simple Pseudo-Color – Gradient Map**

	<b><u>Instructions</u></b>	<b><u>Shortcuts</u></b>
1	<b>Open</b> Image	<b>Ctrl + O</b>
2	<b>Duplicate</b> Image	<b>Image &gt; Duplicate</b>
3	Close original image	<b>Ctrl + F4</b>
4	Switch Image to <b>RGB Color</b>	<b>Image &gt; Mode &gt; RGB Color</b>
5	Open the <b>Gradient Map</b>	<b>Image &gt; Adjustments &gt; Gradient Map</b>
6	Click the color gradient bar	-
7	Select a preset to start from	-
8	Select a <b>Color Stop</b>	-
9	Click <b>Color</b> box to open Color Picker	-
10	Choose color. <b>OK</b>	-
11	Repeat Steps 8, 9, 10 until all colors are as desired	-
12	Gradient Editor, <b>OK</b>	-
13	Gradient Map, <b>OK</b>	-
14	<b>Levels</b>	<b>Ctrl + L</b>
15	<b>Save</b>	<b>Ctrl + S</b>

## V. Complex Pseudo Color – Color Table

	<u>Instructions</u>	<u>Shortcuts</u>
1	<b>Open</b> image	<b>Ctrl + O</b>
2	<b>Duplicate</b> image	<b>Image &gt; Duplicate</b>
3	Close Original	<b>Ctrl + F4</b>
4	Adjust <b>Levels</b>	<b>Ctrl + L</b>
5	<b>Duplicate</b> image until you have 1 image for each color you'll be using	<b>Image &gt; Duplicate</b>
6	Set all images to <b>Indexed Color</b> Mode	<b>Image &gt; Mode &gt; Indexed Color</b>
7	Open the <b>Color Table</b>	<b>Image &gt; Mode &gt; Color Table</b>
8	Click and drag cursor from top left box to bottom right	-
9	Select dark color in the Color Picker. <b>OK</b>	-
10	Select light color. <b>OK</b>	-
11	Preview Image. <b>OK</b>	-
12	Repeat 7 thru 11 for each image	-
13	Choose which images will have areas cut out and which image will be the background	-
14	Enter <b>Quick Mask</b> Mode on any image other than the background image	<b>Q</b>
15	Select <b>Pencil</b> tool & brush diameter	-
16	Draw over bulk areas, avoid edges	-
17	Select <b>Brush</b> tool & brush diameter	-
18	Draw over edges	-
19	Exit <b>Quick Mask</b> mode	<b>Q</b>
20	<b>Invert</b> selection	<b>Ctrl + Shift + I</b>
21	<b>Copy</b>	<b>Ctrl + C</b>
22	Switch Background image to <b>RGB Color</b>	<b>Image &gt; Mode &gt; RGB Color</b>
23	<b>Paste</b>	<b>Ctrl + V</b>
24	Repeat <b>Steps 14</b> thru <b>23</b> on remaining images	-
25	Select background image with all layers copied into it	-
26	Select layer to align	-
27	Set layer <b>Opacity</b> to 50%	-
28	Select <b>Move Tool</b>	<b>V</b>
29	Left click on layer to drag and align	-
30	Use arrow keys for fine alignment	-
31	Set layer Opacity to 100%	-
32	Repeat <b>Steps 26</b> thru <b>31</b> for remaining layers	-

33	Arrange layer order	-
34	<b>Save</b>	<b>Ctrl + S</b>

## VI. Complex Pseudo-Color – The Easy Way

	<u>Instructions</u>	<u>Shortcuts</u>
1	<b>Open</b> image	<b>Ctrl + O</b>
2	<b>Duplicate</b> image	<b>Image &gt; Duplicate</b>
3	Close Original	<b>Ctrl + F4</b>
4	Adjust <b>Levels</b>	<b>Ctrl + L</b>
5	Switch image to <b>RGB Color</b>	<b>Image &gt; Mode &gt; RGB Color</b>
6	Create <b>New Layer</b>	<b>Ctrl + Shift + N</b>
7	Choose <b>Pencil</b> or <b>Brush</b>	<b>B</b>
8	Choose <b>Color</b>	-
9	<b>Paint</b> areas to be colored	-
10	Repeat 6 thru 9 if desired	-
11	Create <b>New Layer</b>	<b>Ctrl + Shift + N</b>
12	Choose <b>Color</b>	-
13	<b>Fill</b> layer	<b>G</b>
14	Adjust <b>Opacity</b> of colored layers	-
15	<b>Save</b>	<b>Ctrl + S</b>

**VII. Calculations**



	<b><u>Instructions</u></b>	<b><u>Shortcuts</u></b>
1	<b>Open</b> Image	<b>Ctrl + O</b>
2	Make 2 <b>Duplicate</b> images & close the original	<b>Image &gt; Duplicate</b>
3	Run a 100 pixel <b>Median Filter</b> on copy_2	<b>Filter &gt; Noise &gt; Median</b>
4	Open <b>Calculations</b>	<b>Image &gt; Calculations</b>
5	Set Source 1 as copy_2 and <b>Invert</b>	-
6	Set Source 2 as the first copy	-
7	Choose <b>Blending</b> option	-
8	Set <b>Result</b> to <b>New Document</b>	-
9	Adjust <b>Levels</b>	<b>Ctrl + L</b>
10	<b>Save</b>	<b>Ctrl + S</b>



## VIII. Combining EDS Color Maps

	<u>Instructions</u>	<u>Shortcuts</u>
1	<b>Open</b> all your color map images	<b>Ctrl + O</b>
2	<b>Select All</b> of one of the images	<b>Ctrl + A</b>
3	<b>Copy</b> it	<b>Ctrl + C</b>
4	Create a <b>New Image</b>	<b>Ctrl + N</b>
5	Set Foreground Color to black	<b>D</b>
6	Use <b>Paint Bucket</b> to fill new image black	-
7	Select one of your color map images	-
8	Select <b>Magic Wand</b> and set <b>Tolerance</b> to <b>0</b>	-
9	Use <b>Magic Wand</b> to select black areas of your image	-
10	<b>Invert</b> the selection	<b>Ctrl + Shift + I</b>
11	<b>Copy</b> the selection	<b>Ctrl + C</b>
12	<b>Paste</b> it into your New Image	<b>Ctrl + P</b>
13	Repeat 7-12 on your other color map images	-
14	Arrange the layers in the New Image	-
15	<b>Save</b>	<b>Ctrl + S</b>

**IX. Multi-Detector Color**

	<b><u>Instructions</u></b>	<b><u>Shortcuts</u></b>
1	<b>Open</b> all your images	<b>Ctrl + O</b>
2	<b>Select All</b> of one of your images	<b>Ctrl + A</b>
3	<b>Copy</b> it	<b>Ctrl + C</b>
4	Create a <b>New</b> image	<b>Ctrl + N</b>
5	Make sure new image is set to <b>RGB Color</b>	<b>Image &gt; Mode &gt; RGB Color</b>
6	Set Foreground Color to black	<b>D</b>
7	Use Paint Bucket to fill new image black	-
8	<b>Paste</b> copied image into <b>New</b> image	<b>Ctrl + V</b>
9	<b>Copy</b> and <b>Paste</b> remaining images into new image	<b>Ctrl + A, Ctrl + C, Ctrl + V</b>
10	Make sure master image layer is at the bottom	-
11	Align other images to master image. (50% opacity, make 1 layer visible at a time  )	<b>V</b>
12	Make all layers visible 	-
13	Right click a layer > <b>Blending Options</b>	-
14	Uncheck 1 or 2 color <b>Channels</b>	-
15	Repeat steps 12, 13 on other layers	-
16	Drag master layer (not colored) to top layer	-
17	Adjust master layer <b>Opacity</b> (lower = more color)	-
18	<b>Save As...</b> (.JPG, .TIF, etc.) & discard layers	<b>Shift + Ctrl + S</b>
19	<b>Open</b> just saved image	<b>Ctrl + O</b>
20	<b>Levels</b>	<b>Ctrl + L</b>
21	Auto Color (optional)	<b>Shift + Ctrl + B</b>
22	<b>Save</b>	<b>Ctrl + S</b>

### X. Drawing a Mag Line Marker – Version A (Pre-Existing Mag Marker)

	<b><u>Instructions</u></b>	<b><u>Shortcuts</u></b>
1	<b>Open</b> your image	<b>Ctrl + O</b>
2	Using the existing line as a guide, draw a rectangle the same length	-
3	To change the size, scale the <b>Width</b>	<b>Edit &gt; Transform Path &gt; Scale</b>
4	Use <b>Text</b> Tool to type a label	-
5	Use <b>Move</b> tool to position line and label	<b>V</b>
6	<b>Save As</b>	<b>Shift + Ctrl + S</b>

**XI. Drawing a Mag Line Marker – Version B (NO Pre-Existing Mag Marker)****[NOT Resized & NO .GIF's]**

	<b><u>Instructions</u></b>	<b><u>Shortcuts</u></b>
1	<b>Open</b> your image	<b>Ctrl + O</b>
2	Turn on the <b>Rulers</b>	<b>Ctrl + R</b>
3	Set <b>Rulers</b> to Centimeters	Right-click on Ruler
4	Divide mag by 10,000	-
5	Using the <b>Info</b> box as a guide, draw a rectangle the length of your answer	<b>F8</b>
6	Use the <b>Text</b> Tool to type a label for the line	-
7	Use the <b>Move</b> Tool to position your line and label	<b>V</b>
8	<b>Save As</b>	<b>Shift + Ctrl + S</b>